


MILTEX HI-LIGHT

ENDO MEASURING BLOCK

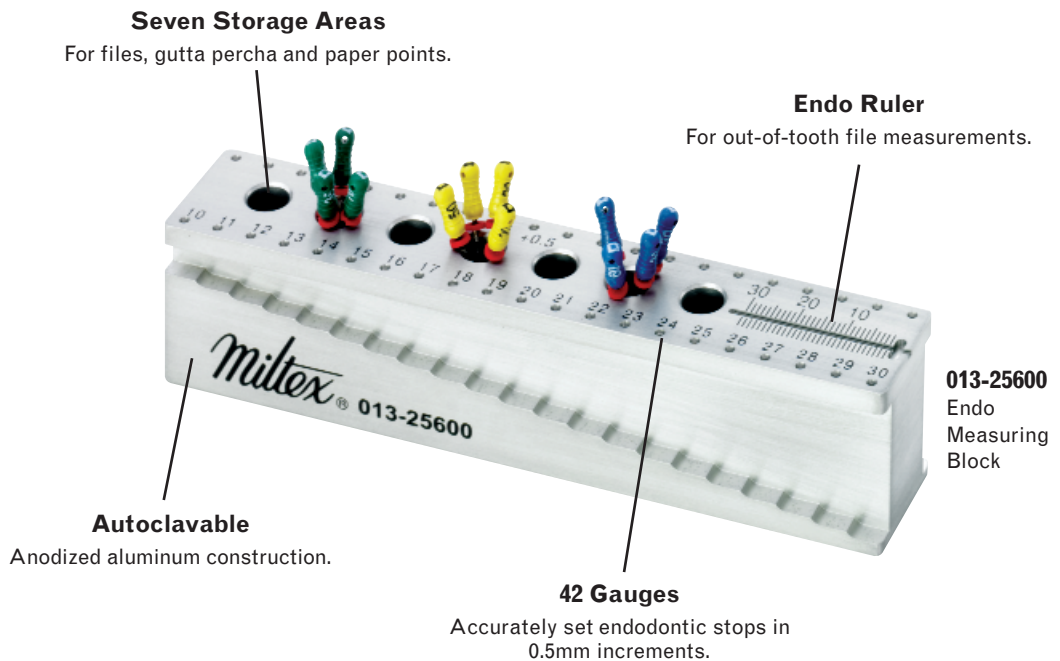

Designed to accurately set stops for all endodontic files, the Endo Measuring Block from Miltex, Inc. is made from anodized aluminum and is autoclavable by any method. The 42 gauges allow precise placement of stops from 10mm to 30mm in 0.5mm increments and can also be used to size gutta percha and paper points. The seven storage areas are ideal for storing files, gutta percha and paper points, while the endo ruler allows accurate out-of-tooth file measurements.

See the Miltex Dental catalog for a complete list of endodontic files, gutta percha, and paper points.

For additional information, an in-office demonstration, or to place an order, please contact your local Miltex, Inc. authorized Distributor or Miltex, Inc. Customer Service at 1-866-854-8300, email customerservice@miltex.com, or visit us at www.miltex.com.

Miltex
Redefining Excellence

589 Davies Drive
York, PA 17402
phone 717 840-9335
toll free 866 854-8300
fax 866 854-8400
www.miltex.com

0304/2/A/1

DENT-ENDOBLK-HI



+H834DENTENDOBLKHC