

Care for  
yourself  
while you  
care for  
others.



**3M™ VFlex™ Health Care  
Particulate Respirator  
and Surgical Masks  
1804/1804S (small)**


Whether you are in the OR or on the floor, you put yourself at risk to care for others every day. The VFlex™ 1804/1804S offers the protection of a NIOSH-approved surgical N95 respirator along with the breathability you desire from a traditional surgical mask.\*

\*See back page for details

- Pleats help increase surface area for ease of breathing
- Respirator expands for a spacious feel
- Tabs designed to help achieve a comfortable seal
- Available in two sizes
- Fits a broad range of faces

# 3M™ VFlex™ Health Care Particulate Respirator and Surgical Masks, 1804/1804S

**Surgical N95 Respirator Differential Pressure (Delta-P) Study**

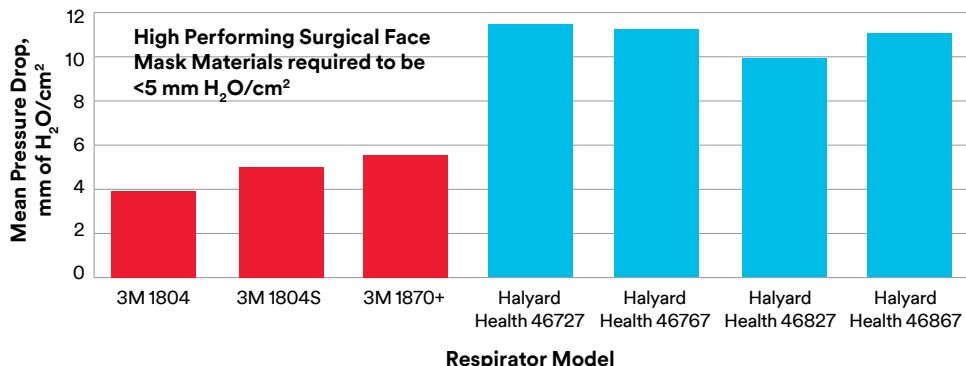

Compared to the Halyard Health models shown, VFlex™ 1804/1804S is a NIOSH-approved surgical N95 respirator that more closely approaches the material breathability of a traditional surgical mask.†\*

†Testing per MIL-SPEC 36945C 4.4.1.2 using 4.9cm<sup>2</sup> sample at 8L/min. Source: Representative pressure drop performance from internal testing of multiple lots of all models, November 2017. Selection of respirator models, testing protocol, data generation and conclusions were reviewed and approved by an expert from the University of Illinois at Chicago. Statistically different per T-tests (P-value <0.0001).

\*Performance requirement ASTM F2100-11 medical facemask materials classification.

## Additional design features

- Respirator expands for a spacious feel
- Folded shaped design helps keep respirator away from your mouth
- Retains comfortable shape

### Tabs

- Uniquely designed for positioning the respirator on the face for convenience and to help achieve a comfortable seal

### Comfortable inner layer

### 3M proprietary filter media

- Proprietary filter media is designed for ease of breathing

## Adjustable noseclip

- Helps provide a custom and secure seal
- Embedded, metal noseclip

## Elastic headbands

- Not made with natural rubber latex

## Pleats

- Flex with mouth movement while talking
- Help increase surface area for ease of breathing

NIOSH-Approved N95  
Fluid resistant • FDA cleared

## Ordering Information

Cat. Number	Description	Item/Box	Box/Case
1804	3M™ VFlex™ Health Care Particulate Respirator and Surgical Mask	50	8
1804S	3M™ VFlex™ Health Care Particulate Respirator and Surgical Mask, Small	50	8

For additional information on all 3M™ Health Care Particulate Respirator and Surgical Masks, visit [3M.com/healthcaremasks](https://www.3M.com/healthcaremasks)



**3M Health Care**  
**Infection Prevention Division**  
 2510 Conway Avenue  
 St. Paul, MN 55144-1000  
 USA  
 1-800-228-3957  
[3M.com/medical](https://www.3M.com/medical)



## WARNING

These respirators help protect against certain airborne particulates. Before use, the wearer must read and understand the User Instructions provided as a part of the product packaging. A written respiratory protection program must be implemented meeting all the requirements of OSHA 29 CFR 1910.134 including training, fit testing and medical evaluation. In Canada, CSA standards Z94.4 requirements must be met and/or requirements of the applicable jurisdiction, as appropriate. Misuse may result in sickness or death. For correct use, see supervisor and package instructions, or call 3M Health Care Helpline in USA at 1-800-228-3957 and in Canada at 1-800-563-2921.

3M and VFlex are trademarks of 3M.  
 Please recycle. Printed in U.S.A.  
 © 3M 2018. All rights reserved.  
 70-0716-9747-1