


## Characteristics

COATS™ barrier technology is a totally powder free, all-natural physical barrier between the glove surface and the wearer's skin, utilizing the anti-irritant and anti-inflammatory properties of colloidal oatmeal. In our COATS™ gloves, the colloidal oatmeal barrier not only shields the skin surface, it also creates a moisturizing and pH balancing environment that helps dry, damaged skin to repair itself. And healthy skin is the body's best defense against infection. What a naturally great idea!

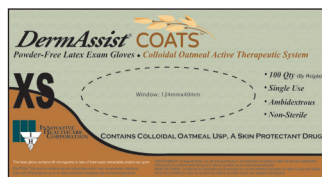

**Exam Glove  
Non-Sterile**

**DermAssist®  
COATS™**

**Latex  
Series 124**



## PRODUCT DETAILS

SIZE	ITEM NO.	PACKAGING	DESCRIPTION
XS	124050	100 Gloves/box, 10 boxes/case	Gloves, Exam, Latex, Non Sterile, Powder-Free, Colloidal Oatmeal, Therapeutic
S	124100	100 Gloves/box, 10 boxes/case	
M	124200	100 Gloves/box, 10 boxes/case	
L	124300	100 Gloves/box, 10 boxes/case	
XL	124350	100 Gloves/box, 10 boxes/case	


## COATS™ - Colloidal Oatmeal Active Therapeutic System

**- A glove that can help protect you and your skin.**

- COATS™ is an innovative new barrier coating comprised of all natural colloidal oatmeal USP which is an FDA recognized skin protectant drug.
- Applied to the don side of the glove through a proprietary process, COATS™ is a powder-free process.
- COATS™ technology reduces irritation and hydration dermatitis by creating a physical barrier between glove film and skin; absorbing perspiration, urea, and salt excretions; providing anti-irritant and anti-inflammatory benefits; effectively pH buffering the skin; and naturally moisturizing the surface of the skin.

*Product Solutions You Trust*

## *DermAssist*® **COATS**™ Powder-Free Latex Exam Gloves • Colloidal Oatmeal Active Therapeutic System

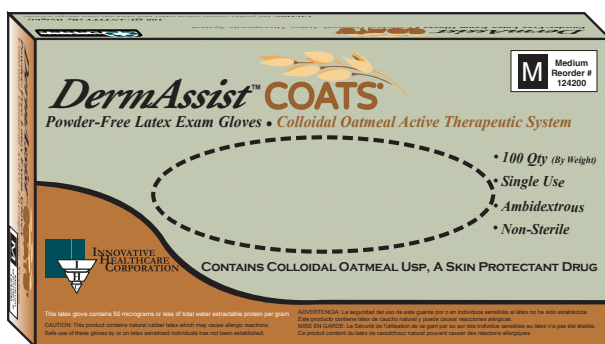

**CONTAINS COLLOIDAL OATMEAL USP, A SKIN PROTECTANT DRUG**

### Product Attributes

- Textured
- Low Modulus
- Low Protein
- Non-Detectable Residual Chemical Level

### Benefits

- Improved Wet/Dry Grip
- Soft, More Comfortable Fit
- Reduced Chances of Latex Allergic Reactions
- Reduced Chance of Dermatitis and Type IV Reactions

*DermAssist® COATS™ is manufactured in compliance with multiple international standards, including the following:*

<b>Designation</b>	<b>Standard</b>
<b>ASTM D3578</b>	Standard Specification for Rubber Examination Gloves
<b>ASTM D5151</b>	Standard Test Method for Detection of Holes in Medical Gloves
<b>ASTM F1671</b>	Standard Test Method for Resistance of Materials Used in Protective Clothing to Penetration by Blood-Borne Pathogens
<b>ASTM D5712</b>	Standard Test Method for Analysis of Aqueous Extractable Protein in Natural Rubber

<b>Average Length</b>	<b>Average Palm Thickness</b>	<b>Average Finger Thickness</b>
9.5 in ± 240 mm	4.0 mil ± 0.10 mm	4.5 mil ± 0.12 mm

<b>Tensile Strength &amp; Elongation</b>	<b>Before Aging</b>	<b>After Accelerated Aging</b>
Tensile Strength (Mpa)	28	24
ASTM Requirement Min. (Mpa)	18	14
Elongation (%)	850	820
ASTM Requirement Min. (%)	650	500



**Intertek**

Innovative Healthcare Corporation is certified to both ISO 9001:2008 QMS and ISO 13485:2003 QMS for medical devices.



**Intertek**



**INNOVATIVE  
HEALTHCARE  
CORPORATION**

*Product Solutions You Trust*

**1-800-272-1533**