



SURE-CHECK® Sterilization Pouches

Giving You 3x The Confidence!



The only way to be SURE about sterilization.

CDC infection control guidelines make it clear: the only way to be sure an instrument has been sterilized is to expose it to **all 3** sterilization criteria: **time, temperature and steam**. Yet traditional sterilization pouches have an internal and/or external process indicator that reacts only to temperature.

Sure-Check® Sterilization Pouches, feature internal/external multi-variable indicator technology that confirms all 3 criteria for sterilization have been met.

For CDC compliance, today more than ever, you need to be sure about sterilization – you need Sure-Check!

Sure-Check™ Sterilization Pouches




Sure-Check®:

- Meets CDC Guidelines for Disinfection and Sterilization in Healthcare Facilities, 2008.
- Tested and certified by an independent laboratory.
- When used in conjunction with weekly biological monitoring, provides the highest level of sterility assurance, unlike traditional sterilization pouches with process indicators that measure only temperature.*
- Meets ANSI/AAMI ST79 Comprehensive Guide to Steam Sterilization and Assurance in Health Care Facilities.
- Eliminates the extra step of adding separate multi-variable internal indicators to each pouch.
- US Food and Drug Administration cleared.

* Multi-Variable Indicators: a multi-variable indicator shall be designed to react to two or more of the critical variables and is intended to indicate exposure to a sterilization cycle at a stated value of the chosen variable.

Sure-Check Sterilization Pouches are not considered hazardous waste and can be disposed of in regular trash (40 CFR 261.24 tested using the TCLP, EPA Publication SW-846).

Patent pending

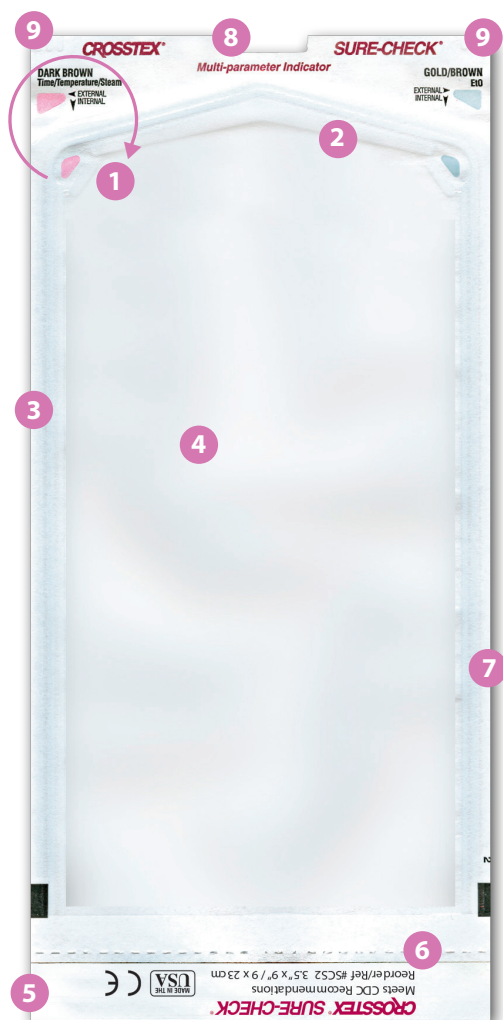

Because You Can't See Sterile, Be Sure with Sure-Check:

- 1 **Internal/External Multi-Variable Indicators** - only change color when all 3 criteria for sterilization have been met: **time, temperature and steam.**
- 2 **Secure Seal** - Extra wide chevron to protect against instrument breaching.
- 3 **Medical Grade Paper** - Wet strength treated, virgin medical grade paper that meets medical pouch standards.
- 4 **Transparent Film** - Allows a clear view of the contents as well as the internal indicators.
- 5 **Wide Self-Seal Adhesive Strip** - Wider sealing area of strip securely seals pouch and provides sterilization integrity of contents during storage. Contents will remain sterile for three years provided the integrity of the pouch is not compromised.
- 6 **Perforated Fold** - Allows for accurate folding of adhesive strip to form an airtight seal, maintaining the integrity of the pouch during storage after sterilization.
- 7 **Complete Film/Paper Separation** - The high quality of the materials and manufacturing process assures a consistently clean separation (peeling) of the film from the paper, maintaining the sterile integrity of the contents.
- 8 **Thumb Notch** - Assures easy opening of pouch.
- 9 **Corner Tack Seals** - Three tack seals on each corner prevent curling and possible separation after sterilization. In addition they prevent the collection of dust and other contaminants in the sealed area.

Ordering Information



AVAILABLE IN 14 SIZES!

Description	Quantity	Ref. #
2.25" x 4" / 6 x 10 cm	200/Box	SCXX2
2.75" x 9" / 7 x 23 cm	200/Box	SCX2
3.5" x 5.25" / 9 x 13 cm	200/Box	SCXS2
3.5" x 9" / 9 x 23 cm	200/Box	SCS2
3.5" x 9" / 9 x 23 cm	500/Box	SCS25
3.5" x 22" / 9 x 56 cm	100/Box	SCL3222
5.25" x 6.5" / 13 x 17 cm	200/Box	SCW2
5.25" x 10" / 13 x 25 cm	200/Box	SCM2
5.25" x 15" / 13 x 38 cm	200/Box	SCL1512
7.5" x 13" / 19 x 33 cm	200/Box	SCL2
8" x 16" / 20 x 41 cm	200/Box	SCL8162
10" x 15" / 25 x 38 cm	100/Box	SCL10152
10.5" x 17" / 27 x 43 cm	100/Box	SCL10172
12" x 15" / 30 x 38 cm	100/Box	SCL12152
12" x 18" / 30 x 46 cm	100/Box	SCL12182

For Questions Relating to Sure-Check® Sterilization Pouches
Please Contact Crosstex Customer Care @ 800-819-3336 or
Email customer@crossstex.com

CROSSTEX
A CANTEL MEDICAL COMPANY

PROTECTS

Manufactured by Crosstex
10 Ranick Road
Hauppauge, NY 11788-4209 USA

Tel: 631-582-6777 | 888-276-7783
Fax: 631-582-1726 | Email: crosstex@crossstex.com

crosstex.com

California • Colorado • Georgia • New York • Pennsylvania • Argentina • Japan • The Netherlands

**Visit CrosstexBMS.com
for Educational Videos**

