

Miltex[®]
Redefining Excellence



DENTAL SYRINGES AND ACCESSORIES

GripRite™ Syringes

REF	Description
76-30	Syringe – Aspirating CW , Blue Silicone Grips
76-31	Syringe – Aspirating CW , Pink Silicone Grips
76-32	Syringe – Aspirating CW , Red Silicone Grips
76-33	Syringe – Aspirating CW , Purple Silicone Grips
76-40	Syringe – Petite Aspirating CW , Blue Silicone Grips
76-41	Syringe – Petite Aspirating CW , Pink Silicone Grips
76-42	Syringe – Petite Aspirating CW , Red Silicone Grips
76-43	Syringe – Petite Aspirating CW , Purple Silicone Grips



76-30



76-40



76-31



76-41



76-32



76-42



76-33



76-43

Silicone covered handles and wings offer a comfortable grip on syringe.

Key Product Features & Benefits

Features

- Silicone covered handles & wings
- 4 color options
- Available in standard & petite styles

Benefits

- Comfortable grip on syringe
- Ability to color code
- Petite style's shortened harpoon rod facilitates aspiration for smaller hands



Aspirating Syringes


Miltex offers a wide variety of aspirating syringes to accommodate the way you practice.



REF	Description
76-70	Aspirating Syringe, C-W Type
76-80	Aspirating Syringe, Type "A"
76-105	Articulating Barrel Syringe, 1.8 cc
77-70	Petite Aspirating Syringe, C-W Type



76-70

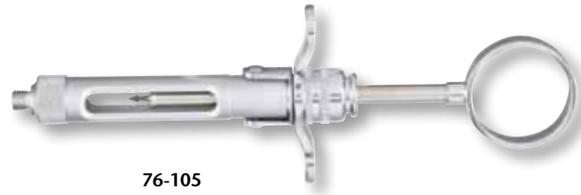

76-105



76-80



77-70

Self-Aspirating Syringes

Self-Aspirating Syringes feature a blunt harpoon rod that allows the syringe to aspirate by means of counter pressure.

REF	Description
76-90	Lightweight Self-Aspirating Syringe, 1.8 cc
76-100	Self-Aspirating Syringe, 1.8 cc



76-90



76-100

Petite Syringe

The Miltex® Petite Syringe was specifically designed to meet the needs of dental professionals who prefer a more comfortable, form-fitting syringe. Ergonomic by design, the Miltex® Petite Syringe features:

- *Reduced thumb ring handle diameter to eliminate the excess gap, allowing just enough room for comfort and ease.*
- *New oval shaped ring handle that fits naturally to the shape of your thumb.*
- *Rounded edges on ring handle to create a smooth and soft feel.*
- *Slightly flared base on ring handle to ease pressure while dispensing.*
- *Shortened harpoon rod to reduce strain during aspiration.*
- *Wide wings allow for a better grip with greater control.*



REF	Description
77-70	Petite Aspirating Syringe, Cook-Waite (CW)

FACT:

Most Miltex® syringe hubs are permanently affixed to the syringe. The hubs do not need to be removed for sterilization, preventing them from being misplaced.

Articulating Barrel Syringe

The Articulating Barrel Syringe features a short barrel and a pivoting barrel that helps facilitate easy loading of cartridges. This pivoting barrel is especially helpful during multiple cartridge injections because the cartridges are loaded and unloaded by placing them directly into the barrel, preventing needles from bending and eliminating the use of multiple needles. The shortened barrel is great for anterior injections and hard to access areas of the mouth. The Articulating Barrel Syringe has a removable hub which allows for the use of the 30° adapter hub for intraligamentary injections.



REF	Description
76-105	Articulating Barrel Syringe 1.8cc
76-56	30 degree adapter hub

N-Tralig® Syringe

The Miltex® N-Tralig® Syringe is one of the only stainless steel syringes on the market that delivers precisely 0.2 cc of anesthetic per ratchet. The positive ratchet action helps protect against the breakage of anesthesia cartridges. The Miltex® N-Tralig® syringe also features a removable hub which allows the use of the 30° adapter hub to help with the most difficult to reach injections. Miltex® N-Tralig® syringe is available with one solid barrel and one window barrel.

REF	Description
76-55	N-Tralig Intraligamentary Syringe
76-56	30 degree adapter for N-Tralig Syringe
PN2	Cartridge barrel solid wall for N-Tralig Syringe
PN3	Cartridge barrel shielded windows for N-Tralig Syringe
PN4	Barrel cap for N-Tralig Syringe



76-55



76-56



PN2



PN4



PN3



www.miltex.com
toll-free phone 866-854-8300
toll-free fax 866-854-8400
phone 717-840-9335
fax 717-840-9347

Please call Customer Service today for additional information, an in-office demonstration, to place an order, or to receive other Miltex catalogs.

Miltex, Miltex Redefining Excellence and N-Trilig are registered trademarks of Miltex, Inc., a subsidiary of Integra LifeSciences Corporation.
GripRite is a trademark of Miltex, Inc., a subsidiary of Integra LifeSciences Corporation.
The Integra wave logo is a trademark of Integra LifeSciences Corporation.
Copyright© 2010 Integra LifeSciences Corporation

