



S-SERIES ULTRASONIC SCALING TIPS

PIEZO TECHNOLOGY



SCALING



S10 UNIVERSAL
USS10

Standard diameter tip used to remove heavy supragingival deposits.

POWER LEVEL MED—HIGH



P UNIVERSAL
USP

Standard diameter tip used to remove moderate to heavy deposits located supragingivally and in shallow pockets.

POWER LEVEL MED—HIGH



P10 UNIVERSAL
USP10

Standard diameter tip used to remove light deposits and plaque. Tip features smooth edges for safe removal of soft and hard deposits.

POWER LEVEL MED



100 UNIVERSAL
US100

Probe-sized tip used to remove light deposits in shallow to moderate pockets and for periodontal maintenance.

POWER LEVEL LOW—MED

AFTER FIVE®

These probe sized tips are designed to enhance access into periodontal pockets and improve adaptation in interproximal areas, allowing the clinician to more thoroughly debride difficult to reach areas to achieve better clinical results.



STRAIGHT
USAFS

POWER LEVEL MED



RIGHT
USAFR

POWER LEVEL MED



LEFT
USAFL

POWER LEVEL MED

S-SERIES TIPS are compatible with Hu-Friedy Symmetry, Satelec/Acteon and NSK Piezo Generators.

How the best perform



PERIO



UNIVERSAL PERIO USU

Thin tip used to remove light deposits and plaque. Tip features smooth edges for safe removal of soft and hard deposits.


POWER LEVEL MED



LEFT PERIO US4L

Standard diameter contra-angle tip used to remove calculus on all surfaces of premolars and molars, especially difficult to reach interproximal surfaces.

POWER LEVEL LOW

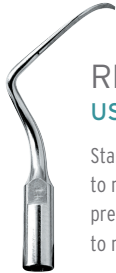

RIGHT PERIO US4R

Standard diameter contra-angle tip used to remove calculus on all surfaces of premolars and molars, especially difficult to reach interproximal surfaces.

LOW




3 PERIO USP3

Thin-bladed curette-like tip used to remove heavy deposits and necrotic cementum on anteriors and premolars.

(Also available for Hu-Friedy EMS devices – UEP3)

LOW

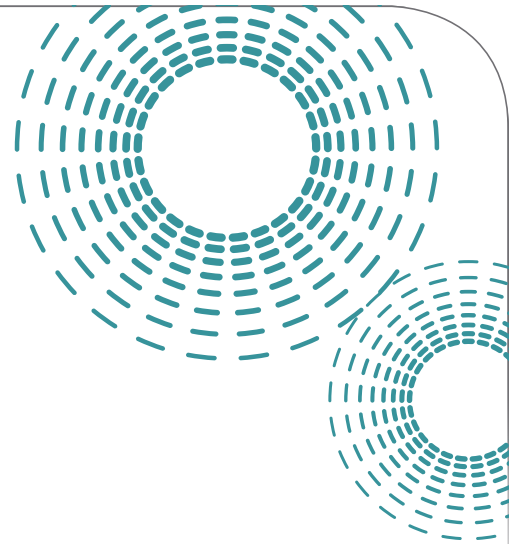

DIAMOND-COATED



LEFT DIAMOND US2LD

Thin, contra-angle diamond-coated finishing instrument used for removal of residual debris and designed to enhance access to deep periodontal pockets and furcations.

POWER LEVEL LOW



RIGHT DIAMOND US2RD

Thin, contra-angle diamond-coated finishing instrument used for removal of residual debris and designed to enhance access to deep periodontal pockets and furcations.

LOW



STRAIGHT DIAMOND US1D

Thin, straight diamond-coated finishing instrument used for removal of residual debris and designed to enhance access to deep periodontal pockets primarily in anterior areas.

LOW


GUARDIAN WRENCH/TIP CARRIERS

S-Series Piezo Tips come with the Guardian™ Wrench Tip Carrier—an innovative, patented torque-limiting wrench/tip carrier designed for remarkably safe and easy connection to the handpiece.

USAGE:

- Tips contained in an ORANGE wrench should be used primarily on medium to high power.
- Tips contained in a BLUE wrench should be used primarily on low to medium power (in Scaling mode if applicable).
- Tips contained in a GREEN wrench should be used primarily on low power only (in Perio mode if applicable).

Color of the Guardian indicates equipment application mode (if applicable). Refer to your equipment for Scaling and Perio application mode.

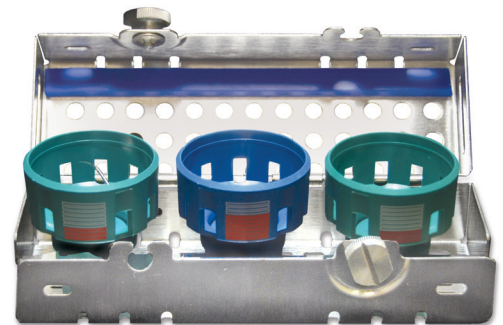

Visit us online at Hu-Friedy.com

Call 1-800-Hu-Friedy or contact your authorized Hu-Friedy representative for more information.

©2017 Hu-Friedy Mfg. Co., LLC. Satelec, NSK, and EMS are trademarks of Acteon, Inc., Nakanishi Inc., and Electro Medical Systems respectively. All rights reserved. HF-507/0817.

