




SEPTOJECT®

SEPTOJECT NEEDLES
SEPTOJECT XL
SEPTOJECT EVOLUTION
ASPIRATING SYRINGES



If you think a needle is just a needle, think again.

**YOUR PATIENTS FEAR NEEDLES.
SHOULDN'T YOU MINIMIZE THEIR
DISCOMFORT BY CHOOSING SEPTODONT?**

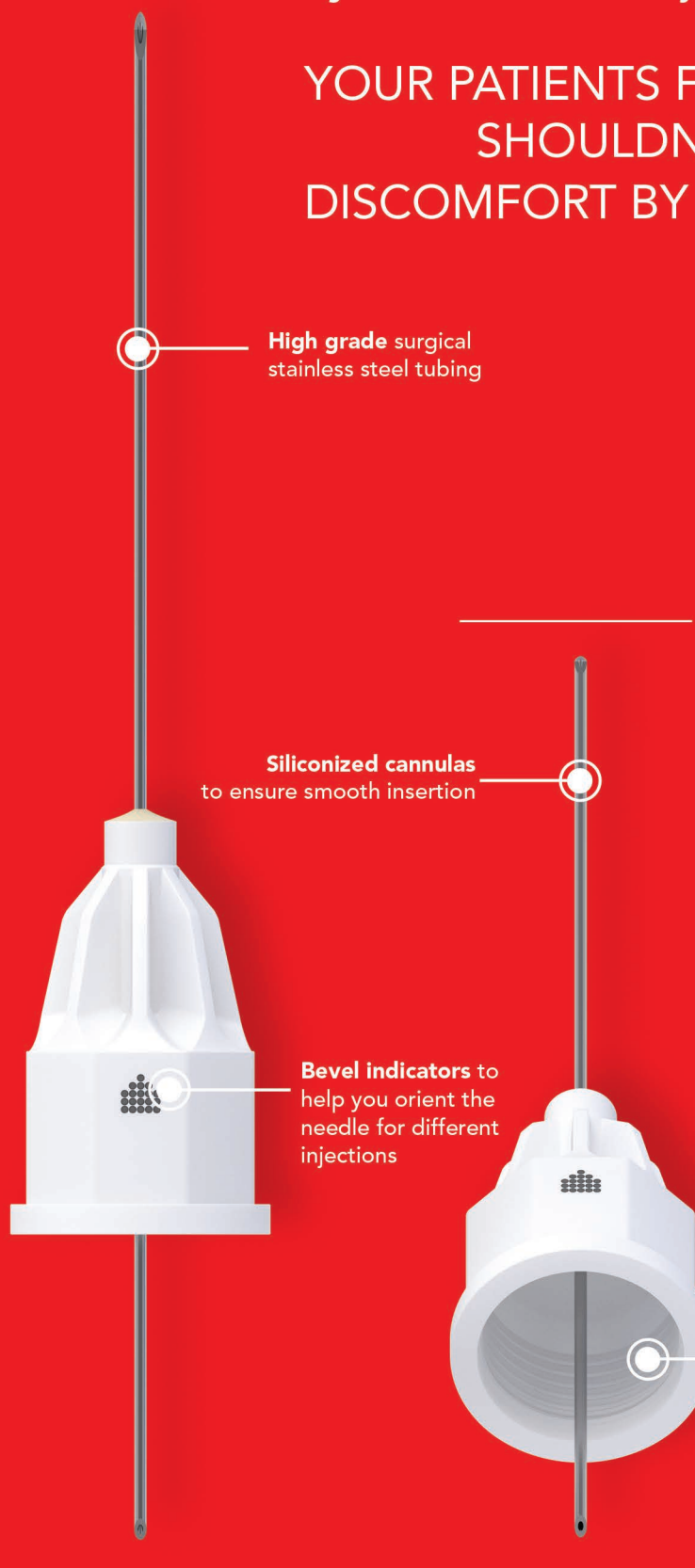

No one offers you more choices.

All our needles offer:

- Pre-threaded plastic hubs for easier threading and no metal to metal contact
- Siliconized cannulas to ensure smooth insertion
- Bevel indicators to help you orient the needle for different injections
- High grade surgical stainless steel tubing

We are the only manufacturer offering you three options:


- Septoject Needles – sharpness you can see
- Septoject XL – less pressure for the same flow
- Septoject Evolution – greater control, less discomfort

SEPTOJECT NEEDLES

TRIPLE BEVEL ▶

PLASTIC HUB — EASIER THREADING

Our **Septoject Standard needles** have a triple bevel design. Many other needles on the market offer only a double bevel design. The added bevel provides for less tissue displacement and sharpness you can see.

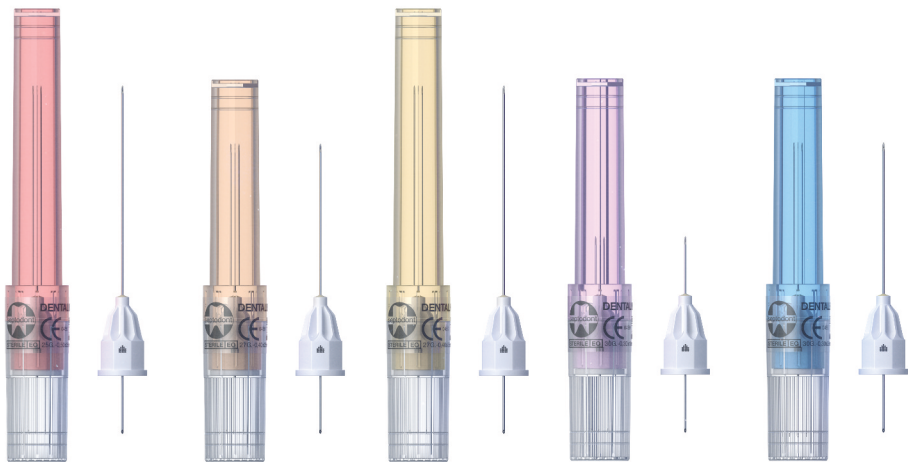

Septoject Standard

Item No.			
01-N1252	25 GAUGE LONG	35MM	●
01-N1271	25 GAUGE SHORT	25MM	●
01-N1272	27 GAUGE LONG	35MM	●
01-N1300	30 GAUGE EXTRA SHORT	10MM	●
01-N1301	30 GAUGE SHORT	25MM	●

— All Gauges: 100 Pre-threaded Plastic Hub Needles per Box. 10 Boxes per Case.


SEPTODONT ASPIRATING SYRINGES

Septodont dental aspirating syringes are lightweight, durable and contain no removable parts. Best of all, they are available in three distinct sizes to fit your specific injection style. They offer lightweight aluminum construction of thumb rings, wings, and a surgical grade stainless-steel body.


Septodont Standard (gold) syringe is equipped with a larger thumb ring and offers a standard length plunger for improved comfort.



01-N2110

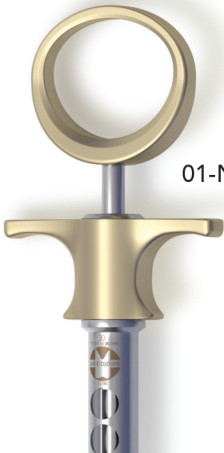

Septodont Fusion (titanium) has a medium thumb ring and a standard length plunger.

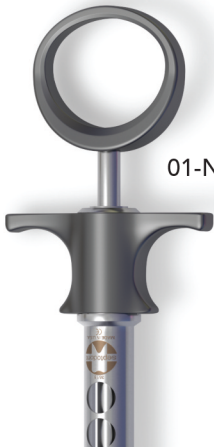

01-N2105

Septodont Petite (blue) syringe is equipped with a compact design, finger grip and a standard length plunger for better control and accuracy.




SEPTOJECT XL NEEDLES

LARGER LUMEN ►

SIMPLIFY YOUR CHOICES — BLOCK OR INFILTRATION

Our **Septoject XL needles** offer a large bore design – 43% wider than a standard needle. It requires less pressure for injection.

Choose block or infiltration and simplify your selection.

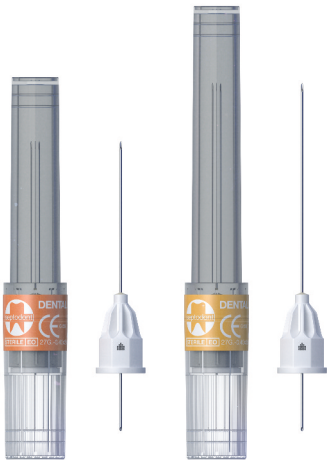

Septoject XL



Item No.			
01-N1500	SEPTOJECT XL - INFILTRATION	25MM	●
01-N1550	SEPTOJECT XL - BLOCK	35MM	●

All Gauges: 100 Pre-threaded Plastic Hub Needles per Box. 10 Boxes per Case.

SEPTODONT ECO SYRINGES

Our Eco syringes share the same basic design of our Septodont syringes but replace the aluminum thumb ring and wings with stainless steel making them more cost-effective.

ite is designed for small hands.
sign provides a smaller thumb ring,
a shorter plunger length. This provides
nd added comfort when aspirating.

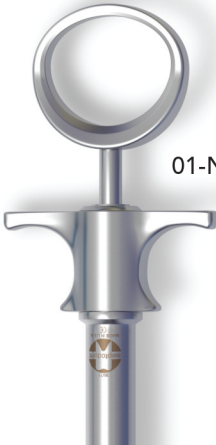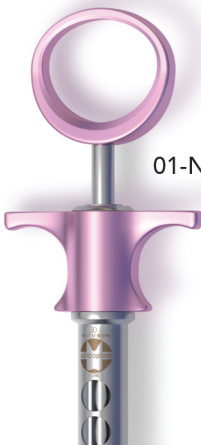

01-N2100

01-N2120

01-N2200


01-N2210

SEPTOJECT EVOLUTION NEEDLES

GREATER CONTROL — REDUCED DISCOMFORT


Our **Septoject Evolution needle** has an innovative patented scalpel-designed bevel. It is engineered for a smoother penetration with less tissue displacement. It requires less force, provides reduced deflection and brings you better control and accuracy. Not intended for block injections.

REDUCED DEFLECTION BETTER ACCURACY

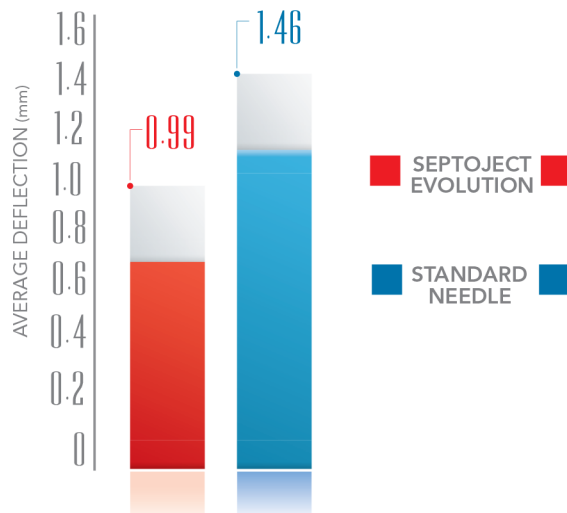

Comparative Deflection (mm) of Needle After Penetration Through a 10mm Block of Silicone²

EASIER INSERTION - 29% LESS FORCE REQUIRED

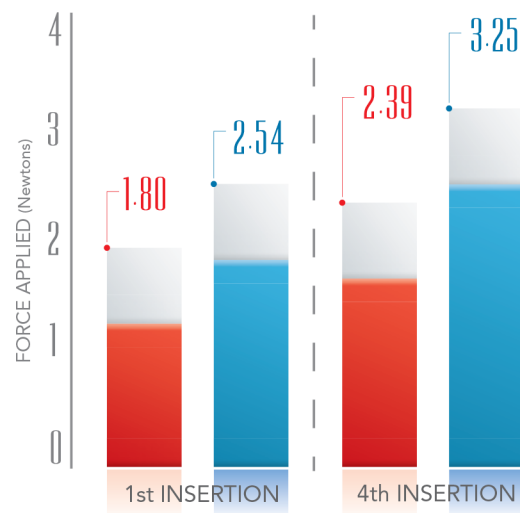

Force required (Newtons) for a 10mm Penetration of a Silicone Block¹



INFILTRATIONS

- Infiltrations (periapical)
- PSA - Posterior Superior Alveolar nerve block
- MSA - Middle Superior Alveolar nerve block
- Palatal Anesthesia



INTRALIGAMENTARY INJECTIONS

- Intraligamentary Injections
- Intraseptal Injections

Septoject Evolution

Item No.				
01N1600	25G	SHORT	25MM	●
01N1630	30G	SHORT	25MM	●

All Gauges: 100 Pre-threaded Plastic Hub Needles per Box. 10 Boxes per Case.

Septoject Evolution

Item No.			
01N1650	30G X-SORT	9MM	●



205 GRANITE RUN DRIVE
SUITE 150
LANCASTER, PA 17601

800.872.8305

www.septodontusa.com

SEPTOJECT NEEDLES
SEPTOJECT XL
SEPTOJECT EVOLUTION
ASPIRATING SYRINGES

- | ANESTHETICS
- | PAIN CONTROL
- | ENDODONTICS
- | RESTORATIVES
- | INFECTION CONTROL
- | PERIO & SURGICAL
- | PROSTHETIC