



FlowTM

DENTAL

Product Catalog

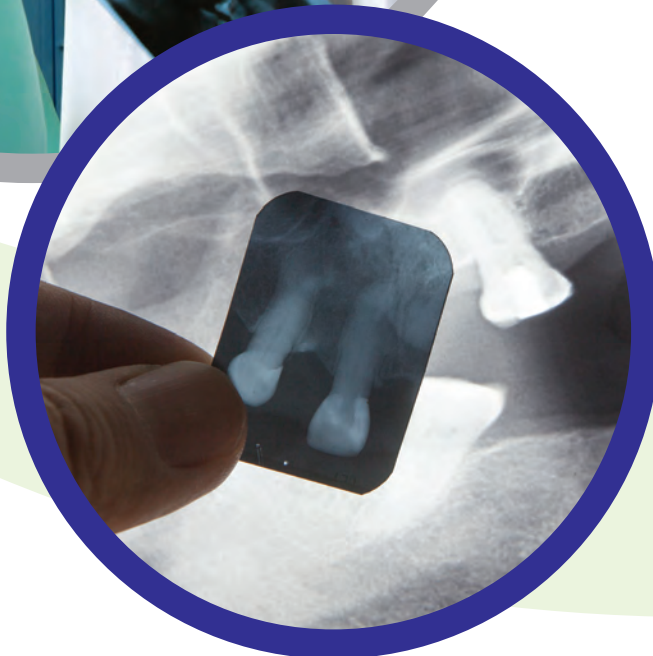


Table of Contents

A

Aprons	22-23
Apron Hangers	23
Auto ID Printer	27
Automatic Processor - Developer, Fixer.....	28
Automatic Processor - Cleaner	28

B

Bib Apron	22-23
Bite Block Baggies	4, 15, 20
Bite Chex DX-42, 100 pieces, Size 3 Film	12
Bite Wing Loops "Loop EZ"	5, 14
Bite Wing Tabs, Self Adhesive "Stickems" ...	5, 14
Butterfly D Speed, Size 2, Single Films	12
Butterfly D Speed, Size 2, Double Films	12
Butterfly D Speed, Size 0, Single Films	12
Butterfly D Speed, Size 0, Double Films	12
Butterfly F Speed, Size 2, Single Films	12
Butterfly F Speed, Size 2, Double Films	12
Butterfly F Speed, Size 0, Single Films.....	12
Butterfly F Speed, Size 0, Double Films	12

C

Cable Saver	20
Cards, Printer Identification	27
Cassettes	6, 13
Cephalometric Tracing Box	25
Chemistry - Developer, Fixer	28
Child's Protective Apron	22-23
Cleaner, Processor	28
Clean Ups PSP Cleaning Wipes.....	3
Clearing Sheets - Processor Cleaner	28
Comfees	21
Cushies Film Cushions	15

D

Darkroom Safelights	26
Darkroom Thermometers	28
Darkroom Timers	27
Deluxe Safe 'N' Sure	2
Dental Film	8-13
Dental Film Hangers	28
Dental Film Hanger Identification Tabs	28
Dental Film Mounts.....	16-18
Dental Illuminators	25
Developer, Liquid	28
Digi Loops Digital Sensor Holders	20
Digital Timer	27
Dosimeter Badge Service	24
Duplicating Film	13
Duplicators	26
DV-58 Econo Pak 600, Size 2 Films	9
DV-58 Value Pak 150, Size 2 Films	9
DV-58 Fresh Pak 100, Size 2 Films	9
DV-57 Value Pak 150, Size 2, Double Films.....	9
DV-57 Fresh Pak 100, Size 2, Double Films.....	9
DV-54 Pedo Pak 100, Size 0 Films	9
DX-58 Paper Pak 150, Size, Double Films	12
DX-54 Paper Pak 100, Size 0 Films	12
DX-42 Bite Chex 100 Size 3 Films	12

E

Econo Comfees	21
Econo Pak, D Speed 600 Films	9
Econo Pak, F Speed 600 Films	11
Econo Safe 'N' Sure	2
Endodontic X-Ray Positioner	7
Exposure Articulating paper Strips.....	4, 15

F

Film Duplicators	26
Film Hangers, Dental	28
Film Mounts	16-18
Film, Panoramic	13
Fixer, Liquid	28
Fresh Pak, D Speed 100 Films	9
Fresh Pak, F Speed 100 Films	11
FV-58 Econo Pak 600, Size 2 Films	11
FV-58 Value Pak 150, Size 2 Films	11
FV-58 Fresh Pak 100, Size 2 Films	11
FV-57 Value Pak 150, Size 2, Double Films	11
FV-57 Fresh Pak 100, Size 2 Films	11
FV-54 Pedo Pak 100, Size 0 Films	11
FV-53 Pedo Pak 100, Size 0, Double Films	11

G

Greylight Mounts	16, 18
------------------------	--------

H

Hangers, Apron	23
Hangers, Dental Film	28
Hangers, Film Identification Tabs	28

I

Identification Printer	27
Identification Printer Cards.....	27
Identification Tabs, Dental Hanger	28
Illuminators	25
Index, Numerical	29
Intensifying Screens.....	25

K

Klear-View Cardboard Mounts	16, 18
Kopy Kat Film Duplicator	26

L

Lead Free Protective Aprons	23
Loops, Bite-Wing "Loop-EZ"	5, 14

M

Mounts:	
Greylight Mounts	16, 18
Grip Right Plastic Mounts	16, 18
Klear-View Cardboard Mounts	16, 18
Perfect Pocket Mounts	16, 17, 18
Tear Off Perfect Pocket Mounts	17, 18

O

OPT, Safe 'N' Sure	3
--------------------------	---

P

Panoramic Film	13
Paper Pak, DX-58 150, Size 2 Films	12
Paper Pak, DX-54 100, Size 0 Films	12

Perfect Pocket Mounts	16, 17, 18
Plastic Illuminator - Slimline	25

Positioners:

Cushies	15
Digi Loops Sensor Holders	20
Loop-EZ Bite Wing Loops	5, 14
RAPiD Positioning Device	7
Sensibles Universal Sensor Holders	19
Stickems Self Adhesive Bite-Wing Tabs	5, 14
Sticky Bites	5, 14
SUPA Disposable Bite Blocks	5, 14
Printer, Identification	27
Processor Solutions	28
Processor, Cleaner - Automatic	28
Protective Aprons	22-23

R

Radiation Badge Service.....	24
RAPiD Positioning Device	7
RAPiD Shorty	7

S

Safelights	26
Safe 'N' Sure	2-3
Sensibles.....	19
Self-Adhesive Bite Wing Tabs "Stickems"	5, 14
Silver D, D Speed Film	8-9
Cassettes	13
Slimline Illuminator	25
SMART Phosphor Plate Cassettes	6
SMART Positioning System	6
Snuggies Adjustable Sensor Sleeves	21
Stickems	5, 14
Sticky Bites	5, 14
SUPA	5, 14

T

Tabs, Identification	28
Tabs, Self Adhesive Bite Wing "Stickems"	5, 14
Tear Off "Perfect Pocket" Mounts	17
Thermometers, Darkroom	28
Thyroid Collar	22-23
Timers, Darkroom	27
Tracing Box, Cephalometric	25

X

Xpress F Speed Film	11-12
X-Ray Duplicator	26
X-Ray Film, Dental	8-13



For over 40 years Dental Professionals have relied on Flow products. Our rich history in dental film and film accessories has led us to develop many unique imaging products for today's phosphor plate and sensor based practitioners.

Regardless of your choice of media; film, phosphor plates, or sensors, Flow Dental is your source for innovative accessories that help to make your practice run more efficiently, safely and economically!

This catalog is divided into sections devoted to each imaging media for your added convenience.

For additional information please log on to FlowDental.com. Thanks for choosing Flow!

William Winters

William Winters
President

Save Money with Flow Dough!


Select Flow products contain valuable **Flow Dough Dollars!** Redeem them for free film or any one of hundreds of imaging accessories in this catalog.

Sign up and redeem your Flow Dough today at:
www.FlowDough.com



PSP Users

Barrier Envelopes	2
Positioning Accessories	4
Positioners	6

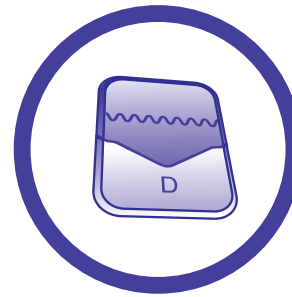

Film Users

Positioners	7
D Speed Films	8
F Speed Films	10
Barrier Packet Films	12
Panoramic Films	13
Single Use Positioning	14
X-Ray Mounts	16



Sensor Users

Sensor Positioning	19
Positioning Accessories	20
Sensor Protection	21



Protection

Lead-Free Aprons	22
Lead Aprons	23
Dosimeter	
Badge Service	24



Darkroom

Illuminators	25
Cassettes, Screens	25
Darkroom Accessories	26



Phosphor Plate Users

Barrier Envelopes for PSP Plates

Safe 'n' Sure barrier envelopes are the best phosphor plate envelopes you can buy. Safe 'n' Sure makes loading and unloading quick and easy. Reliable seals prevent cross-contamination. Rounded corners and soft edges provide added patient comfort. There are three Safe 'n' Sure versions to choose from.



Safe 'n' Sure Deluxe™

Soft clear vinyl face allows clinician to visually confirm the direction of the plate. Yellow EZ Glide Tab facilitates loading and sealing the envelope. Easy-Tear center seam allows for quick and easy plate removal.

Safe 'n' Sure Deluxe

Size 0: 80101 Size 1: 80111 Size 2: 80121 Size 3: 80131

300/box



LOAD



SEAL



UNLOAD



NEW THUMB NOTCH!

Safe 'n' Sure Econo™

New Thumb Notch on Easy-Tear Seam allows for fast removal of plate after exposure. Black on black vinyl for complete protection against light exposure.

Safe 'n' Sure Econo

Size 0: 80100 Size 1: 80110 Size 2: 80120 Size 3: 80130

300/box



LOAD



SEAL



UNLOAD



Barrier Envelopes for PSP Plates

Safe 'n' Sure OPT™

Pre-loaded cardboard inserts for plates with magnets.

Soft clear vinyl face allows clinician to visually confirm the direction of the plate. EZ Glide Tab facilitates loading and sealing the envelope. Easy-Tear center seam allows for quick and easy plate removal.

Safe 'n' Sure OPT

Size 0: 80141 Size 1: 80151 Size 2: 80161 Size 3: 80171

250/box



LOAD



SEAL



UNLOAD

Clean Ups™

PSP Cleaning Wipes. Individually wrapped, single-use wipes remove fingerprints, dirt and artifacts from your phosphor plates prior to exposure.

Simply tear open and wipe your plates!

Clean Ups 70107

50/box



Sticky Bites™

Sticky Bites are our most comfortable bite block! Made of tasteless polyethylene foam, Sticky Bites easily affix to film or plate to help you capture the perfect image. Disposable and 100% Latex free.

Sticky Bites 15780

500/box






Diagnostic & Positioning Accessories



NEW PRODUCT!

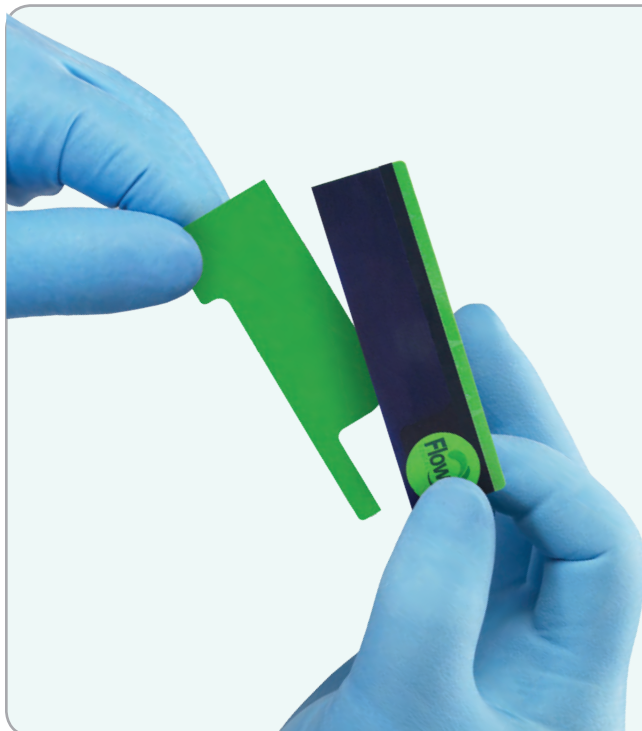

Exposure™ Articulating Film

Forget forceps...Exposure Articulating Film makes checking contacts a “snap!”

Exposure's built-in handle eliminates the need for forceps. Simply grasp, snap and bite! Easily exposes high occlusal contacts on both arches. Use for one tooth or an entire quadrant, on wet or dry teeth, even metal and non-metallic restorations!

Super thin EXPOSURE Articulating Film makes it easy to identify and adjust high contacts. Each film measures 4 ½ x 2 ¾”.

Exposure Articulating Film 

24300

200/box




Cushies™

Cushies are foam tubes that create a soft edge for added patient comfort. Pre-slit to easily slip onto the long or short sides of the barrier envelope or film edge. Great for lower anteriors!

Cushies 

15800

50/box

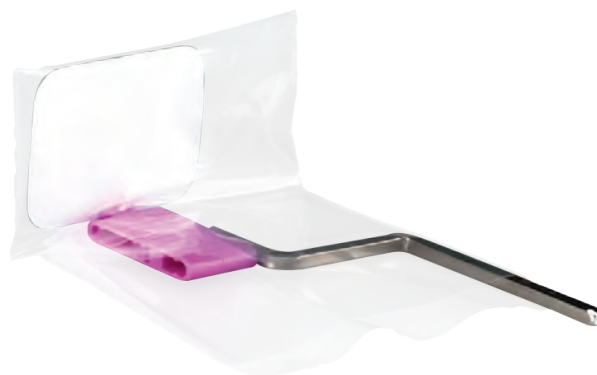

Bite Block Baggies™

Bite Block Baggies™ slip over the bite block and alignment arm to prevent cross contamination! After use, simply re-bag for the next patient. Bite Block Baggies are economically priced and come in a handy self-dispenser box.

Bite Block Baggies 

70104

500/box



Single-Use Positioning Accessories

Everything you need to protect your phosphor plates or film and ensure perfect alignment.



Stickems™

Adhesive bitewing tabs for use with either PSP barrier envelopes or film. For anterior, posterior and bite wings!

Stickems™

15720

500/box



SUPAs™

SUPA foam bite blocks hold your PSP envelopes at a 90° angle. For perfect alignment, use SUPA with our RAPiD Positioning Rings and Arms.

SUPAs

16087

100/box



Loop EZ™

Made of heavy duty paper board, Loop-EZ are made to accommodate adult or child size films or PSP envelopes. Great for taking bite wings!

Loop-EZ

Adult Econo Pack 15725

1,500/box

Adult: 15700

Child: 15710

500/box



Sticky Bites™

Sticky Bites are our most comfortable bite block! Made of tasteless polyethylene foam, Sticky Bites easily affix to film or plate to help you capture the perfect image. Disposable and 100% Latex free.

Sticky Bites

15780

500/box



Phosphor Plate Users

Positioning Systems

Everything you need to protect your film or phosphor plates and ensure perfect alignment.



LOAD



SNAP CASSETTE IN
BITE-BLOCK & EXPOSE



PROCESS

SMART™ Positioner

Flow's patented SMART system makes using phosphor plates just as easy, familiar and convenient as film.

Simply load your phosphor plate into the SMART Cassette and snap to close. The disposable cassette keeps phosphor plates safe from damage and protects them from moisture and light exposure. Place the cassette into the Universal Bite Block, attach an arm and an aiming ring, then position in the mouth.

SMART ensures perfect plate/cone alignment for every exposure. When you're done, take the plate out, process it and discard the cassette. The SMART system, including the bite blocks, arms and rings, is completely autoclavable.

The SMART Kit 40854

Includes:


- 1 Anterior Aiming Arm
- 2 Anterior Aiming Rings
- 2 Posterior Aiming Arms
- 2 Posterior Aiming Rings
- 8 Universal Bite Blocks
- 20 Size 1 Cassettes
- 50 Size 2 Cassettes

Reorder Number

- 408570 (1 ea)
- 408590 (1 ea)
- 408580 (1 ea)
- 408600 (1 ea)
- 408700 (1 ea)
- 408500 (50 ea)
- 408520 (100 ea)

Additional parts can be ordered separately.



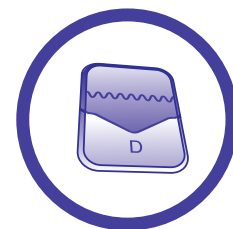

Positioning System for both Film and PSP Plates

RAPiD™ Positioner

For accurate x-rays, nothing beats the RAPiD system! Our RAPiD is the original color-coded film positioner. RAPiD bite blocks feature rounded corners and are flexible for easy loading and unloading. And our tabs ensure your film or PSP plates won't fall out during x-rays. Color-coding makes assembling the correct positioner quick and easy and RAPiD is autoclavable for over 100 uses.*

RAPiD is a right angle positioning device that assures correct film/cone alignment and prevents "cone cutting". It can be used for anterior, posterior, and bite-wing x-ray exams using any brand and size (3,2,1, and 0) of dental film.

RAPiD Kit with Bite Wing 408560



RAPiD Kits

All Kits available with short positioning arms.

Add a -S to the end of the item numbers below if you want to order a kit with short arms. Example: 408560-S

RAPiD Paralleling Kit with Bite-Wings \$ 408560

- 4 Aiming Rings: 2 Anterior/Bite Wing, 2 Posterior
- 4 Indicator Arms: 1 Anterior, 2 Posterior, 1 Bite Wing
- 16 Bite Blocks: 4 Anterior, 8 Posterior, 4 Bite Wings

RAPiD Paralleling Kit without Bite Wings \$ 409280

- 3 Aiming Rings: 1 Anterior, 2 Posterior
- 3 Indicator Arms: 1 Anterior, 2 Posterior
- 18 Bite Blocks: 6 Anterior, 12 Posterior

RAPiD Bite Wing Kit with #2 Vertical Bite Blocks 409250

- 1 Aiming Ring, 1 Indicator Arm, 4 Bite Blocks

RAPiD Bite Wing Kit with #2 Horizontal Bite Blocks 409260

- 1 Aiming Ring, 1 Indicator Arm, 4 Bite Blocks

RAPiD Bite Wing Kit with #3 Horizontal Bite Blocks 409300

- 1 Aiming Ring, 1 Indicator Arm, 4 Bite Blocks

RAPiD Bite Wing Kit with #0/1 Horizontal Bite Blocks 409320

- 1 Aiming Ring, 1 Indicator Arm, 4 Bite Blocks

RAPiD Anterior Kit with Size 0/1 Anterior Bite Blocks 408540

RAPiD Posterior Kit with Size 2 Horizontal Bite Blocks 408550



Anterior Part Numbers

1 Ant. & Bite Wing Ring	408590
25 Anterior Bite Blocks	408610
1 Anterior Arm	408570
1 Shorty Arm	408570-S



Posterior Part Numbers

1 Posterior Aiming Ring	408600
25 Posterior Bite Blocks	408620
1 Posterior Arm	408580
1 Shorty Arm	408580-S



Anterior & Posterior

Bite Blocks 10 of each	408630
------------------------	--------

Bite Wing Part Numbers

6 #2 Hor.	409290
6 #2 Vert.	420410
6 #0/1 Hor.	409330
6 #3 Hor.	409310
1 Bite Wing Arm	409270
1 Shorty Arm	409270-S




Endodontics

Endodontic X-Ray Positioner Kit	40303
1 Aiming Ring	40310
1 Indicator Arm	40311
3 Bite Blocks	40312



*Note: You cannot chemiclave or dry heat sterilize or use autoclave drying cycle.

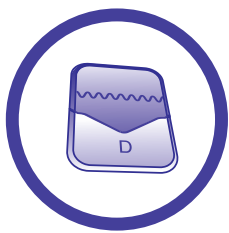

Film Users

Silver D: D Speed X-Ray Film

Flow Silver D is a high contrast film especially designed for today's most discerning dentist. This higher contrast delivers greater resolution and detail.

Now That's an Image!

In independent test results, Silver D delivered as much as 21% greater contrast and over 40% better resolution than other competitor's films.

What's more, our unique formulation stands up to oxidized chemistry better than other films. That means you'll get more films out of each gallon of developer used.

That's what on the inside. Outside, our vinyl packets are thinner, softer, and more flexible. Your patients will love our softer edges, your assistants will enjoy the added flexibility and easy-opening tabs!

Silver D has the same exposure settings as other D speed films and can be processed manually or automatically.



Econo Pak
600/Box

Value Pak
150/Box

Fresh Pak
100/Box



Size 2

1 1/4" x 1 5/8"
31 x 41 mm

DV-57 DV-58
FV-57 FV-58
DX-57 DX-58

Size 0

7/8" x 1 3/8"
22 x 35 mm

DV-53 DV-54
FV-53 FV-54
DX-53 DX-54

Speeds

Flow Dental's dental films are available in D speed (similar to Kodak Ultra-Speed*) or F speed (similar to Kodak Insight*) which requires roughly half the exposure time of D speed.


Specifications

Each Flow Intraoral film packet contains films sandwiched between black paper and lead foil. The film is embossed with a dot and the packets are color coded to indicate speed type and to ensure the correct side of the film is placed to the tube. Flow film is available in five intraoral and two extraoral sizes, and can be processed manually or automatically.

Flow Dental Silver D Film (D Speed) Exposure Time In Seconds

REGION		Exposure Factors	60kV 7mA	65kV 7.5mA	65kV 8mA	70kV 7mA	70kV 8mA	80kV 10mA
REGION	Maxillary	Anterior	0.55	0.32	0.30	0.27	0.24	0.10
		Premolar	0.73	0.43	0.40	0.37	0.32	0.13
		Molar	0.82	0.48	0.45	0.41	0.36	0.14
		Bite Wing	0.55	0.32	0.30	0.27	0.24	0.10
	Mandibular	Anterior	0.46	0.27	0.25	0.23	0.20	0.08
		Premolar	0.50	0.29	0.28	0.25	0.22	0.09
		Molar	0.55	0.32	0.30	0.27	0.24	0.10
		Bite Wing	0.55	0.32	0.30	0.27	0.24	0.10

Exposure times shown are for an 8 inch Source-image distance. For a 12 inch cone, multiply the exposure time by 2.25. For a 16 inch cone, multiply the exposure time by 4. To convert times to impulses, multiply the exposure times by 60. NOTE: For large patients, increase time by approximately 25%. For children, reduce exposure time by approximately 30%. For edentulous patients, reduce exposure time by approximately 25%. For DC x-ray generators, reduce exposure time by about one-third.

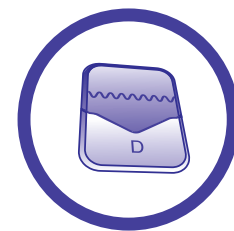

Silver D: D Speed X-Ray Film



DV-58 Film

D Speed Size 2 Adult Film

DV-58 D-speed, size 2 intraoral film is packaged in soft, supple vinyl packets for greater patient comfort. Perfect for a broad range of x-ray exams, including anterior, posterior, and bite-wing exams. DV-58 is available in three package sizes.

Econo Pak® \$
18204

600 films per box

Value Pak® \$
18200

150 films per box

Fresh Pak®
18600

100 films per box



DV-57 Double Film

D Speed Size 2 Adult Film - 2 Film Packets

All the great qualities of Silver D Film in a double film packet. One exposure = Two films! Perfect for referral situations or when a film duplicate is required for insurance processing. Twin D has the same exposure settings as regular Silver D, and can be processed manually or automatically.

Value Pak® \$
18201

130 films per box

Fresh Pak®
18601

100 films per box



DV-54 & DV-53 Double Film

D Speed Size 0 Child Film

All the great image qualities of Flow's Silver D in a size 0, child film size. Pedo Pak is perfect for kids ages 3 through 8. Pedo Pak is available in single film packets and double film packets, (Double D Pedo Pak). Flow D speed child film has 30% less exposure than adult Silver D. Pedo Paks can be processed automatically or manually.

Pedo Econo Pak® \$
18206


400 single films per box

Pedo Pak® \$
18235

100 single films per box

Double D Pedo Pak® \$
18236

100 double films per box

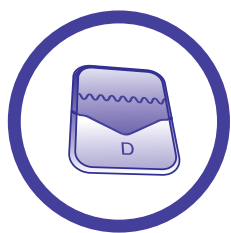

Xpress: F Speed X-Ray Film

As fast as any intraoral film on the market, Flow F speed films combine reduced radiation exposure with superb contrast and image resolution.

Don't Sacrifice Image Quality for Value, Speed & Comfort!

Available in sizes 2 and 0, in single and double film packets, Flow F speed takes conventional dental radiography to a new level in terms of dose minimization.

Inside: A fast film that doesn't sacrifice image quality for reduced radiation exposure! Outside: softer, thinner, more flexible vinyl packets your patients will love! Flow F speed films require approximately half the exposure time of D speed (see the exposure table below) and can be processed manually or automatically. Flow F speed has the same exposure settings as other F speed films making it easy to switch.



Econo Pak
600 per box

Value Pak
150 per box

Fresh Pak
100 per box

Size 2

1 1/4" x 1 5/8"
31 x 41 mm

DV-57 DV-58
FV-57 FV-58
DX-57 DX-58

Size 0

7/8" x 1 3/8"
22 x 35 mm

DV-53 DV-54
FV-53 FV-54
DX-53 DX-54

Speeds

Flow Dental's dental films are available in D speed (similar to Kodak Ultra-Speed*) or F speed (similar to Kodak Insight*) which requires roughly half the exposure time of D speed.

Specifications


Each Flow Intraoral film packet contains films sandwiched between black paper and lead foil. The film is embossed with a dot and the packets are color coded to indicate speed type and to ensure the correct side of the film is placed to the tube. Flow film is available in five intraoral and two extraoral sizes, and can be processed manually or automatically.



Flow Dental Xpress Film (F Speed) Exposure Time In Seconds

REGION		Exposure Factors	60kV 7mA	65kV 7.5mA	65kV 8mA	70kV 7mA	70kV 8mA	80kV 10mA
	Maxillary	Anterior	0.25	0.14	0.14	0.12	0.11	0.04
	Mandibular	Premolar	0.33	0.19	0.18	0.16	0.14	0.06
		Molar	0.37	0.22	0.20	0.19	0.16	0.06
		Bite Wing	0.25	0.14	0.14	0.12	0.11	0.04
		Anterior	0.21	0.12	0.11	0.10	0.09	0.04
	Mandibular	Premolar	0.23	0.13	0.12	0.11	0.10	0.04
		Molar	0.25	0.14	0.14	0.12	0.11	0.04
		Bite Wing	0.25	0.14	0.14	0.12	0.11	0.04

Exposure times shown are for an 8 inch Source-image distance. For a 12 inch cone, multiply the exposure time by 2.25. For a 16 inch cone, multiply the exposure time by 4. To convert times to impulses, multiply the exposure times by 60. NOTE: For large patients, increase time by approximately 25%. For children, reduce exposure time by approximately 30%. For edentulous patients, reduce exposure time by approximately 25%. For DC x-ray generators, reduce exposure time by about one-third.

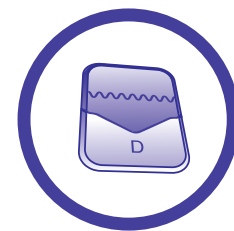

Xpress: F Speed X-Ray Film



FV-58 Film

F Speed Size 2 Adult Film

Flow's FV-58 is an F speed, size 2, intraoral dental film in soft vinyl packets. FV-58 requires approximately half the exposure time of DV-58.

FV-58 can be used for anterior, posterior and bite-wing x-ray exams in adults, and posterior and bite-wing exams in children 9 years or older and can be processed manually or automatically.

Econo Pak® \$

18215

600 films per box

Value Pak® \$

18310

150 films per box

Fresh Pak®

18620

100 films per box



FV-57 Double Film

F Speed Size 2 Adult Film - 2 Film Packets

All the great image qualities of our super fast F speed Adult dental x-ray film in a double film packet. One exposure = two films!

Perfect for referral situations or when a duplicate film is required for insurance processing. Flow F Speed film uses half the exposure settings of our Silver D, and can be processed manually or automatically.

Value Pak® \$

18410

130 films per box

Fresh Pak®

18621

100 films per box



FV-54 & FV-53 Double Film

F Speed Size 0 Child Film

All the great image qualities of Flow's Xpress F in a size 0, child film size. Pedo Pak is perfect for kids ages 3 through 8. Pedo Pak is available in single film packets and double film packets, (Double F Pedo Pak). Flow F speed child film has 30% less exposure than adult Xpress F Pedo Paks can be processed automatically or manually.

Pedo Econo Pak® \$

18207

400 single films per box

Pedo Pak® \$


18635

100 single films per box

Double F Pedo Pak® \$

18636

100 double films per box

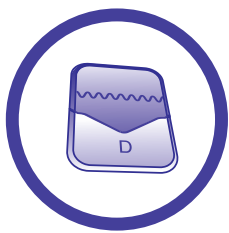

Barrier Packet Dental X-Ray Films

Butterfly Barrier Packet Films

Butterfly Barrier Packet Films are your preferred alternative to Kodak's Clinasept® film packets! Our easy-opening Butterfly packets come pre-loaded with your favorite Flow film, and protect you and the film from oral fluids. Enjoy the added safety and convenience while reducing your cost!

Butterfly's unique "wings" at the top of the envelope allow you to open the packet with minimum contact, keeping fingers free of blood and saliva. As an extra assurance of minimum contact, our envelopes open by pulling the "wings" apart from left-to-right as opposed to front-to-back. Once cleanly removed from its envelope, the remaining piece of film is a standard, easy-to-open, Flow film packet.

Butterfly films are available in the most popular sizes, styles and speeds: D and F Speeds, in Adult and Child sizes, and in single and double film packets.



	Speed	Films per Packet	Butterfly Name	Pieces per Box	Stock Number	Type of Film	Flow Dough
Size 2	D-Speed	Single	DSS-58	100	19600	DV-58	\$
		Double	DSS-57	100	19601	DV57	\$
	F-Speed	Single	FSS-58	100	19620	FV-58	\$
		Double	FSS-57	100	19610	FV-57	\$
Size 0	D-Speed	Single	DSS-54	75	19235	DV-54	\$
		Double	DSS-53	75	19236	DV-53	
	F-Speed	Single	FSS-54	75	19635	FV-54	\$
		Double	FSS-53	75	19636	FV-53	

*Clinasept is a registered trademark of Carestream Dental/Eastman Kodak Company

Paper Packet Intraoral Film

DX-58 Paper Pak® is a D speed, size 2 intraoral dental film in paper packets. DX-58 has the same exposure settings and developing techniques (manual or automatic) as all other D speed films.

18001

150 films per box

DX-42 Bite Chex® is a size 3, D speed film used to x-ray the coronal portion of contacting teeth.

18060


100 films per box with tabs for film positioning

DX-54 Paper Pak® is a D speed, size 0 child film in paper packets. DX-54 is for anterior, posterior and bite-wing exams in children from 3-8. DX-54 has the same exposure settings as any other D speed film and can be manually or automatically processed.

18035

100 films per box



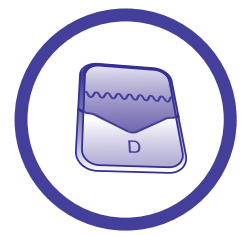

Panoramic Film & Accessories



Green Team® 512 & 612

Full speed, high contrast, orthochromatic (green sensitive) dental x-ray film available in 5 x 12" and 6 x 12" sheets.

Flow's Green Team® is ideally suited for all panoramic and TMJ imaging. The spectral sensitivity of Green Team is perfectly matched for use with all brands of green emitting, rare earth, intensifying screens. Green Team has wide surface area silver halide crystals to yield enhanced image sharpness.

Green Team® 512 & 612 \$

Green Team 512 - 5 x 12" Green Sensitive Pano Film **18095-G**
 Green Team 612 - 6 x 12" Green Sensitive Pano Film **18096-G**
50 films per box **See page 25 for Intensifying Screens.**



True Blue® 512 & 612

Full speed, blue sensitive dental x-ray film available in 5 x 12" and 6 x 12" sheets.

Flow's True Blue® Panoramic film yields excellent detail with high contrast. The blue spectral sensitivity of True Blue is perfectly matched for use with all brands of calcium tungstate (blue emitting) intensifying screens.

True Blue® 512 & 612 \$

True Blue 512 - 5 x 12" Blue Sensitive Pano Film **18095**
 True Blue 612 - 6 x 12" Blue Sensitive Pano Film **18096**
50 films per box **See page 25 for Intensifying Screens.**

Duplicating Films

Kopy Kat 512® is a 5 x 12" duplicating film. One sheet will reproduce one pano film or up to 27 intraoral x-rays. Packaged 50 sheets per box. Kopy Kat 512 can be exposed on any duplicator. Note that Kopy Kat film is notched to identify which side should face the film being exposed.

See our Kopy Kat X-Ray Film Duplicators on Page 26

Kopy Kat® 512 \$

5 x 12" *50 sheets per box* **18070**
 6 x 12" *100 sheets per box* **18071**

Flexible Panoramic Film Cassette

Flow Dental's Flexible Panoramic Film Cassettes are the perfect accessory for use with our Panoramic film. Fabricated of sturdy, heavy gauge vinyl for decades of reliable use, this flexible cassette curves to any desired position.

The cassette is easy to load and closes by two snaps on the side. One side is labeled "front" so you are always assured of placing the correct side to the x-ray source.

The cassette can be used with either 5 x 12" or 6 x 12" film and of course any brand film. As with any x-ray film the cassette must be loaded using standard film precautions in a darkroom setting. **See page 25 for Intensifying Screens.**

Flexible Panoramic Film Cassette


For Morita and Yoshida machines (15 x 30cm)
 For Panoramic Corp Machines (15 x 30 cm) comes with Tongue flap
 For Gendex machines (5 x 12 inches)

1 per box

11530
11532
11531



Internal screens needed. Sold separately. Call for details.

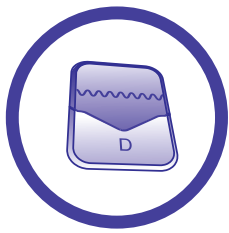

Film Users

Single Use Positioning Accessories

Everything you need to protect your film or phosphor plates and ensure perfect alignment.



Stickems™

Adhesive bite wing tabs for use with either PSP barrier envelopes or film. For anterior, posterior and bite wings!

Stickems™ 15720

500/box



SUPAs™

SUPA foam bite blocks hold your film or PSP envelopes at a 90° angle. For perfect alignment, use SUPA with our RAPiD Positioning Rings and Arms.

SUPAs 16087

100/box



Loop EZ™

Made of heavy duty paper board, Loop-EZ are made to accommodate adult or child size films or PSP envelopes. Great for taking bite wings!

Loop-EZ

Adult Econo Pack 15725

1,500/box

Adult: 15700

Child: 15710

500/box




Sticky Bites™

Sticky Bites are our most comfortable bite block! Made of tasteless polyethylene foam, Sticky Bites easily affix to film or PSP plate to help you capture the perfect image. Disposable and 100% Latex free.

Sticky Bites 15780

500/box

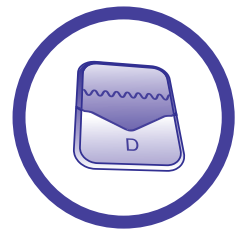

Diagnostic & Positioning Accessories

NEW PRODUCT!



Exposure™ Articulating Film

Forget forceps...Exposure Articulating Film makes checking contacts a “snap!”

Exposure's built-in handle eliminates the need for forceps. Simply grasp, snap and bite! Easily exposes high occlusal contacts on both arches. Use for one tooth or an entire quadrant, on wet or dry teeth, even metal and non-metallic restorations!

Super thin Exposure Articulating Film makes it easy to identify and adjust high contacts. Each film measures 4 ½ x 2 ¾”.


Exposure Articulating Film 💰 24300
200/box



Cushies™

Cushies are foam tubes that create a soft edge for added patient comfort. Pre-slit to easily slip onto the long or short sides of the barrier envelope or film edge. Great for lower anteriors!


Cushies 💰 15800
50/box

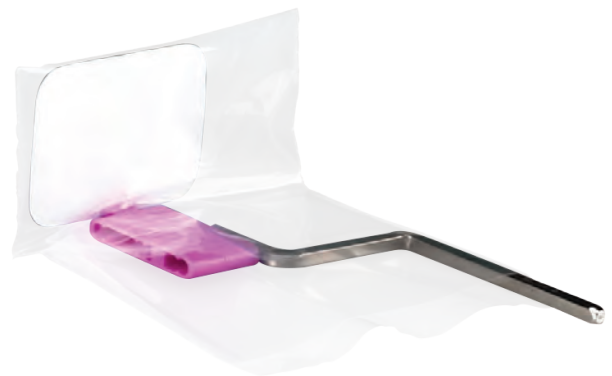

Bite Block Baggies™

Bite Block Baggies™ slip over the bite block and alignment arm to prevent cross contamination! After use, simply re-bag for the next patient. Bite Block Baggies are economically priced and come in a handy self-dispenser box.

Bite Block Baggies 💰 70104
500/box

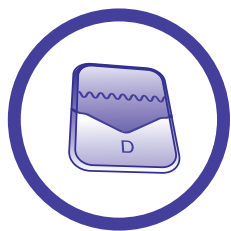

X-Ray Mounts

All Flow Dental X-Ray mounts are manufactured of the finest archival materials. Flow makes four types of mounts, each available in 19 styles. Whichever your preference, Flow mounts will give you decades of use.

Perfect Pocket™ Mounts

Masked pocket mounts were never this good! A great alternative to Rinn's E-Z View* masked mounts, Perfect Pocket mounts feature an innovative design loaded with benefits others don't have. First our material is thinner than others, making it easier to load films. Next our vinyl is clearer than others, so you have a better view of your films. Every Perfect Pocket mount features rounded corners—no more finger cuts while loading x-rays! We also added additional space for writing patient information, too!

[Perfect Pocket Mounts](#) 💰

See chart on page 18 for item #s 100/box



Klear-Vue™ Cardboard Mounts

Three layers of cardboard form a recessed pocket and lip that hold film in place. Flow mounts are printed with room for patient data (name, date, notations). Mounts can be imprinted with personalized doctor information. Klear-Vue mounts are comparable to Rinn's EeZee mounts*.

[Klear-Vue Cardboard Mounts](#) 💰

See chart on page 18 for item #s 100/box




Greylight™ Mounts

Similar to the Rinn EZ-Tab* and ADA* style mounts, our Greylight Mounts feature open window viewing. Films are held in place via embossed frames with raised borders that load easily and lock films firmly into place. Greylight mounts can be imprinted with custom doctor information.

[Greylight Mounts](#) 💰

See chart on page 18 for item #s 100/box

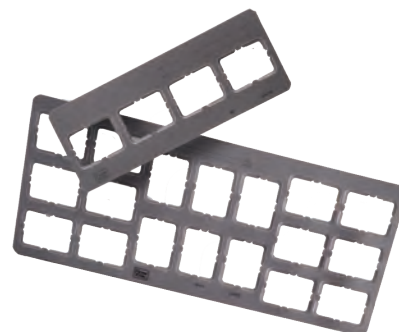

Grip Right™ Plastic Mounts


Similar to the Rinn EeZee Mount* plastic and Greene* TransView Mounts. Grip Right Plastic mounts are made from a thin semi-rigid vinyl with die-cut tabs. Thinner than other mount styles, Grip Right mounts allow for compact storage of your patient records. Can be imprinted with custom doctor information.

[Grip Right Plastic Mounts](#) 💰

See chart on page 18 for item #s 100/box



*ADA Mounts are a Registered Trademark of ADA Dental Products. *Rinn EeZee Mount, EZ-View, EZ Tab and Greene Transview are Registered Trademarks of Dentsply International. *Trollmount is a Registered Trademark of Trollplast.

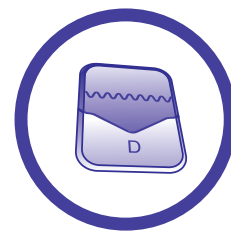

Tear-Away X-Ray Mounts

Mounts that allow you to custom-fit your x-ray records to your needs




Perfect Pocket™ Tear-Away Mounts

The first ever masked pocket tear-away mounts you can buy! Now, with Flow's new Tear-Away Perfect Pocket Mounts you can use as many windows as you need, no more... no less!

Our Perfect Pocket Tear Away Mounts are the same great mounts as our standard Perfect Pocket Mounts (see page 16) but in a format that allows you to tear off two mounts at a time. And each one-over-one window section is a complete mount, with space for patient info and dates to be written in. The windows accommodate a 2 H film. The dark blue of the mount background and the rounded corners of the windows make a perfect masking device, assuring you maximum contrast and illumination of your x-rays.


Perfect Pocket Tear Away Mounts come in a box of 80 mount sheets, 12 windows per sheet. That gives you a total of 960 windows per box. Each sheet is 11 1/2" long, 4 1/8" high making storage of whole sheets easy in a standard file draw set to accommodate letter size pages. Individual one-over-one sections are 1 7/8" long. Like our standard Perfect Pocket Mounts, our Tear Aways have rounded corners to protect you from sharp edges and their rigid construction allows for easy x-ray insertion.

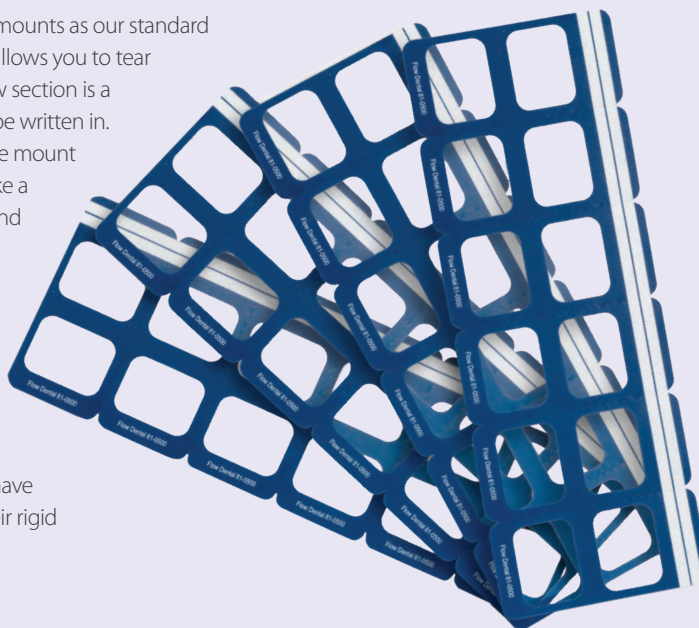

Perfect Pocket Tear-Away Mounts \$

81-0500

80 pieces per box (960 Windows)

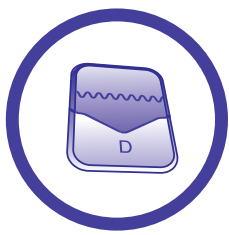

X-Ray Mounts Reference Chart


Layout	Description		Perfect Pocket Mounts	Graylight Plastic Mounts	Klear-View Cardboard Mounts	Grip Right Plastic Mounts	
	1H	1 Wind. Horiz Size 2, Size 0	81-6011	601-200G	25200	23200 23607	
	2H	2 wind. Horiz. Size 2, Size 0	81-6021 81-4021	602-203G 402-608G	25203 25608	23203 23608	
	2HBT	2 wind. Horiz. B Wing	81-7021	702-206G	25206	23206	
	4H 4H "0"	4 Wind. Horiz. Size 2 4 Wind. Horiz. Size 0	81-6040 81-4040	604-209G 404-609G	25209 25609	23209 23609	
	4V	4 Wind. Vert. Size 2		604B-510G			
	4H-2R	4 Wind. Horiz. 2 Rows, Size 2	81-6045	6045-239G	25239	23259	
	6H	6 Wind. Horiz. 2 Rows, Size 2	81-6061	606-215G	25215	23215	
	6-2V	6 Wind. 4 Horiz. 2 Vert., Size 2		606B-216G	25216	23216	
	7-3V	7 Wind. 4 Horiz. 3 Vert., Size 2		607-217G	25217	23217	
	8H	8 Wind. Horiz. 2 Rows, Size 2	81-6081	608-219G	25219	23219	
	10H	10 Wind. Horiz. 2 Rows, Size 2		610-220G			
	12H	12 Wind. Horiz. 2 Rows, Size 2	Perfect Pocket Tear-Away Mounts - 80 pcs/box, 960 windows/box 81-0500				
	10BT - "0" 10BT - "2" 10BT - "3"	10 wind B wing Size 0 Size 2 Size 3	81-2345 81-3456	610-261G 610-262G	25260 25261 25262	23260 23261 23262	
	10-6V	10 Wind. 4 Horiz. 6 Vert. Size 2			25240	23240	
	14-6V	14 Wind. 8 Horiz. 6 Vert. Size 2 Size 0		614-225G 614B-625G	25225 25625	23225 23625	
	14-6V	14 Wind. 8 Horiz. 6 Vert. Size 2, Size 0	81-8142				
	16-6V	16 Wind. 10 Horiz. 6 Vert., Size 2	81-6162	616-226G	25226	23226	
	16-6V 2BT	16 Wind. 8 Horiz. 6 Vert. Size 2, B Wing		716-227G	25227	23227	
	18-6V	18 Wind. 12 Horiz. 6 Vert. Size 2	81-6180	618-228G	25228	23228	
	20H	20 Wind. Horiz. Size 2	81-4567	618-230G			
	20-12H "2" BV	20 Wind. 12 Horiz. Size 2, 8 Vert. Size 0	81-6200	610-270G			



Perfect Sensor Positioning

Imagine... just one bite block for every patient!

Sensibles™ - Universal Sensor Holder The Sensible Way to X-Ray!

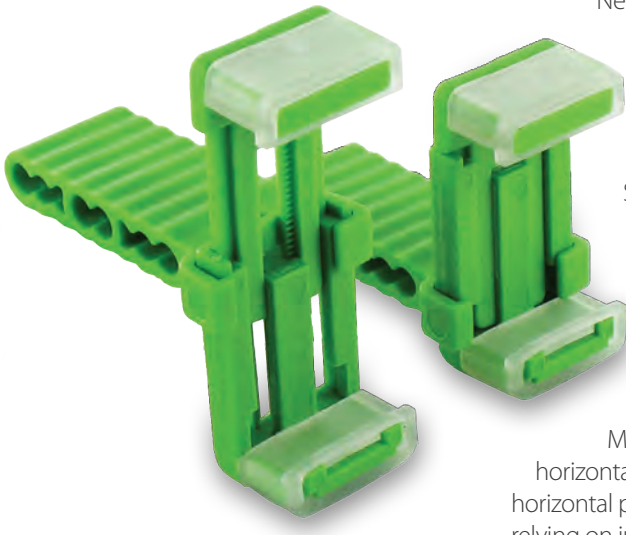

New, improved Sensibles is a "universal" sensor positioning system that does it all. Our patent-pending ratcheting bite block locks down tight to form a custom-fit around any size sensor in either the Horizontal or Vertical position. After use, simply press the quick release button and autoclave your Sensible Bite Block.


Sensibles now feature a unique "floating" bite block. When taking a bite wing you simply slide the bite block to the approximate mid-line of your sensor! It's quick, easy, and one bite block does it all!

Sensibles come in two sizes. Large is suited for most size 1 & 2 sensors in the vertical position. Medium is suited for most size 1 & 2 sensors in the horizontal position, and size 0 sensors in either the vertical or horizontal positions. Once you try them, you'll quickly find yourself relying on just one size Sensible for the majority of your x-rays!



*Patent Pending

Color coded sizes. Quick and easy to use. Versatile and autoclavable. Economically priced too. What a Sensible way to take x-rays!

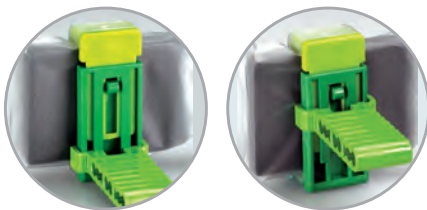

Using Sensibles is Easy As 1-2-3:



Sensibles Deluxe w/Bite Wing

6 Large & 10 Medium Bite Blocks
2 Anterior Positioning Arms
2 Anterior/Bite Wing Aiming Rings
2 Bite Wing Positioning Arms
2 Posterior Positioning Arms
2 Posterior Aiming Rings
50 Deluxe Comfees Samples
40915

Sensibles Deluxe Starter Kit

4 Large & 8 Medium Bite Blocks
2 Anterior Positioning Arms
2 Anterior Aiming Rings
2 Posterior Positioning Arms
2 Posterior Aiming Rings
50 Deluxe Comfees Samples
40900 \$

Sensibles Econo Starter Kit

2 Large & 2 Medium Bite Blocks
1 Anterior Positioning Arm
1 Anterior Aiming Ring
1 Posterior Positioning Arm
1 Posterior Aiming Ring
25 Deluxe Comfees Samples
40901 \$

Sensibles Bite Wing Kit

2 Large & 2 Medium Bite Wing Bite Blocks
2 Anterior/Bite Wing Aiming Rings
2 Bite Wing Alignment Arms
40912 \$

Sensibles Refill Kit 40910
12 Large (Two-Tone) Bite Blocks

Sensibles Refill Kit 40909
12 Medium Bite Blocks

1 Anterior Positioning Arm 408570

1 Anterior/Bite Wing Aiming Ring 408590

1 Posterior Positioning Arm 408580

1 Posterior Aiming Ring 408600

1 Bite Wing Positioning Arm 409150



Sensor Users

Perfect Sensor Positioning



Sticky Bites™

Sticky Bites™ make sensor positioning fast and easy. Simply peel and stick. The soft polyethylene foam is comfortable to bite down on, and can be used for Anterior, Posterior, or Bite-Wing images. Also protects sensor wires from patient biting! Quick, easy, affordable, and 100% Latex free, New Sticky Bites help you take a great image!


Sticky Bites 15780 500/box



Stickems™

Stickems™ self-adhesive bite wing tabs affix to the sensor sleeve. Packed in an automatic dispenser box, simply pull, fold, and stick!


Stickems 15720 500/box

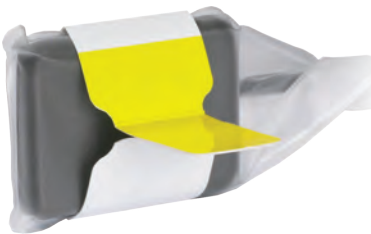

Digi Loops™

Digi Loops™ are made of heavy duty paperboard that simply 'loop' around the sensor. They're color-coded so you can easily pick out the proper size. Large and Medium hold most size 2 sensors. Small is suitable for most size 1 and size 0 sensors. Packaged 300 per box. Also available in an Assortment with 100 of each size.


Digi Loops Large: 15732 Medium: 15731 Small: 15730 Assortment 15733
300/box

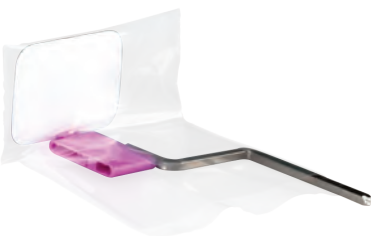

Bite Block Baggies™

Bite Block Baggies™ slip over the bite block and alignment arm to prevent cross contamination! After use, simply re-bag for the next patient. Bite Block Baggies are economically priced and come in a handy self-dispenser box.

Bite Block Baggies 70104 500/box

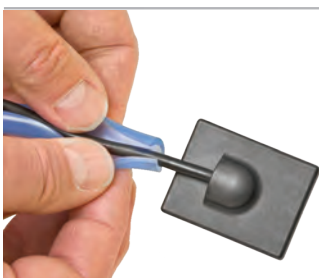

Cable Saver™

Cable Saver™ protects your sensor cables from being stepped on, rolled over, twisted, and tangled. Easy to install and easy to keep clean, Cable Savers will add years to the life of your sensor cables. Available in 8' lengths. Simply cut to size!

Cable Saver 80500 8' length



Perfect Sensor Protection

NEW PRODUCT!

Snuggies™ Custom-Fit Sensor Sleeves

Introducing new Snuggies, the Snug-Fit Sensor Sleeve. Now you can protect your equipment and your patients with a snug, custom-fit sensor sleeve—regardless of sensor size. Simply peel, stick, and go!

Snuggies works with sensors of any size, shape, and thickness. Patients will love Snuggies because there's less loose plastic to gag on. Practitioners will love Snuggies because one-size-fits-all, eliminating the need to buy and stock sleeves of many sizes. Plus, with Snuggies, you don't have to twist the sleeve around delicate sensor cords to keep the sensor from sliding around in the sleeve.

To custom-fit your Snuggies, just peel off the clear tab, then pull the white tab snug against the sensor. That's it!

Snuggies are economically priced and packaged 500 per box. Order a box today and see for yourself!

Snuggies Custom-Fit Sensor Sleeves 80222
500/box



Comfees™ Comfortable Barrier Sleeves

Comfees™ are simply the best barrier sleeve you can buy, making sensor imaging easier and quicker for the technician, safer and more comfortable for the patient.

Comfees Deluxe™ feature a super-soft vinyl material that makes loading the sensors effortless. The added strength of the material also protects against tears.

Both the Deluxe and Econo versions feature rounded corners for better fit and added patient comfort. Both also feature our unique tear-away seams, so you can easily peel back the Comfee instead of tugging on the sensor wire.

Comfees Deluxe

Large: 80210

Medium: 80200

500/box


Comfees Econo

Large: 80211

Medium: 80201


500/box

Deluxe

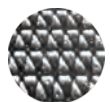

Textured Vinyl

Econo

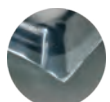

Smooth Seams

Large: Fits most size 2 sensors. Medium: Fits most size 0 and size 1 sensors.



Protection

Green, Lead-Free Protective Aprons

Flow Dental's new No-Lead Aprons make it easy to be Green!

Our lead-substitute material provides 0.25mm lead equivalency protection while making our aprons 25% lighter than traditional aprons. And No-Lead aprons don't require special disposal.

They are better for you, better for your patients, and better for the environment! The safety of lead protection without the weight... or the lead!

Our new aprons are made of environmentally-friendly, lead-free, non-toxic materials. They're lighter and more flexible than conventional lead aprons, but have the same protective qualities!

All Flow No-Lead Protective Aprons provide 0.25mm lead equivalency protection. Aprons are backed with a tear-resistant cling fabric. Aprons are available with and without attached thyroid collars, and thyroid collars can also be ordered separately.



Panoramic Apron in Camouflage

Bib Apron with Attached Collar in Royal Blue

Separate Thyroid Collar in Paw Prints



Uni-Hanger Item #16401




Bib Apron in Tye Dye

How to order your Flow No-Lead Apron

To order, select from any of the patterns shown on the right. Use the stock number below and add a fabric number after it. **Example:** If you want a Bib Apron, without a thyroid collar, in the Tye Dye fabric, the correct product number is 75050NL-30.

Description	Size	Stock Number
Adult Bib Apron without Collar	24" x 26"	75050NL (+ color #) \$
Adult Bib Apron with Attached Collar	24" x 26"	75051NL (+ color #) \$
Child Size Bib Apron with Attached Collar	20" x 20"	75053NL (+ color #) \$
Adult Panoramic Apron	24" x 27"	75061NL (+ color #) \$
Child Panoramic Apron	20" x 20"	75062NL (+ color #) \$
Thyroid Collar	5" x 27"	75056NL (+ color #) \$
Apron Uni-Hanger		16401 \$

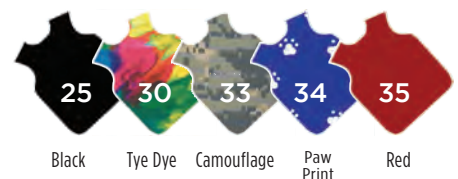



Traditional Lead Protective Aprons

Traditional lead protection for your patients in fashionable styles and colors.

Flow lead protective aprons are the industry standard. Manufactured from lead rubber, sandwiched between an attractive vinyl front and a tear-resistant cling back, our aprons clean easily with soap and water.

Made for years of heavy use, all our lead aprons provide 0.3mm lead equivalency protection, and all styles are also available in our No Lead option. When ordering, just write "NL" after the stock number and before the color number.

And for the perfect accessory to store your apron, Flow carries the UniHanger which allows you to hang your apron properly, on a hook or in a closet, to prevent folds which will damage the lead lining.

Aprons are available with and without attached thyroid collars, and thyroid collars can also be ordered separately. Select from any of the fabric colors shown below.



Bib Apron with Collar in Teal



Bib Apron in Purple



Child Bib Apron



Panoramic Apron in Burgundy



How to order your Flow Lead Apron

To order, select from any of the colors shown on the left. Use the stock number below and add a fabric number after it. **Example:** If you want a Bib Apron, without a thyroid collar, in the Teal fabric, the correct product number is 75050-14.

Description	Size	Stock Number
Adult Bib Apron without Collar	24" x 26"	75050 (+ color #) \$
Adult Bib Apron with Attached Collar	24" x 26"	75051 (+ color #) \$
Child Size Bib Apron with Attached Collar	20" x 20"	75053 (+ color #) \$
Adult Panoramic Apron	24" x 27"	75061 (+ color #) \$
Child Panoramic Apron	20" x 20"	75062 (+ color #) \$
Thyroid Collar	5" x 27"	75056 (+ color #) \$
Apron Uni-Hanger		16401 \$



Protection

Dosimeter Badge Service

Keep you and your staff safe

Working around x-rays is a concern for everyone; dentists and staff. X-ray exposure is cumulative. It's vital you maintain an accurate history and monitor your daily exposure. And for that, you want to use the most accurate service you can get. And you want it at an affordable price.

The answer: Flow's new Thermoluminescent Dosimeter TLD Service.

TLD is state-of-the-art technology. Four elements are incorporated into our Dosimeters, so we measure down through three layers of tissue. And our processing software eliminates any guesswork or errors.

Here's how to keep you and your staff safe:

Everyone receives his or her own individual badge, to be worn at all times while at work. These lightweight badges measure 1"x2" and clip easily onto clothing. Each badge is assigned a unique identification code.

Every 3 months a new replacement badge arrives in the mail, along with a return envelope. You send us back your old badge. We generate a detailed report and it's sent to you for your records. And if the readings ever show exposure in excess of published NRC safety and standards, your Flow service staff will contact you immediately by phone!

In addition to personal safety concerns, OSHA and many state agencies require that personnel use x-ray monitoring devices when operating, or in proximity to, any radiation-emitting equipment. Failure to comply may result in fines. And having individual documentation on the history of your staff's exposure to x-rays may prove vital in case of litigation any time in the future.

For only \$122 a year per person you get your own individual badge and quarterly processing and reporting. That's it! No extra fees! No set-up charges, no charge for replacement reports, or even for a changeover in staff personnel. And if you're currently enrolled in another service, you can sign up now and won't start billing you until your present service is finished!

Enrollment is quick and easy:

- 1) Select who at your office will be the main contact and have them fill out the form below. If you need to list more names, copy the form.
- 2) Fax the form to (631) 242-1001
- 3) You will receive a return fax confirmation within 72 hours. Your first delivery of badges will arrive within 10 business days.

Your subscription includes automatic badge replacement every 3 months, for one year, reply envelopes for returning badges for processing, and copies of our quarterly reports (including reprints at no charge).



X-RAY BADGE SERVICE ORDER FORM Refer to Stock Number 16102 Fax this form to 631-242-1001

Account Number: _____ Clinic: _____ Dr.: _____

Address: _____ City: _____ State: _____ Zip: _____

Contact: _____ Phone: _____ Email: _____

☐ Yes, I want to switch my service. My contract ends: ____/____/____

☐ Yes, I want Automatic Annual Renewal

1. _____

Employee ID Number: _____

2. _____

Employee ID Number: _____

3. _____

Employee ID Number: _____

4. _____

Employee ID Number: _____

5. _____

Employee ID Number: _____

6. _____

Employee ID Number: _____



Illuminators



Slimline® Illuminator

Made for desk or wall mounting, Flow's Slimline Illuminator series gives you the maximum viewing possibilities with an economy of space and price. The Slimline's 4 ¾" x 13 ¾" viewing area accommodates any mount up to a full series.

The Slimline comes with bulbs, 6 foot line cord and keyholes in back for easy wall mounting.

Specifications: 6 ½" High x 15 ¼" Wide x 2 ½" Deep
Approx. shipping weight: 3lbs.

Dental Illuminators	Viewing Area	Stock Number
Heavy Duty White Plastic	13 ½ x 5	21502
White Enamel Finish	13 ½ x 6	21730
Beige Enamel Finish	13 ½ x 6	21730-BF
Gray Enamel Finish	13 ½ x 6	21730-GF
Gray Enamel Finish	13 ½ x 10 ½	21740
Cephalometric Tracing Box	10 ½ x 10 ½	21750

Slimline Cassettes

Slimline cassettes are easy to open and close, thin for easier storage, and lightweight for easier handling. They feature a 'snap' type open and close latch.

5 x 12" **11365**
6 x 12" **11366**

Intensifying Screens


Helps reduce exposure to patients and eliminate scatter radiation.

BG Green: 400 Speed **Blue: 400 Speed**
5 x 12" **10736** 5 X 12" **10727**
6 X 12" **10737** 6 X 12" **10768**

Cassette and Screen Combinations are also available. Please call for details and ordering information.



Darkroom Accessories

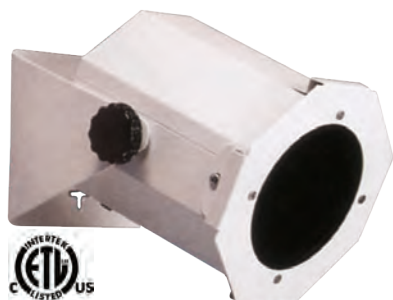

Swivel Safelight

Perfect for direct spotlighting, this wall mounted unit swivels to any position. Easy to install via 2 keyholes, the Swivel Safelight includes a GBX type filter, light bulb, and a 6 foot cord with a three prong plug.

Specifications: 7 1/2" H x 6" W.
Approx. ship weight 5lbs.

Swivel Safelight 19198



Angular Safelight

The Angular Safelight is ideal where a direct light is wanted. Easy to install via 2 keyholes, the Angular Safelight includes a GBX type filter, light bulb, and a 6 foot cord with a three prong plug.

Specifications: 10" H x 7 1/2" W x 5 1/2" D.
Approx. ship weight 5lbs.

Angular Safelight 19103



Combination Safelight

This unique two-compartment unit features a safelight at the bottom and a normal incandescent room light at the top. It is filtered by our GBX type filter and comes with a 6 foot cord with a three prong plug.

Specifications: 15 1/2" H x 10" W x 2 1/2" D.
Approx. ship weight 12lbs.

Combination Safelight 19106

The Swivel and Angular safelights are also available with an OPTICAL SENSOR that turns the safelight on automatically when the regular light is shut off (the unit can still be shut off with the on/off switch). The stock number for the Swivel safelight with the sensor feature is 19298. The stock number for the Angular is 19203.



Circular Safelight

Our economy model screws into any existing outlet. Made of aluminum, filtration is from a GBX type filter that is removable for bulb replacement. Specifications: 5 1/4" L x 4 1/8" D. Approximate shipping weight 5 lbs.

Circular Safelight 19104

Kopy Kat X-ray Film Duplicators

Flow X-ray offers two duplicators. The 10x12 inch surface of the 21522 is perfect for duplicating full series or panoramic x-rays, the 21520 will duplicate all films up to 6x12 inches.

Our duplicators can be used with any duplicating film. To operate the duplicators, place the film to be duplicated on the plexiglass screen with a sheet of duplicating film on top*, close the cover and set the exposure time. A red light signals you when exposure is complete.



Screen Size	Unit Dimensions			Shipping Weight	Stock Number
	L	W	H		
6 x 12"	18.25"	8"	6.25"	12lbs.	21520
10 x 12"	18.25"	12"	6.25"	22lbs.	21522

*Don't forget to order Kopy Kat X-Ray Film... see Page 13!



Darkroom Accessories

Timers

The Flow Big Timer - Our analog timer features a giant 8" diameter dial that gives you timing from 1 second to 60 minutes. It automatically shuts-off a duplicator or other electrical appliances, and has a clear buzzer signal. **18900**

The Electronic Timer - Our digital timer is a pocket sized 2 5/16"W x 2 1/16"H x 1/2"D. It attaches to anything: it has a spring-loaded clip that converts into a stand, a magnet for metal surfaces, and a hanging hole. Times from 100 minutes down to 1 second and emits a tone when timing period is completed. **18906**



ID Printers

ID Printers allow you to expose patient information onto extraoral film.

Auto ID Printer - has an adjustable exposure timer. Set the timer, insert ID Card, and press the exposure plate. **17105**

Manual ID Printer - Just insert the film and ID Card, hold down the exposure plate, and release. Compatible with Kodak and Wolf ID Cards. **17104**

ID Printer Cards - 3"x5" 1,000 per box:

Standard ID Cards - With blank spaces for patient name and address. **15100**

Personalized Standard ID Cards - Same as above with custom imprint, 3 line maximum. **15102**

Double Printed Personalized ID Cards - Identification blocks printed top and bottom. Custom Imprinted **15151**



Econo-Top Film Bin

Now you can safely store your Pano Film without having to worry about accidental exposure. Flow's new Econo-Top Film Bins are your economical choice!

The Econo-Top Film Bin is made of durable plastic so it's lightweight and easy to handle. It features an easy-to-use lid that secures with Velcor fasteners, and comes fully assembled.

25130 Holds up to 300 5" x 12" or 6" x 12" films





Darkroom Accessories


Thermometers

A) Digital Thermometer Registers temperatures from -40° to 300° F. Just place it in the solution and press the button. The temperature is digitally displayed on an LCD. Runs on battery (included) and comes with a lifetime warranty. **19107**

B) Dial Thermometer Calibrated from 0° to 140° F. Includes floating ring. **18903**

C) Glass Thermometer Floats in the developer tank. Calibrated from 0° to 100° F. **18904**

D) Stainless Steel Thermometer Calibrated for accuracy between 20° and 150° F. **18905**

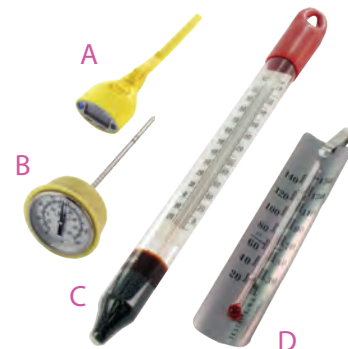

Hangers

E-Z Dip Film Developing Hangers Consist of a stainless steel strip, curved at the top, to which spring clips and identification tabs are attached.

Single Clip	16210	12 clip	16207	4 clip (long)	16203
Box of 12		14 clip	16208	10 clip (long)	16215
6 clip	16218	16 clip	16209	12 clip (long)	16216
8 clip	16219	18 clip	16205	14 clip (long)	16217
10 clip	16206	20 clip	16204		

Replacement Dental Hanger ID Tabs - **14000**
Box of 100

Small Single Clip - **14001**
Box of 12



Value Chem® Concentrated Developer and Fixer

Flow Developer has a Phenidone base for a higher concentration of developer, yielding better radiographs and slower oxidation of the solution. Flow Fixer is designed to fix a film in three minutes. A bottle of concentrate makes one gallon of working solution.

Developer & Fixer **24220**
1 each



Flow Clearing Sheets

Our 7 X 8 1/2" clearing sheets fit all intraoral and extraoral processors. Sheets help remove silver residue and algae off the roller racks. Simply feed a sheet through in the morning, after machine warm-up, and as many times during the day as necessary. Each sheet can be used up to 3 times.

Flow Clearing Sheets **70106**
50 sheets/box



Stock Number Index

STOCK #	PAGE #	STOCK #	PAGE #	STOCK #	PAGE #	STOCK #	PAGE #	STOCK #	PAGE #
10727	28	18200	9	23217	18	408540	7	75050	22
10728	28	18201	9	23219	18	408550	7	75051	22
10736	28	18204	9	23225	18	408560	7	75051NL	23
10737	28	18206	9	23226	18	408570	6, 7, 19	75053	22
10738	28	18207	11	23227	18	408570-S	7, 19	75053NL	23
11365	28	18215	11	23228	18	408580	7, 19	75056	22
11366	28	18235	9	23239	18	408580-S	7, 19	75056NL	23
11530	13	18236	9	23240	18	408590	7, 19	75061	22
11531	13	18310	11	23259	18	408600	7, 19	75061NL	23
11532	13	18410	11	23260	18	408610	7	75062	22
14000	28	18425	11	23261	18	408620	7	75062NL	23
14001	28	18600	9	23262	18	408630	7	80100	2
15100	27	18601	9	23608	18	408700	7	80101	2
15102	27	18620	11	23609	18	40900	19	80110	2
15151	27	18621	11	24220	28	40901	19	80111	5
15700	5, 14	18635	11	24300	4, 15	40909	19	80120	2
15710	5	18636	11	24301	15	40910	19	80121	2
15720	5, 14, 20	18900	27	24302	15	40912	19	80130	2
15725	5, 14	18903	28	24303	15	40915	19	80131	2
15730	20	18904	28	24304	15	409250	9	80141	3
15731	20	18905	28	24220	28	409260	9	80151	3
15732	20	18906	27	25130	27	409270	9	80161	3
15733	20	19103	26	25200	18	409270-S	9	80171	3
15780	3, 5, 14, 20	19104	26	25203	18	409280	9	80200	3
15800	4, 15	19106	26	25206	18	409290	9	80201	3
16087	5, 14	19107	28	25209	18	409300	9	80210	21
16102	24	19198	26	25215	18	409310	9	80211	21
16203	28	19203	26	25216	18	409330	9	80222	21
16204	28	19235	12	25217	18	420410	9	80500	20
16205	28	19236	12	25219	18	601-200G	18	81-0500	17
16206	28	19298	26	25225	18	602-203G	18	81-2345	18
16207	28	19600	12	25226	18	604-209G	18	81-3456	18
16208	28	19601	12	25227	18	6045-239G	18	81-4021	18
16209	28	19610	12	25228	18	604B-510G	18	81-4567	18
16210	28	19620	12	25239	18	606-215G	18	81-6011	18
16215	28	19635	12	25240	18	606B-216G	18	81-6021	18
16216	28	19636	12	25260	18	607-217G	18	81-6040	18
16217	28	21502	25	25261	18	608-219G	18	81-6045	18
16218	28	21520	26	25262	18	610-220G	18	81-6061	18
16219	29	21522	26	25608	18	610-261G	18	81-6081	18
16401	23	21730	25	25609	18	610-262G	18	81-6162	18
17104	27	21730BF	25	25625	18	610-270G	18	81-6180	18
17105	27	21730GF	25	402-608G	18	614-225G	18	81-6200	18
18001	12	21731	25	40303	7	614B-625G	18	81-7021	18
18035	12	21740	25	40310	7	616-226G	18	81-8142	18
18060	12	21750	25	40311	7	618-228G	18		
18070	13	23200	18	40312	7	618-230G	18		
18071	13	23203	18	408500	6	70104	4, 15, 20		
18095	15	23206	18	40850-10	7	70106	28		
18096	15	23209	18	408520	6	70107	3		
18095-G	13	23215	18	40852-10	7	702-206G	18		
18096-G	13	23216	18	40854	6	716-227G	18		



100 West Industry Court
Deer Park, New York 11729

Phone: 800.356.9729 | 631.242.9729

Fax: 631.242.1001

www.flowdental.com

info@flowdental.com