



Everything For Imaging And More



Illuminators
 Mammography Illuminators
 Podiatry Illuminators
 Chiropractic Illuminators
 Film Digitizers
 Monitors
 Adjustable PACS Workstations
 Protective Apparel
 Lead Free Protective Apparel
 Protective Eyewear
 Protective Barriers
 Medical X-ray Film
 Cassettes, Screens and Grids
 CR Phosphor Plates
 CR Cassettes, Screens and Grids
 Cassette Markers
 Digital And Mammo Markers
 Cassette Holders
 Positioning Devices
 Utility, CR and Cassette Carts
 Film Filing Cabinets
 Film Storage Bins
 MRI Carts
 X-ray Accessories
 Darkroom Accessories
 Ultrasound Gels And Disinfectants



Wolf XRay



The inSight X-ray Illuminator

State Of The Art
Illuminating
Technology
Applied To
Medical X-ray
Viewing

The new inSight illuminator from Wolf X-ray utilizes LCD technology to bring you an x-ray viewbox with a uniquely even light surface in a super-thin cabinet. Just 2 inches deep, it gives you the look of a recessed unit without the hassle of cutting holes in you wall.

Available as a single unit, two in one and four in one; each panel accommodates up to one 14 x 17" x-ray. All units come fully equipped with auto switches that Activate the viewing area when films are inserted, as well as a manual on/off switch.

The inSight comes with a one year warranty for defects in parts or manufacturing. Sleek modern, affordable film viewing; the new inSight is the perfect addition to you reading room.



	Stock Number	Viewing Area	Dimensions Height Length Depth	Approx. Shipping Weight	Power Consumption
Single Unit	27001	14"x17" 35.6x43.2	20 1/4"x21"x1 5/8" 515x535x40	6.5 lbs.	Max. 70W
2 Bank Unit	27002	28"x17" 71.2x43.2	20 1/4"x35"x1 5/8" 515x937x40	10.5 lbs.	Max. 140W
4 Bank Unit	27004	56"x17" 142.4x43.2	20 1/4"x54"x1 5/8" 515x1595x40	20.0 lbs.	Max. 280W

Illuminators - Definition of Terms and Summary Tables

Wolf X-Ray manufactures the most comprehensive range of illuminators available. Familiarizing yourself with "The Definition of Terms" and "Summary Tables" presented below will enable you to pick the Wolf Illuminator best suited to your needs.

Definition of Terms

Bulbs per 14" x 17" Viewing area: All illuminators come with fluorescent bulbs, the number of bulbs per 14" x 17" panel are listed below. Actual candela per square meter per panel is listed on the individual descriptive pages.

Cluster Switch: Switches include a "master switch" to shut off the entire unit. These are standard for 2 over 2 & up.

Reflectors: All Wolf illuminators are equipped with advanced reflectors which guarantee a diffused evenly lit viewing area.

Auto Switch: Switches on when film is inserted, off when it's removed.

Individual Switch: Individual "on/off" switch per 14" x 17" section.

Hospital Grade Plug: A spark proof hospital approved 3 prong plug.

Roller Grip: A roller bar grip extending full viewing length.

Magic Grip: Built in tension grip, no moving parts.

Wolf's E-"Z" Bar™ allows you to install Wall Hung units simply, quickly and without opening the unit!. Only critical dimension is height. Simply screw the provided bar in the desired spot, and mount the illuminator on it.

Recessed Mounting: Cabinet is installed in the wall. Units ordered for recessed mounting are provided with a kit that easily adapts standard units for recessed mounting. (See page 8 for details)

Mobile Stand: Mobile desk console unit for multiple banks.

Summary Table - Standard Illuminators

Illuminator Style	Bulbs per 14" x 17" Viewing Area				Switches			Film Retainers			Standard Hospital Grade Plug	Options					Catalog Page Number
	6	4	3	2	Cluster Switches	Auto Switch	Ind. Switches	Roller Grip	Magic Grip	"Z" Bar		Recessed Mounting	Hospital Grade Plug	Additional Cluster Switches	Mobile Stand	Handle	
Liberator		■			■	■	■	■		■	■	■		■	■		3
Liberator HiLo	■			Dims to 4 Bulbs	■	■	■	■		■	■	■		■	■		3
TrimlinePlus HiLo	■			Dims to 4 Bulbs	■	■	■	■		■	■	■		■	■		5
Trimline Plus		■			■	■	■	■		■	■	■		■	■		6
Trimline			■		■	■	■	■		■		■	■		■		7
Trimline Basic			■		■		■	■		■		■	■	■	■		7
Trimline III			■		■	■	■	■		■		■	■	■	■		6
MG-7			■				■		■				■			■	11
Econoline			■				■		■			■	■			■	11

The OR Series of Illuminators have been designed to meet the requirements of more restrictive areas such as operating rooms and as such are listed to a more stringent UL standard.

Summary Table - Specialty Illuminators – The OR Series

Illuminator Style	Bulbs per 14" x 17" Viewing Area				Switches			Film Retainers			Standard Hospital Grade Plug	Options				Catalog Page Number
	6	4	3	2	Cluster Switches	Auto Switch	Ind. Switches	Roller Grip	Magic Grip	"Z" Bar		Recessed Mounting	Additional Cluster Switches	Mobile Stand	Handle	
Liberator OR		■			■	■	■	■		■	■	■	■	■		4
Liberator HiLo OR	■			Dims to 4 Bulbs	■	■	■	■		■	■	■	■	■		4
TrimlinePlus HiLo OR	■			Dims to 4 Bulbs	■	■	■	■		■	■	■	■	■		15
Trimline Plus OR		■			■	■	■	■		■	■	■	■	■		15



Summary Table - Specialty Illuminators

Illuminator Style	Specialty	Film Capacity #	Light Intensity	Features	Catalog Page #
inSight	Super Thin LCD TV Technology Illuminator				1
Mammoline	Mammography	24 x 30 cm films 18 x 24cm films	Dimmer controlled from 270 to 7,600 candela per square meter	Four independently moving masking shutters and an E-Z Bar for easy wall mounting	14
Podolite	Podiatry & 10" x 12"	Accommodates 8"x10"	2 Bulb	Wall hung or desk mounted.	17
Trimline 36 & 51	Orthopedic	F m 1-2 14 x 36; film	3 or 4 Bulbs	Can be wall mounted or recessed	18
Vu Plus	Mammography	24x30 - 18x24cm	5 Bulb Hi/Lo	Wall hung or desk mounted.	14



The most unique
viewbox ever
built

The Liberator Illuminators

-  The Plus
-  The Hi-Lo


Never lose a film again! Wolf's Liberator series ingeniously attaches the film grip to the plexiglass to form a perfect film holding assembly that never goes out of alignment. The entire plexiglass panel & grip assembly come off quickly and effortlessly with a simple pull. Our Magna-Lock assembly holds the panels in place until you want to get inside. Changing a light bulb or cleaning inside the panel was never this simple. Whether it's installation, service, or operation the *Liberator* is at the top of the class.

The *Liberator* is available in two styles. The Liberator Plus is a 4 bulb illuminator, available from a single 14" X 17" viewing area up to 2 tiers of 84" X 17" each. Our Hi Lo version, available in the same sizes as the standard, features 6 bulbs per 14" x 17" viewing area. Each 14" x 17" area of the Hi Lo can be dimmed to 4 bulbs at the press of a button.

Any model can be wall mounted in minutes with our E-Z bar or beautifully recessed in a wall with our recessing kit. Specifications and ordering information for the Liberator are on the next page.

Wolf
Wolf
Wolf
Wolf
Wolf
Wolf
Wolf
Wolf
Wolf
Wolf
Wolf

Liberator And Liberator Hi Lo Illuminator Series Specifications

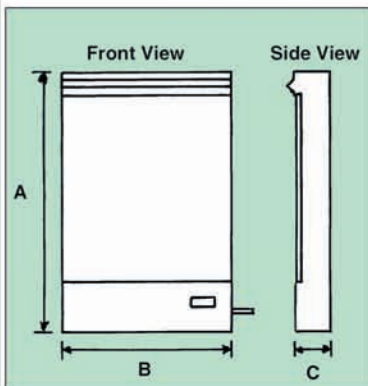

Letters reference dimensions shown in the right hand column.

Standard Features...

- Plexiglass viewing area and film retainer lift off for easy interior access.
- Film activated auto switches.
- Cluster switches and a master switch on 2 over 2 units and up.
- Gravity roller grip film retainer.
- Wolf's E-"Z" Bar™ for easy surface wall mounting. Simply mount the bar on the wall and slip the illuminator over it.
- Hospital grade plug.

Light Intensities

Standard Liberator

- 4 bulbs per 14" x 17" viewing area.

Hi Lo Liberator

- Variable Light intensities (6 bulbs to 4 bulbs) per each 14" x 17" viewing area.

Optional Features

- A frame for recessed mounting.
- A mobile stand.
- Custom room installations.

Ordering Information


- Use stock numbers to place your order
- If you want a recessed illuminator, order the appropriate Recessing Kit.

Description


Liberator Stock Numbers Standard HI Lo Recess Kit

Viewing Area

Wall Mounted Dimensions A B C

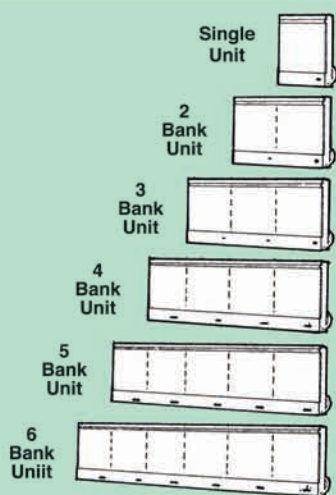

**NOTE: TO ORDER
ANY OF THESE
SIZES AS AN OR
(Described on the
previous page)
add an OR after
the stock number**

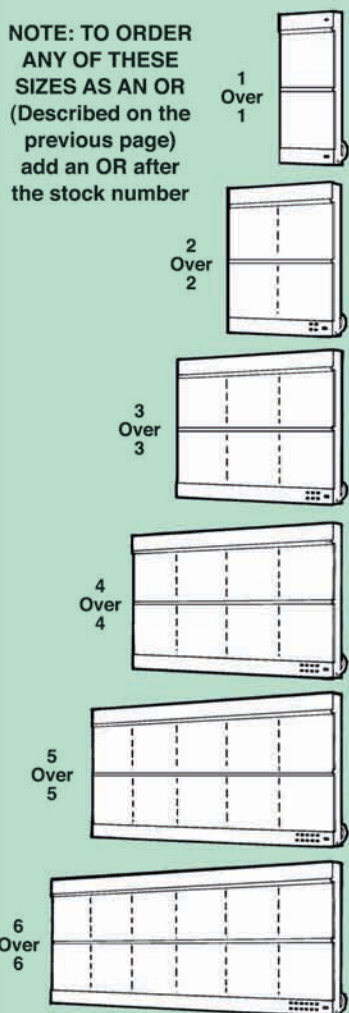

Standard	HI	Lo	Recess Kit	Viewing Area	Wall Mounted Dimensions A B C
22801	22901	21701		14" x 17" 35.6 x 43.2	22 14 3-1/4 55.9 35.6 8.3
22802	22902	21702		28" x 17" 71.2 x 43.2	22 28 3-1/4 55.9 71.2 8.3
22803	22903	21703		42" x 17" 106.8 x 43.2	22 42 3-1/4 55.9 106.8 8.3
22804	22904	21704		56" x 17" 142.4 x 43.2	22 56 3-1/4 55.9 142.4 8.3
22805	22905	21705		70" x 17" 178.0 x 43.2	22 70 3-1/4 55.9 178 8.3
22806	22906	21706		84" x 17" 213.6 x 43.2	22 84 3-1/4 55.9 213.6 8.3
22811	22911	21711		14" x 34" 35.6 x 86.4	44 14 3-1/4 111.8 35.6 8.3
22822	22922	21722		28" x 34" 71.2 x 86.4	44 28 3-1/4 111.8 71.2 8.3
22833	22933	21733		42" x 34" 106.8 x 86.4	44 42 3-1/4 111.8 106.8 8.3
22844	22944	21744		56" x 34" 142.4 x 86.4	44 56 3-1/4 111.8 142.4 8.3
22855	22955	21755		70" x 34" 178.0 x 86.4	44 70 3-1/4 111.8 178 8.3
22866	22966	21766		84" x 34" 213.6 x 86.4	44 84 3-1/4 111.8 213.6 8.3



Wolf X-Ray Illuminators comply with all ACR Recommendations.



THE TRIMLINE ILLUMINATORS

Wolf Manufactures Five Trimline Styles, Each With Unique Features Designed To Satisfy Your Individual Viewing Situation. Whatever Your Needs, There Is A Trimline Style To Fill Them.



Trimline Plus Hi Lo



The TrimlinePlus Hi-Lo provides two light intensities per each 14 x 17" viewing area. You can go from 5,150 to 7,250 Candela per square meter at the push of a button. Illumination is from six lamps in high intensity, four lamps in low intensity. Shown here is a 2 in 1 unit with the left hand unit switched to "Hi" to illustrate the difference between Hi and Lo.

The Trimline Plus Hi-Lo features a 3 1/4" deep cabinet with a baked on white finish, a roller grip film retainer, film activated auto switches (the illuminator goes on automatically when film is inserted), and cluster switches, on 2 over 2 units and up.

See page 9 for dimensions and stock numbers.

Trimline Plus



The TrimlinePlus offers four lamps (5,150 Candela per square meter) per each 14 x 17" viewing area.

The TrimlinePlus is available in 12 sizes (shown above is a 4 over 4 unit), or may be custom built to your specifications. Multi tier units are specially designed with only 1 1/2" of "dead space" between the top and bottom viewing panels.

The TrimlinePlus comes standard with a baked on white finish, film activated auto switches, a roller grip film retainer and cluster switches on 2 over 2 units and up.

See page 9 for dimensions and stock numbers.

Trimline III

Should your needs require Wolf offers the Trimline in a three bulb version (4,000 Candela per square meter of illumination).

The Trimline III has all the standard Trimline family features: viewing areas from 14 x 17" up to 2 tiers of 84 x 17", film activated auto switches, cluster switches (on 2 over 2 units and up) and a roller grip film retainer all in a sleek cabinet finished in white baked on enamel.

Options include a Mobile Stand, custom installations (see below) and a frame for recessed mounting.

As with all the Trimlines, the Trimline III can be combined with any other Trimline style to create a viewing area uniquely suited to your viewing needs.

See page 9 for dimensions and stock numbers.



All Trimlines Series Illuminators can be recessed flush in the wall for a sleek unobtrusive look. Simply order your Trimline with a recessing kit. The frame attaches easily in minutes with a screwdriver.


Trimline

The Trimline has 2 lamps (2,200 Candela per square meter) per each 14 x 17" viewing area. The Trimline is available in any one of 12 sizes with viewing areas from 14 x 17" up to two tiers of 84 x 17" each (shown is a 4 in 1 which gives you a viewing area of 56 x 17").

The Trimline comes standard with a baked on white finish, film activated auto switches, a roller grip film retainer and cluster switches on 2 over 2 units and up.

Options include an outer frame for recessed installation and mobile stand

See page 9 for dimensions and stock numbers..

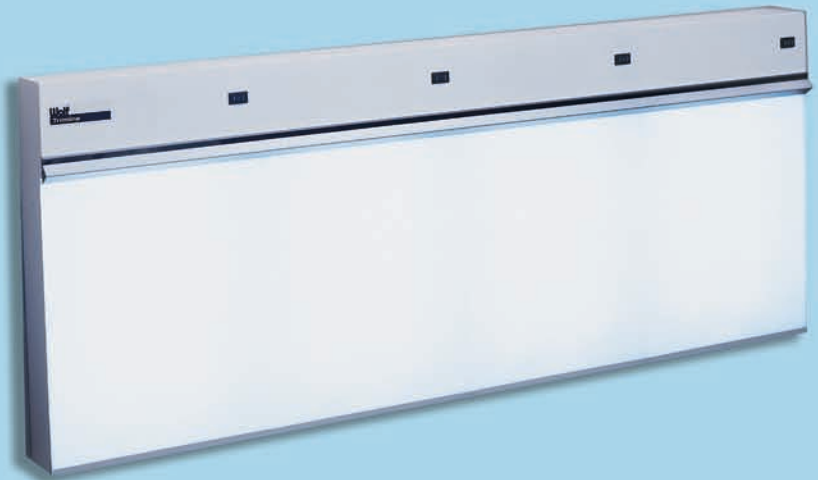

Trimline Basic




Wolf's Trimline Basic Illuminator is the "meat and potatoes" version of our Trimline Series. You don't get an automatic film activator, what you do get is Wolf quality and dependability.

The Trimline Basic is a 2 bulb illuminator (2,200 Candela per square meter) that is available in 12 sizes, in multiple sizes (2 in 1 through 6 over 6) each 14 x 17" viewing area has an individual switch. Master switch and cluster switches are on 2 over 2 units and up.

See page 9 for dimensions and stock numbers.


Wall Mounting Your Illuminator

All Trimline Series Illuminators come with Wolf's E-"Z" bar wall mounting system. Simply follow the steps below to quickly and accurately mount any size illuminator on the wall.

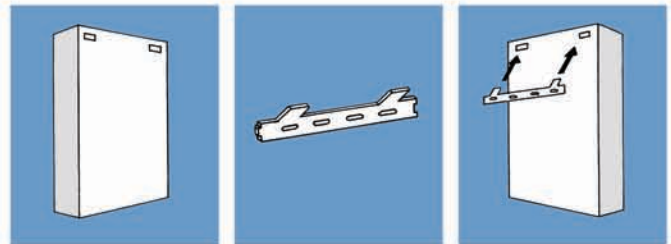

1. The back of your illuminator already has an E-"Z" bar taped to it. Remove the bar and note the 2 holes that E-"Z" Bar came out of.
2. Take the E-"Z" bar and using a spirit level, mount it on the wall with the flanges on the bar facing out as shown.
3. Slip the holes in your illuminator over and down onto the E-"Z" Bar flanges, and you're done.

Recessing Your Illuminator In The Wall

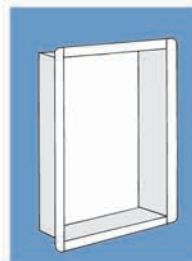

As shown on the picture of the Trimline Plus any illuminator can have a frame attached for installation flush in the wall. The frame attaches in minutes via pre-drilled holes and screws.

Note: Order the recessing kit when ordering your illuminator. Stock numbers for the kits are on page 9 next to the appropriate illuminator stock number.



Mobile stands are available for all Trimlines from 2 bank to 6 over 6. See the next page for Trimline Series stock numbers.

TrimlinePlus Hi-Lo™

Standard Features...

- Variable light intensities.
(6 bulbs dim to 4 bulbs) per viewing area.
- Film activated auto switches.
- Cluster switches and a master switch on 2 over 2 units and up.
- Gravity roller grip film retainer.
- Hospital grade plug.

Optional Features

- A frame for recessed mounting.

TrimlinePlus™

Standard Features...

- 4 bulbs per viewing area.
- Film activated auto switches.
- Cluster switches and a master switch on 2 over 2 units and up.
- Gravity roller grip film retainer.
- Hospital grade plug.

Optional Features

- A frame for recessed mounting.

Trimline III™

Standard Features...

- 3 bulbs per viewing area.
- Film activated auto switches.
- Cluster switches and a master switch on 2 over 2 units and up.
- Gravity roller grip film retainer.

Optional Features

- Hospital grade plug.
- A frame for recessed mounting.

Trimline

Standard Features...

- 2 bulbs per viewing area.
- Film activated auto switches.
- Cluster switches and a master switch on 2 over 2 units and up.
- Gravity roller grip film retainer.

Optional Features

- Hospital grade plug.
- A frame for recessed mounting.

Trimline Basic


Standard Features...

- 2 bulbs per viewing area.
- 3 prong plug.
- Gravity roller grip film retainer.

Optional Features

- Hospital grade plug.
- A frame for recessed mounting.

Trimline Series Specifications

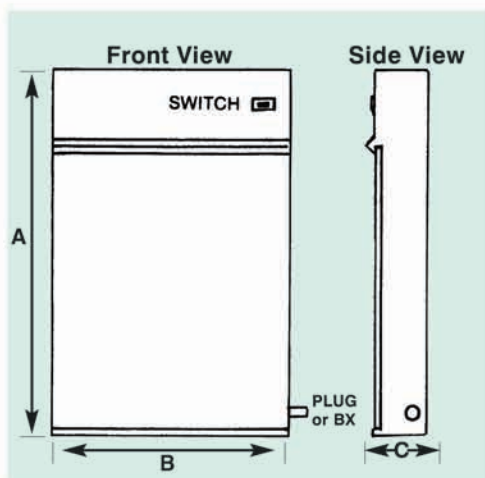


Wolf's Trimline Series Illuminators can be ordered in any of the twelve configurations shown in the right hand column.

The letters "A", "B" and "C" in the drawings above reference the corresponding letters in the "Wall Mounted Dimensions" column on page 9. Pick the size illuminator you want from one of the 12 shown and follow that line to the right to get the dimensions of the unit. Follow the line further right to get the correct stock number for that unit from your choice of one of the five Trimline Series illuminators available.


All five of the Trimline Series Illuminators come standard with Wolf's E-Z Bar for simple surface wall mounting. Simply mount the bar on the wall and slip the illuminator over it.


All the Trimlines shown can be installed flush in the wall by attaching an optional "recessing frame". Stock numbers for these frames are in the last column of the right hand page. As shown below, frames are easy to install.

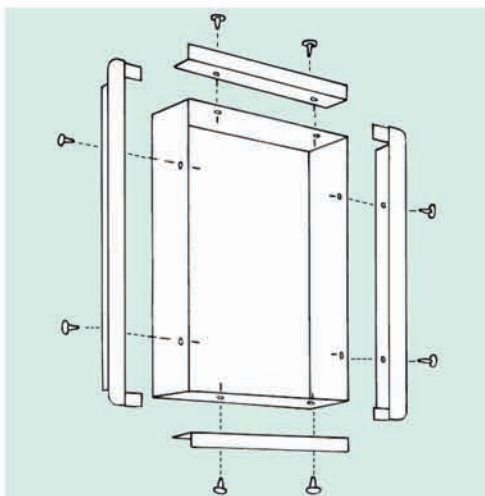

Using the screws and holes provided, screw in the 2 vertical frame sections first. Then screw in the 2 horizontal frame sections (the flange will overlap behind the horizontal sections.)


Trimline Series Available Sizes

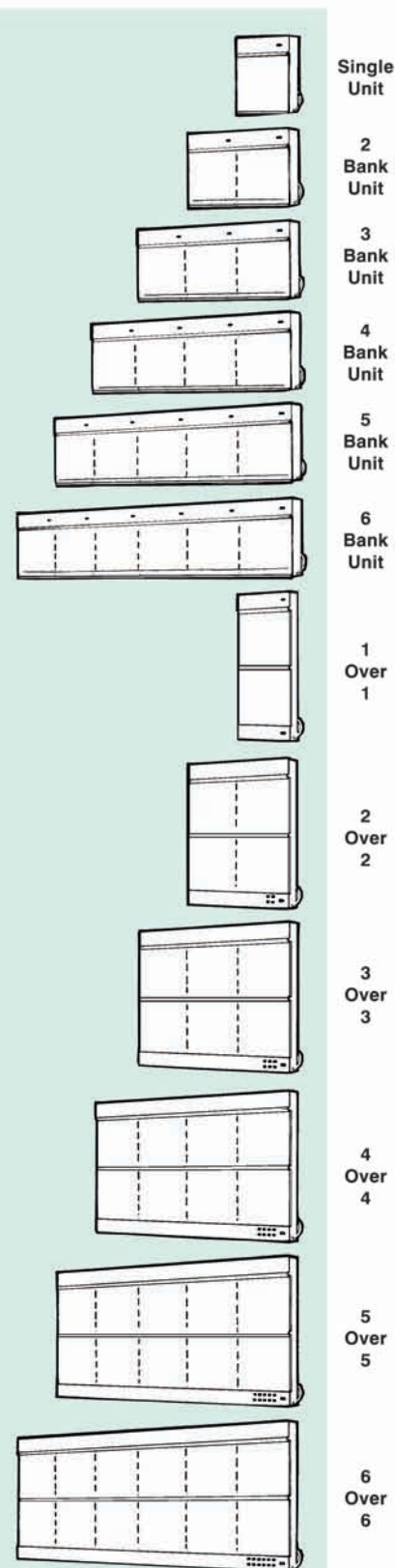

Viewing Area	Wall Mounted Dimensions			Stock Numbers for Trimline Plus HiLo	Stock Numbers for Trimline Plus	Stock Numbers for Trimline III	Stock Numbers for Trimline	Stock Numbers for Trimline Basic	Stock Numbers for Recessing Frames All Trimline
	A	B	C						
17" x 14" 43.2 x 35.6	22 55.9	14 35.6	3-1/4 8.3	29001 19 lbs. <i>Approx. shipping weight</i>	29401 19 lbs.	23401 17 lbs.	21401 17 lbs.	29501 17 lbs.	21701
17" x 28" 43.2 x 71.2	22 55.9	28 71.2	3-1/4 8.3	29002 35 lbs.	29402 35 lbs.	23402 30 lbs.	21402 30 lbs.	29502 30 lbs.	21702
17" x 42" 43.2 x 106.8	22 55.9	42 106.8	3-1/4 8.3	29003 57 lbs.	29403 54 lbs.	23403 54 lbs.	21403 54 lbs.	29503 52 lbs.	21703
17" x 56" 43.2 x 142.4	22 55.9	56 142.4	3-1/4 8.3	29004 70 lbs.	29404 67 lbs.	23404 67 lbs.	21404 67 lbs.	29504 64 lbs.	21704
17" x 70" 43.2 x 178.0	22 55.9	70 178	3-1/4 8.3	29005 120 lbs.	29405 120 lbs.	23405 120 lbs.	21405 120 lbs.	29505 118 lbs.	21705
17" x 84" 43.2 x 213.6	22 55.9	84 213.6	3-1/4 8.3	29006 143 lbs.	29406 143 lbs.	23406 143 lbs.	21406 143 lbs.	29506 141 lbs.	21706
34" x 14" 86.4 x 35.6	44 111.8	14 35.6	3-1/4 8.3	29011 35 lbs.	29411 35 lbs.	23411 35 lbs.	21411 35 lbs.	29511 35 lbs.	21711
34" x 28" 86.4 x 71.2	44 111.8	28 71.2	3-1/4 8.3	29022 72 lbs.	29422 72 lbs.	23422 72 lbs.	21422 72 lbs.	29522 72 lbs.	21722
34" x 42" 86.4 x 106.8	44 111.8	42 106.8	3-1/4 8.3	29033 130 lbs.	29433 130 lbs.	23433 130 lbs.	21433 130 lbs.	29533 128 lbs.	21733
34" x 56" 86.4 x 142.4	44 111.8	56 142.4	3-1/4 8.3	29044 160 lbs.	29444 140 lbs.	23444 140 lbs.	21444 140 lbs.	29544 138 lbs.	21744
34" x 70" 86.4 x 178.0	44 111.8	70 178	3-1/4 8.3	29055 180 lbs.	29455 180 lbs.	23455 180 lbs.	21455 180 lbs.	29555 178 lbs.	21755
34" x 84" 86.4 x 213.6	44 111.8	84 213.6	3-1/4 8.3	29066 210 lbs.	29466 210 lbs.	23466 210 lbs.	21466 210 lbs.	29566 208 lbs.	21766
 Wolf X-Ray Illuminators comply with all ACR Recommendations.									


Trimline Plus Hi-Lo

- 6 lamps 15 Watt Daylight F15T8D. Hi/Lo switch alters illumination level by activating 4 lamps in the low setting, 6 lamps in the high setting. Provides an increase in illumination between high, and low settings of 50% (typical luminous intensities are 7250 to 5150 candela per meter squared), with uniform light distribution.
- 6 Solid State Electronic Ballasts with inherent thermal protection. A sound rating, max. 1.50 Amps, 120VAC, 60/50 Hz, UL, CSA listed. Also available in 220VAC, 60/50 Hz.
- Three-position rocker switch controls individual viewing panels on single through 6 in 1 bank units. All tier units come standard with cluster switch controls for individual Hi-Lo/On-Off Automatic film activated panel viewing and Master On/Off switch.
- Film activated switch (Automatic Model) conserves energy, eliminates distracting light from adjacent panels when not in use.
- Wolf's unique E-"Z" Bar™ wall mounting system makes hanging your illuminator quick and accurate.
- Film held with roller gravity grip. Rollers are coated to prevent scratching the film.
- A hospital grade plug.

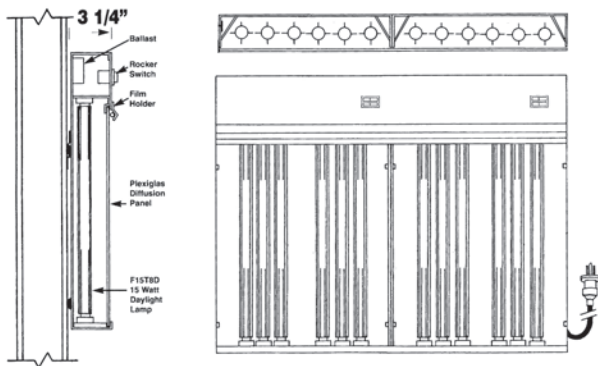


Trimline Plus

- 4 lamps 15 Watt Daylight F15T8D Illuminator.
- 4 Solid State Electronic Ballasts with inherent thermal protection. A sound rating, max. 1.00 Amps, 120VAC, 60/50 Hz, UL, CSA listed. Also available in 220VAC, 60/50 Hz.
- Three-position rocker switch controls individual viewing panels on single through 6 in 1 bank units. All tier units come standard with cluster switch controls for individual On-Off Automatic film activated panel viewing and Master On/Off switch.
- Film activated switch (Automatic Model) conserves energy, eliminates distracting light from adjacent panels when not in use.
- Wolf's unique E-"Z" Bar™ wall mounting system makes hanging your illuminator quick and accurate.
- Film held with roller gravity grip. Rollers are coated to prevent scratching the film.
- A hospital grade plug.

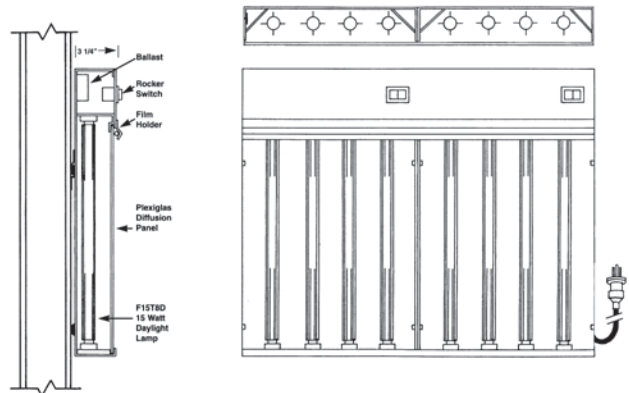


Trimline III

- 3 lamps 15 Watt Daylight F15T8D Illuminator.
- 3 Solid State Electronic Ballasts with inherent thermal protection. A sound rating, max. .75 Amps, 120VAC, 60/50 Hz, UL, CSA listed. Also available in 220VAC, 60/50 Hz.
- Three-position rocker switch controls individual viewing panels on single through 6 in 1 bank units. All tier units come standard with cluster switch controls for individual On-Off Auto film activated panel viewing and Master On/Off switch.
- Film activated switch (Automatic Model) conserves energy, eliminates distracting light from adjacent panels when not in use.
- Wolf's unique E-"Z" Bar™ wall mounting system makes hanging your illuminator quick and accurate.
- Film held with roller gravity grip. Rollers are coated to prevent scratching the film.

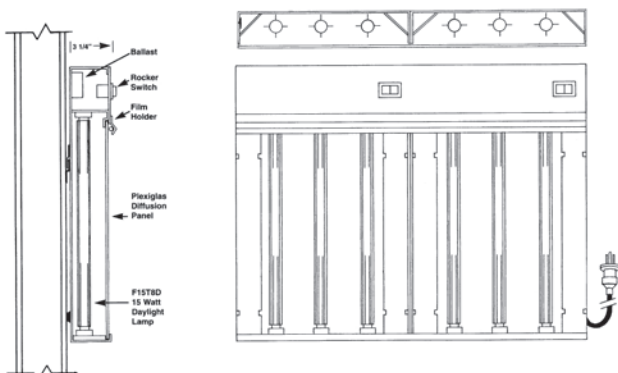

Trimline

- 2 lamps 15 Watt Daylight F15T8D Illuminator.
- 2 Solid State Electronic Ballasts with inherent thermal protection. A sound rating, max. .50 Amps, 120VAC, 60/50 Hz, UL, CSA listed. Also available in 220VAC, 60/50 Hz.
- Three-position rocker switch controls individual viewing panels on single through 6 in 1 bank units. All tier units come standard with cluster switch controls for individual On-Off Automatic film activated panel viewing and Master On/Off switch.
- Film activated switch (Automatic Model) conserves energy, eliminates distracting light from adjacent panels when not in use.
- Wolf's unique E-"Z" Bar™ wall mounting system makes hanging your illuminator quick and accurate.
- Film held with roller gravity grip. Rollers are coated to prevent scratching the film.

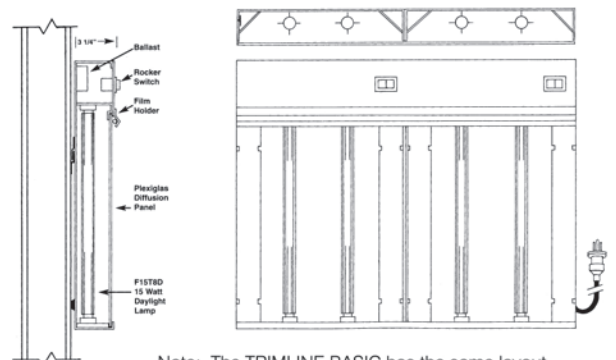

Note: The TRIMLINE BASIC has the same layout as the Trimline and the same "2 lamps" specs, most other specs are different, see page 7 for details.


MG-7

Wolf's workhorse, the MG-7 is designed to give the private practitioner years of reliable use.

Easily wall surface mounted, the MG-7's five inch depth also makes it ideal for desk top use.

Available as a single unit with a 14 x 17" viewing area or up to 4 in 1 (56 x 17" viewing area), the MG-7 features two 15 watt fluorescent lamps (2,200 Candela per square meter), a white baked on enamel finish, magic grip film retainer and a 6' line cord with 3 prong plug.

Specifications and ordering information are on the next page.

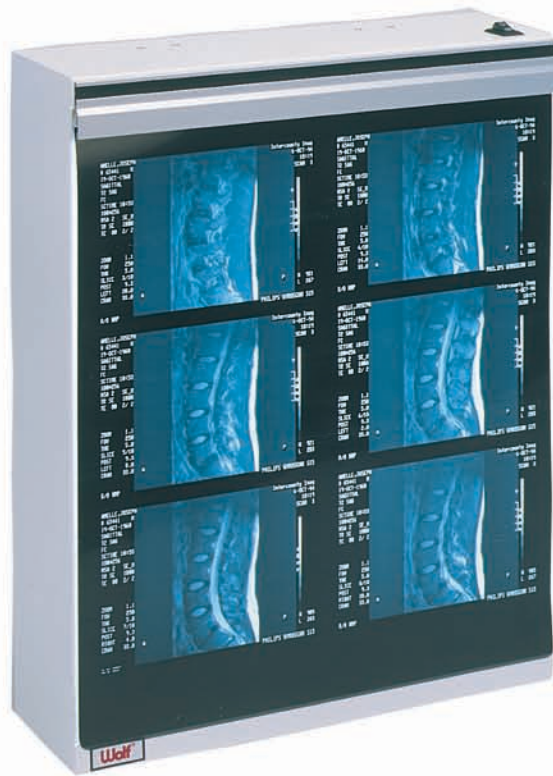

MG-7 and Econoline Illuminators

Wolf's durability
and reliability
in an inexpensive
format

Econoline

The Econoline is Wolf's most economical illuminator. Available in 4 sizes, the Econoline is available with a kit for recessed installation (page 18).


The Econoline features two bulb illumination, (15 watts each – 2,200 Candela per square meter), a 3 1/2" deep cabinet with a white finish, a magic grip film retainer and a 6 foot line cord with 3 prong plug.

Specifications and ordering information are on the next page.



Wolf
Wolf
Wolf
Wolf
Wolf
Wolf
Wolf
Wolf
Wolf
Wolf

MG-7 Series™ Specifications

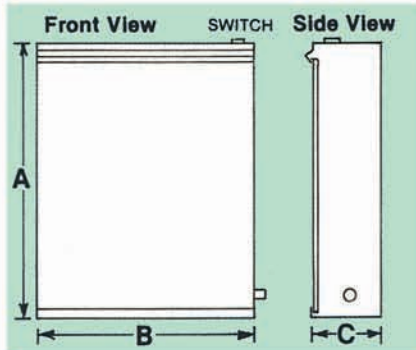

Standard Features...

- 2 bulbs per viewing area.
- A 3 prong plug.
- Magic grip film retainer.

Optional Features

- Handle kit.
- Can be used on a desk top

Description


	Stock Number
Single Unit	21104
2 Bank Unit	21801
3 Bank Unit	21802
4 Bank Unit	21803

**MODEL CANNOT
BE RECESSED**



Viewing Area	Wall Mounted Dimensions			Approx. Shipping Weight
	A	B	C	
14" x 17" 35.6 x 43.2	18-1/8 46.0	14-1/8 35.9	5 12.7	13 lbs.
28" x 17" 71.2 x 43.2	18-1/8 46.0	28-1/4 71.8	5 12.7	29 lbs.
42" x 17" 106.8 x 43.2	18-1/8 46.0	42-3/8 107.7	5 12.7	49 lbs.
56" x 17" 142.4 x 43.2	18-1/8 46.0	56-1/2 143.6	5 12.7	60 lbs.

Econoline Series™ Specifications

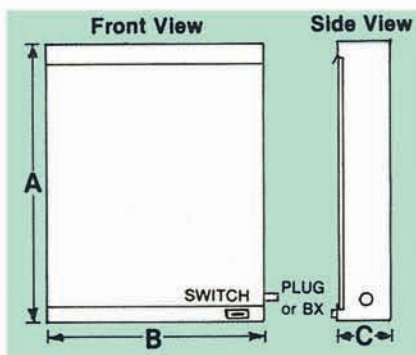

Standard Features...

- 2 bulbs per viewing area.
- A 3 prong plug.
- Magic grip film retainer.

Optional Features

- A frame for recessed mounting.
- A handle mounted on top of the box for easy transport

Ordering Information

- Use stock numbers to place your order
- If you want a recessed illuminator, order the appropriate Recessing Kit. Stock numbers for Recessing Kits are on page 18.

Description

	Stock Number	Recess Kit
Single Unit	29601	29701
2 Bank Unit	29602	29702
3 Bank Unit	29603	29703
4 Bank Unit	29604	29704

Viewing Area	Wall Mounted Dimensions			Approx. Shipping Weight
	A	B	C	
14" x 17" 35.6 x 43.2	18-1/4 46.0	14 35.9	3-1/2 9.0	10 lbs.
28" x 17" 71.2 x 43.2	18-1/4 46.0	28 71.8	3-1/2 9.0	20 lbs.
42" x 17" 106.8 x 43.2	18-1/4 46.0	42 107.7	3-1/2 9.0	30 lbs.
56" x 17" 142.4 x 43.2	18-1/4 46.0	56 143.6	3-1/2 9.0	40 lbs.



**Wolf X-Ray Illuminators
comply with all ACR
Recommendations.**

All Trimline and Liberator Illuminators can be custom designed to fit any area. The unit shown right is a good example of the kind of viewing environment Wolf engineers can create for you. This unit contains 16 illuminators and features ... individual automatic film activated on/off switches, individual 14" x 17" viewing area switches in a console panel, a master switch for each entire wing, and a combination of Trimline and TrimlinePlus units. Call the Wolf (800) number for more information, and a portfolio of viewing areas we've designed.




Mammoline, Trimline Plus Hi-Lo, Trimline Plus, Trimline, Trimline III, Trimline Basic, The Liberator and OR Series

Wolf's custom manufactured mobile stands provide attractive desk-style seated viewing in any location. Manufactured of high grade welded steel, finished in white baked enamel, Wolf's mobile stands ride on easy glide, ball bearing casters with positive locking mechanism. Available for the Mammoline all Trimlines, The Liberator and the OR Series in 2 through 6 bank and 2 over 2 through 6 over 6 sizes. Call for dimensions and weights.

Trimline Desk Stand

If you're using your Trimline or Liberator Viewbox on a desk, then you need our new Trimline Desk Stand! The Trimline Desk Stand provides great stability and prevents rocking during film loading. The Trimline Desk Stand works with all Liberator models and

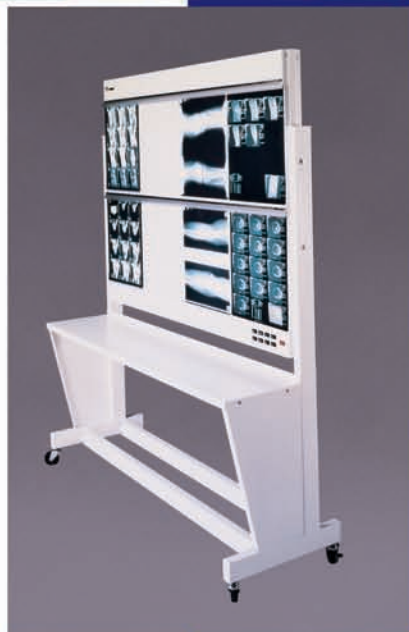

You can Save 50% on replacement bulbs by purchasing Wolf Replacement Bulbs ...and remember Wolf bulbs fit any brand illuminator. #21104-19



Here's a handy innovative option for your Trimline Illuminator. Available for single and double units, any style Trimline can be equipped with a handy, unobtrusive tray to hold pens, pencils or markers. The tray easily attaches to the illuminator using the pre-drilled holes for the frame. Single **#21508**, Double **#21509**



The Mammoline Mammography Illuminator

The Mammoline combines superior lighting technology with ACR mandated masking shutters for optimum film viewing.

At top intensity the Mammoline produces over 7,600 Candela per square meter of illumination! You can also adjust the light output down to 280 Candela via a built in dimmer (the ACR recommends a minimum of 3,400 Candela per square meter per viewing area).

Wolf's unique reflector design disperses light evenly to eliminate "hot spots" found on other viewers. The ballasts are state of the art energy efficient solid state units for flicker free instant on lighting. In addition, the Mammoline has a unique venting design that dissipates heat build-up. The Mammoline's shutters are completely light proof and glide easily on horizontal and vertical tracks. Film grips are gravity roller type for easy film mounting.

The Mammoline is easily wall mounted by means of Wolf's unique E-"Z" Bar mounting system. An optional mobile stand is also available. **(Stock No. 21504MS-28)** The Mammoline cabinet is 48" long, 36" high and 5" deep. The viewing area is two tiers measuring 38 1/2" x 11" each. The finish is baked on enamel. The approximate shipping weight for wall hung units is 146 lbs.

Mammoline Illuminator Stock No. 28601
The Mammoline cannot be recessed.

VuPlus Mammography Viewer



The VuPlus is a mammography viewer with all the special features mammography imaging requires. The viewing area is custom sized to 29 x 48 1/2 cm to minimize light spillover and the unit has two light intensities.

The VuPlus features a white baked on enamel finish, a gravity roller grip film retainer, five fluorescent lamps (3,426 candela per square meter) per viewing area, a rocker switch and an 8' line cord with 3 prong plug.

Stock Number Dimensions

28001 (single unit) 16 3/8" H x 19 1/8" W x 4"D

28002 (double unit) 16 3/8" H x 39 1/8" W x 4"D

Approximate shipping weight (single unit) 19 lbs.

The VuPlus cannot be recessed.



PLEASE NOTE: This Product Has Been Discontinued

The Wolf OR Illuminators

The Wolf OR Illuminator Series have been designed to meet the requirements of more restrictive areas such as operating rooms and as such are listed to a more stringent UL standard.


- The OR Plus, a four bulb illuminator with the same features as the Trimline Plus described on page 6.
- The OR Plus Hi-Lo, an illuminator that goes from 4 to 6 bulb intensity at the press of a button as described on page 5.
- The OR Liberator Plus, all the features of the OR Plus with Wolf's uniquely removable viewing area, described on page 3.
- The OR Liberator Plus Hi-Lo, like the OR Plus Hi-Lo described above (4 to 6 bulb intensity) and on page 5 with our "Liberator" easy-off plexi surface as described on page 3.

All are available in the 12 sizes shown at left and feature your choice of stainless steel or white finish. Standard features include film activated auto switches, a roller grip film retainer, cluster switches on 2 over 2 units and up and a hospital grade plug. Options include an outer frame for recessed installation (as shown), and a mobile stand. Specifications and ordering information for the ORs are on the following page.

	OR Plus	OR Hi-Lo	OR Liberator Plus	OR Liberator Plus Hi-Lo
Single	29401OR	29001OR	22801OR	22901OR
2 Banks	29402OR	28002OR	22802OR	22902OR
3 Banks	29403OR	28003OR	22803OR	22903OR
4 Banks	29404OR	28004OR	22804OR	22904OR
5 Banks	29405OR	28005OR	22805OR	22905OR
6 Banks	29406OR	28006OR	22806OR	22906OR
1 OVER 1	29411OR	29011OR	22811OR	22911OR
2 OVER 2	29422OR	29022OR	22822OR	22922OR
3 OVER 3	29433OR	29033OR	22833OR	22933OR
4 OVER 4	29444OR	29044OR	22844OR	22944OR
5 OVER 5	29455OR	29055OR	22855OR	22955OR
6 OVER 6	29466OR	29066OR	22866OR	22966OR\

Wolf
Wolf
Wolf
Wolf
Wolf

OR Specifications

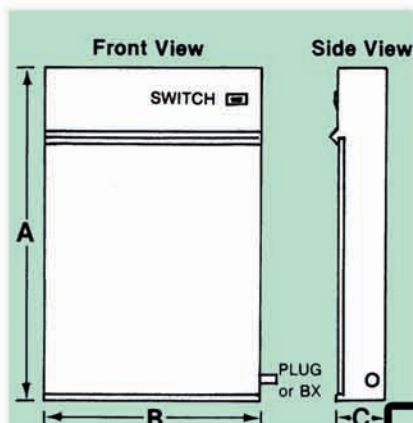

Letters reference dimensions shown in the right hand column.

Standard Features...

- Film activated auto switches.
- Cluster switches and a master switch on 2 over 2 units and up.
- Gravity roller grip film retainer.
- Wolf's E-"Z" Bar™ for easy surface wall mounting. Simply mount the bar on the wall and slip the illuminator over it.
- Hospital grade plug.

Optional Features


- A frame for recessed mounting.
- A mobile stand.
- Custom room installations
- A stainless steel finish.

Ordering Information


- Use stock numbers to place your order
- If you want a recessed illuminator, order the appropriate Recessing Kit. Stock numbers for Recessing Kits are on page 18.
- An "S" should be placed after the stock # for a stainless steel finish.

Call Wolf X-Ray For Pricing on Stainless Steel Units.

Description

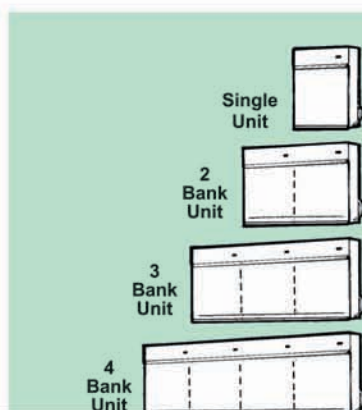

PLEASE NOTE: This Product Has Been Discontinued

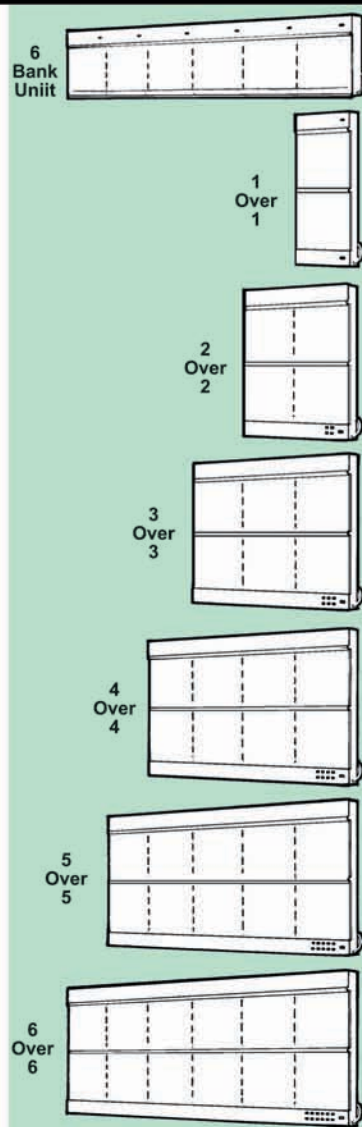

Viewing Area	Wall Mounted Dimensions			Approx. Shipping Weight
	A	B	C	
14" x 17" 35.6 x 43.2	22 55.9	14 35.6	3-1/4 8.3	19 lbs.
28" x 17" 71.2 x 43.2	22 55.9	28 71.2	3-1/4 8.3	32 lbs.
42" x 17" 106.8 x 43.2	22 55.9	42 106.8	3-1/4 8.3	65 lbs.
56" x 17" 142.4 x 43.2	22 55.9	56 142.4	3-1/4 8.3	70 lbs.
70" x 17" 178.0 x 43.2	22 55.9	70 178	3-1/4 8.3	120 lbs.
84" x 17" 213.6 x 43.2	22 55.9	84 213.6	3-1/4 8.3	143 lbs.
14" x 34" 35.6 x 86.4	44 111.8	14 35.6	3-1/4 8.3	30 lbs.
28" x 34" 71.2 x 86.4	44 111.8	28 71.2	3-1/4 8.3	72 lbs.
42" x 34" 106.8 x 86.4	44 111.8	42 106.8	3-1/4 8.3	130 lbs.
56" x 34" 142.4 x 86.4	44 111.8	56 142.4	3-1/4 8.3	140 lbs.
70" x 34" 178.0 x 86.4	44 111.8	70 178	3-1/4 8.3	180 lbs.
84" x 34" 213.6 x 86.4	44 111.8	84 213.6	3-1/4 8.3	210 lbs.



Wolf X-Ray Illuminators comply with all ACR Recommendations.

Podo Lite

The 12 x 24" viewing area makes the Podo Lite perfect for the Podiatrist.

The Podolite can be desk mounted, wall hung or recessed in the wall. The Podo Lite features steel construction, a white baked on enamel finish, Magic Grip film retainer, two fluorescent lamps (20 watts each) a rocker switch and a 6' line cord with 3 prong plug.

Specifications: 24 1/8" long, 13 1/2" high, 3 1/2" deep # **21506-N**. Note make the hole cut in wall for recessed units 24 1/4" long by 13 3/4" high and order the Recessing Kit: #**27506**



Wolf Intensi-Spot

Designed to provide a concentrated light source. The Intensi-Spot utilizes a 75-watt lamp, shielded by an octagon mask for light concentration.

The Intensi-Spot tilts to any angle and includes foot switch operation. The unit comes in a baked enamel finish, with holes in the base for wall mounting, and a 6' line cord with 3 prong plug.

Specifications: 8 1/2"H x 6"D.

Approx. shipping weight: 5 lbs.

Stock Number **21295**



Wolf Cool-Lite

This unique viewing device provides concentrated light like our intensi-spot, but eliminates the heat by the use of a fluorescent bulb. Energy efficient, the Cool-lite can be left on all day.

The Cool-Lite tilts to any angle, has a baked enamel finish, and includes holes in the base for easy wall mounting. The Cool-Lite comes with a 6' line cord with a 3 prong plug.

Specifications: 8 1/2"H x 6"D.

Approx. shipping weight: 5 lbs.

Stock Number **21695**



Magnifying Glasses

Wolf's four inch diameter (2x) magnifying glass (left hand picture) is perfect for close up diagnostic film inspection. It has a polished aluminum rim and a four inch black plastic handle. # **21604**

The 5-inch magnifier (right hand photo) magnifies 2x. It complies with all accepted procedures for mammography and general use. # **21602**

Wolf
Wolf
Wolf
Wolf
Wolf
Wolf
Wolf
Wolf
Wolf
Wolf

36 Inch and 51 Inch Spinal Illuminators

The 36 inch is available as a Trimline 36 (3 bulbs per 14 x 36" viewing area) or a Trimline Plus 36 (4 bulbs per 14 x 36" viewing area), the bulbs are 30 watts each.

Both feature roller grip film retainers, an 8' line cord with 3 prong plug and a rocker switch. An E-"Z" Bar™ is provided for easy wall mounting. Unit dimensions are 14" wide, 42" long, 3 1/4" deep. The unit can be recessed flush in the wall by attaching a recessing frame, stock numbers shown below.

Trimline 36	
Single	#29280
Recessing Frame	#29780
Double	#29281
Recessing Frame	#29781
Trimline Plus 36 *	
Single	#28280
Recessing Frame	#29780
Double	#28281
28" wide, 42" long, 4" deep	
Recessing Frame	#29781

*4 bulb spinals come with an auto switch

The 51 inch is available as a Trimline 51 (3 bulbs per 14 x 51" viewing area) or a Trimline Plus 51 (4 bulbs per 14 x 51" viewing area), the bulbs are 30 watts each, with an auto switch.

Both feature roller grip film retainers, an 8' line cord with 3 prong plug and a rocker switch. An E-"Z" Bar™ is provided for easy wall mounting. Unit dimensions are 14" wide, 55 3/4" long, 3 1/4" deep. The unit can be recessed flush in the wall by attaching a recessing frame.

Trimline 51	
Single	#29282
Recessing Frame	#29782
Trimline 51 Plus *	
Single	#28282
Recessing Frame	#29782



Wolf Illuminator Recessing Kit


To recess your illuminator, you must order the appropriate frame. Stock numbers for the frames are shown below. A single illuminator gets a single frame, a 2 over 2 illuminator gets a 2 over 2 frame, etc.

The following stock numbers are for: The Liberator, Liberator HiLo, Trimline Plus Hi-Lo, Trimline Plus, Trimline III, Trimline Series, Trimline Basic and the OR Series.

Single - 21701	1 over 1 - 21711
2 in 1 - 21702	2 over 2 - 21722
3 in 1 - 21703	3 over 3 - 21733
4 in 1 - 21704	4 over 4 - 21744
5 in 1 - 21705	5 over 5 - 21755
6 in 1 - 21706	6 over 6 - 21766


These numbers are for the Econoline.

Single - 29701
2 in 1 - 29702
3 in 1 - 29703
4 in 1 - 29704

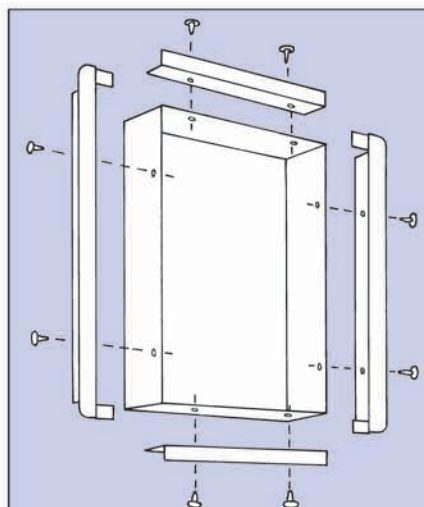

Using the screws and holes provided, screw in the 2 vertical frame sections first. Then screw in the 2 horizontal frame sections (the flange will overlap behind the horizontal sections).

□
□
□

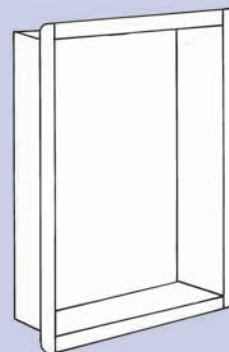

Your illuminator will now look like this, and is ready to install in the wall.

□