



Hu-Friedy's

STERILIZATION CONTAINER SYSTEM

Provides Efficient Sterilization of Cassettes and Instruments

Hu-Friedy's Sterilization Containers are a reusable, rigid system used for the packaging of cassettes and instruments during sterilization. An alternative to the wrap and pouch materials commonly used to package cassettes and instruments, Sterilization Containers provide significant economic and environmental benefits along with enhanced safety.

- Allows sterilization of enclosed contents and maintains sterility during transport and storage, until opened.
- Provides fully enclosed and locked transit of loose instruments and cassettes.
- Offers economic and environmental benefits by providing reusable packaging, low cost of use, a long service life, and minimal use of disposable materials.
- Uniquely sized to fit many Hu-Friedy cassettes, including the Container Cassette that features an integrated handle for easy insertion and removal.
- Includes a full line of Container accessories including filters, indicator cards, and tamper evident seals.
- Suitable for use in pre-vacuum Steam Sterilizers.



POINTS OF PERFORMANCE

Two external display windows provide space for a process indicator card and container identification

The filter retention plates audibly click to lock into place providing easy filter changes

Tamper evident seals on latches provide visual evidence of sterility breach


Lid and bottom are made of anodized aluminum making them lightweight with excellent drying properties. The aluminum alloy with an anodized oxide surface helps resist corrosion.

Latches on the lid act as handles for easy, contamination free lid removal

How the best perform




STERILIZATION CONTAINER SYSTEM

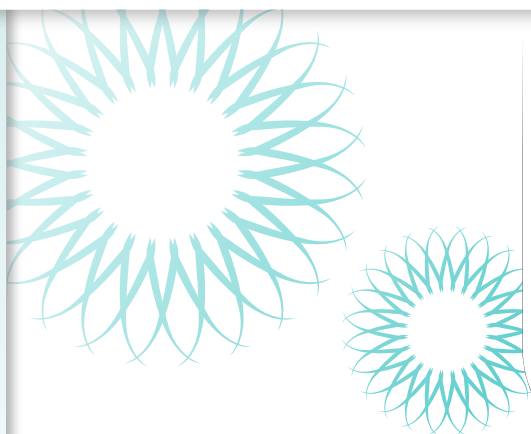

Container Cassette for 12 Instruments

Blue
268 x 176 x 34 mm | **IMECC128**
10-1/2 x 7 x 1-1/3"

Container Cassette with handle on lid and filter outlets.
Fits 12 instruments, with an accessory area.
Fits all sterilization containers

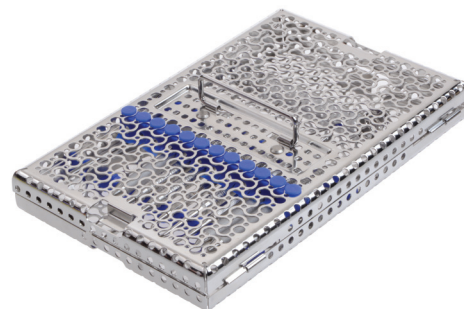

312 x 190 x 50 mm
12-1/4 x 7-1/2 x 2" | **IMDINCO0M**

312 x 190 x 65 mm
12-1/4 x 7-1/2 x 2-1/2" | **IMDINCO1M**

Top and bottom with filter.
Fits one container cassette.



312 x 190 x 100 mm
12-1/4 x 7-1/2 x 4" | **IMDINCO2M**

Top and bottom with filter.
Fits two container cassettes.



312 x 190 x 135 mm
12-1/4 x 7-1/2 x 5-1/3" | **IMDINCO3M**

Top and bottom with filter.
Fits three container cassettes.

CONTAINER ACCESSORIES

Container Indicator Card
250 pieces, single use | **IMCO-INDNA**

Tamper Evident Seal for Latch
100 pieces, single use | **IMCOM-SEAL**

Disposable Paper Filter
100 pieces, blank, single use | **IMCO-FIL3M**



Tamper Evident Seal for Latch



Container Indicator Card

Visit us online at Hu-Friedy.com/SterilizationContainers

Call 1-800-Hu-Friedy or contact your authorized Hu-Friedy representative for more information.
©2015 Hu-Friedy Mfg. Co., LLC. All rights reserved. IMS-372/0315

