

Ketac™ Fil Plus Aplicap™

Glass Ionomer Filling Material



The esthetic glass ionomer solution you can count on for restoring a confident smile.


- Chemical bond to enamel and dentin guarantees a tight marginal seal
- Tooth-like coefficient of thermal expansion
- Aplicap capsule delivery for consistent handling characteristics and reproducible results
- Minimal preparation
- High fluoride release
- Radiopaque
- Available in eight esthetic shades

Choose Ketac Fil Plus material to closely match natural tooth structure.

Ketac Fil Plus filling material is a conventional glass ionomer that helps you achieve virtually invisible anterior restorations.

There are eight shades to choose from, making it easy to match natural tooth color.

And, with the chemical bond to enamel and dentin, less tooth structure needs to be removed during preparation—which supports a healthy smile.

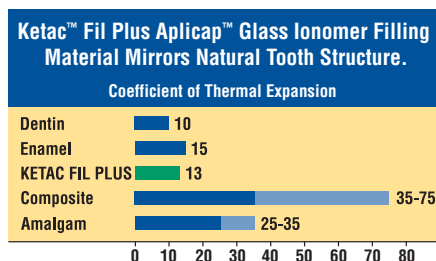

Source: Internal data 3M ESPE AG

Clinical advantages:

- Ideal for pediatric and geriatric patients
- Can be used to treat lost or cervical restorations
- Less shrinkage than traditional composites—leads to better marginal integrity
- Unique capsule design assures proper activation—so time-wasting, costly “duds” are not a problem

Indications:

- Class III and V restorations
- Wedge-shaped defects
- Small Class I fillings
- Core build-ups
- Minimally Invasive Dentistry (MID)

Glass Ionomer solutions from 3M ESPE produce the results you depend on.

After more than 30 years developing glass ionomer products, 3M ESPE has you covered. Use this handy reference to quickly determine which material is best suited to meet your procedural needs.

Indications									
	Small Class I	Class III	Class V	Minimal Invasive Dentistry (MID)	Primary Teeth	Base Under Composite	Temporary Treatment	Core Build-up	Delivery
RESIN MODIFIED									
Ketac™ Nano	●	●	●	●	●	●	●	●	
Photac™ Fil	●	●	●	●	●		●		
Vitremer™	●	●	●	●	●	●	●	●	
CONVENTIONAL									
Ketac™ Molar Quick	●		●	●	●	●	●	●	
Ketac™ Molar	●		●		●	●	●	●	
Ketac™ Fil Plus	●	●	●	●	●			●	
METAL REINFORCED									
Ketac™ Silver	●		●		●		●	●	



Quick Mix Capsule



Clicker™ Dispenser



Capsule



Powder Liquid

Ordering Information



Ketac™ Fil Plus Aplicap™ Introductory Kit (55000)

3M™ ESPE™ RotoMix™ Capsule Mixing Unit (76310)



Item Number	Product Information
55000	Ketac™ Fil Plus Aplicap™ GI Anterior Restorative Intro Pack Contains: 50 Capsules, asst. shades, 20 each shade A2, 10 each shades A3, A3.5, A4; 1 Activator; 1 Applier; 10mL Ketac™ Conditioner; 2.5mL Ketac™ Glaze.
55010	Ketac™ Fil Plus Aplicap™ GI Anterior Restorative — 50 Capsules, Assorted Contains: 50 Capsules, asst. shades, 20 each shade A2, 10 each shades A3, A3.5, A4
55020	Ketac™ Fil Plus Aplicap™ GI Anterior Restorative — 50 Capsules, Shade A1
55030	Ketac™ Fil Plus Aplicap™ GI Anterior Restorative — 50 Capsules, Shade A2
55040	Ketac™ Fil Plus Aplicap™ GI Anterior Restorative — 50 Capsules, Shade A3
55050	Ketac™ Fil Plus Aplicap™ GI Anterior Restorative — 50 Capsules, Shade A3.5
55060	Ketac™ Fil Plus Aplicap™ GI Anterior Restorative — 50 Capsules, Shade A4
55070	Ketac™ Fil Plus Aplicap™ GI Anterior Restorative — 50 Capsules, Shade C2
55080	Ketac™ Fil Plus Aplicap™ GI Anterior Restorative — 50 Capsules, Shade C3
55090	Ketac™ Fil Plus Aplicap™ GI Anterior Restorative — 50 Capsules, Shade D60
37190	Ketac™ Glaze, 2.5g Bottle
37470	Ketac™ Conditioner, 10mL
76310	3M™ ESPE™ RotoMix™ Capsule Mixing Unit At the touch of a button, the RotoMix Capsule Mixing Unit provides a void-free, completely homogenous mix which is easier to handle so you do your best work. The unique design of 3M™ ESPE™ Aplicap™ and Maxicap™ capsules prevents "duds" — saving your practice time and money.

Customer Care Center: 1-800-634-2249 www.3MESPE.com

3M ESPE

Dental Products

3M Center
Building 275-2SE-03
St. Paul, MN 55144-1000
U.S.A.
1-800-634-2249

3M Canada

Post Office Box 5757
London, Ontario N6A 4T1
Canada
1-888-363-3685

Please recycle. Printed in U.S.A.
© 3M 2011. All rights reserved.
70-2009-3714-5

3M, ESPE, Aplicap, Ketac, Maxicap, Photac and Vitremer are trademarks of 3M or 3M ESPE AG. Used under license in Canada.