


# ProCart II™

## Mobile Treatment Console

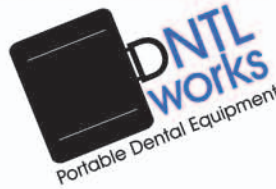

**T**HE INNOVATIVE PROCART II COMBINES THE BEST ATTRIBUTES OF SELF-CONTAINED MOBILITY AND A HIGH-TECH DESIGN THAT PERFECTLY SUITED FOR PROFESSIONAL IN-OFFICE USE.

### FEATURES INCLUDE:

- ✿ *Sophisticated, one-piece design: completely self-contained with built-in air compressor, separate vacuum pump, water and waste management.*
- ✿ *High-strength aluminum chassis with innovative WhisperQuiet sound reduction.*
- ✿ *Separate oil-free adjustable vacuum pump for powerful suction.*
- ✿ *Large, four liter waste container with automatic overflow shutoff and electrical waste purge.*
- ✿ *Ample, 1.5 liter water reservoir with additional external water connection.*



*Designed specifically for use in the modern dental office, nursing home, hospital and medical center, the ProCart II combines superior performance with state-of-the-art design.*



# PRO CART II™

## Mobile Treatment Console



### Specifications

#### *DIMENSIONS/WEIGHT*

Width: 20" (50.8 cm)

Depth: 21" (53.3 cm)

Height: 30" to 34" adjustable  
(76.2 cm 86.4 cm)

Weight: 100 lbs. (45.5 kg)

#### *ELECTRICAL*

Power Supply: 115 Vac, 60 Hz or  
220 Vac, 50 Hz

Circuit breaker protection: 15 amps  
Hospital Grade, NEMA type  
5-15 Plug

#### *COMPRESSOR*

Oil-free 1/3 hp motor  
Working pressure: 70 to 100 psi  
Flow rate: 2.5 cfm  
Sound level 48 dB  
Air reservoir: 4 liters  
Air tank purge

#### *VACUUM MOTOR*

Oil-free 1/3 hp motor  
Vacuum: 11" Hg  
Vacuum flow: 4.6 cfm  
Waste receptacle: 4 liters  
Sound level: 48 dB  
Automatic overflow shutoff  
Aseptic filtered exhaust  
Large volume solids trap  
Automatic waste removal

#### *WATER MODULE*

1.5 liter reservoir  
Pressurized system: 30 to 40 lbs  
Aseptic non-retracting HPC coolant  
Manual handpiece flush system  
External quick connects for water and air



### FEATURES INCLUDE

- ♥ Wet/dry foot control
- ♥ Automatic controls for one, two and three handpieces
- ♥ Accomodates most brands of handpieces and fiber optic systems
- ♥ Three-way syringe
- ♥ Adjustable H.V.E. system
- ♥ Adjustable saliva ejector
- ♥ Asepsis designed chassis
- ♥ Recessed controllrs for easy cleaning
- ♥ Air filter moisture seperator
- ♥ Repair/Tool kit with operators manual
- ♥ Side maneuvering handles and high-grade casters for easy movement of the cart
- ♥ Optional arm-mounted trays and system cover

DNTLworks Equipment Corp.  
7300 South Tucson Way  
Centennial, CO 80112 USA  
Toll-free: 1-800-847-0694  
Phone: +303-693-1410  
Fax: +303-693-6189  
E-mail: [info@DNTLworks.com](mailto:info@DNTLworks.com)

ProCart II is a registered trademark of  
DNTLworks Equipment Corporation ©2011

