



## BETTER TECHNOLOGY BETTER VALUE



707 Executive Blvd. Unit D • Valley Cottage, NY 10989 Phone (845) 589-0210 • Fax (845) 589-0211 info@johnsonpromident.com www.deldent.com • www.johnsonpromident.com

# Better Airpolishing Technology

The **JetPolisher 2000™** (airpolishing system) and the **JetSonic 2000M™** (combination airpolishing and scaling system) feature Deldent's exclusive, patented Homogenous Stream Technology (H.S.T.)™.

#### The Homogenous Stream

**Technology (H.S.T.)** mixes air, water, and sodium bicarbonate powder together within a single nozzle to produce a homogenous stream, unlike other airpolishing systems in which there are separate nozzles for water and powder streams.

Benefits of the H.S.T. system include:

- Virtual elimination of nozzle clogging.
- A more efficient polishing stream, so only 35 psi of air pressure is needed.
- A gentler, softer polishing action because sodium bicarbonate particles begin to soften within the spray head
  - The JetPolisher 2000 and the JetSonic 2000M's polishing unit are harmless to titanium implant surfaces and are the ONLY polishers cleared by the FDA for use around titanium implants.\*

\* for cleaning around implants under dentures

#### Clinical Applications.



#### General prophylaxis

- Plaque & stain removal
- Cleaning abutments
- · Prior to fissure sealing



#### Periodontics

- Plaque & stain removal
- Flutings & furcations
- GTR root debridement



#### Orthodontics

- Around fixed appliances
- · Cleaning prior to bonding
- Post-treatment cleaning



**Implantology** 

- Implant maintenance
- Cleaning superstructures



- Non-clogging nozzle due to exclusive H.S.T. technology.
- Harmless to titanium implants.
- No installation required, completely portable, and can be moved easily from room to room (measures just 3.5" x 3.5" x 7.5").
- Connects directly to existing handpiece tubing (30-35 psi) via 2-3 hole or 4 hole connectors (Kavo, W&H, Sirona connectors also available).
- · Uses existing dental unit's foot pedal.
- Lightweight and easy on the wrist.
- Fully detachable autoclavable handpiece and water control knob to prevent cross-infection.
- · Large, easy-view powder reservoir.
- 2 year warranty.

\* for cleaning around implants under dentures

#519310	JetPolisher 2000 with 4-Hole Connector	#509120	Regular Handpiece
#519900	JetPolisher 2000 with Kavo-Type Connector	#509111	Booster Handpiece (for heavy stain)
#519500	JetPolisher 2000 with Sirona-Type Connector	#509550	Spray Head Tubes (bag of 50)
#519600	JetPolisher 2000 with W&H-Type Connector	#500004	Jetstream Mint Powder
			(4 x 9 oz. tubs per carton)

Recommended to clean tooth surfaces prior to bonding procedures, placing sealants, and bleaching.

### JetSonic 2000M™

With Cavitron-Type Scaler



The airpolishing portion of the JetSonic 2000M has all the features of the JetPolisher 2000 with the addition of warm water flow. The scaler side of this unique unit also offers advanced magnetorestrictive technology, including:

- Unique sensor technology recognizes and activates the selected handpiece (scaling or polishing).
- Advanced patented digital circuitry automatically tunes any "stack-type" insert, whether 25k or 30k.
- Advanced microchip feedback system maintains optimum scaling power for all settings.
- Power-Surge Mode<sup>™</sup> technology increases power by 15-20% with a "double-click" on the foot switch for removal of tenacious deposits. The pre-selected power setting returns automatically when the foot control is released.
- Large, easy-view powder reservoir (no need to stop to add powder).
- Convenient rinse mode available.
- · Water heating feature.
- 2 year warranty.

Comes with two scaling inserts, polishing powder, and sprayhead tubes to get you started.

#533000 #509120 #509111 #509550 #500004

JetSonic 2000M 115v Regular Handpiece Booster Handpiece (for heavy stain) Spray Head Tubes (bag of 50) Jetstream Mint Powder (4 x 9 oz. tubs per carton)

#### Piezo-Type Scaler


#### DelSonic<sup>™</sup> 2000



High-Tech Piezo Ultrasonic Scaler

Advantages of linear tip vibration, compared with elliptical vibration (Cavitron-Type):

- Smoother root surface
- More efficient tartar removal
- Easier amalgam overhang removal
- Sleek, compact, quiet, powerful, and portable.
- Linear tip vibration, for more efficient and more comfortable scaling.
- High frequency 30 kHz design provides greater patient comfort than the competition (less vibration against teeth).
- Cool running, high torque handpiece (minimal tip warmth).
- Fully autoclavable handpiece sleeves, so purchasing additional handpieces is not necessary.
- Range of tips available for all clinical needs.
- Tips are designed to introduce spray even into deep pockets.
- May be used with any independent fluid supply system dispensing purified water or medicaments.
- 2 year warranty (compare with competitive brands' one year warranty on handpiece).

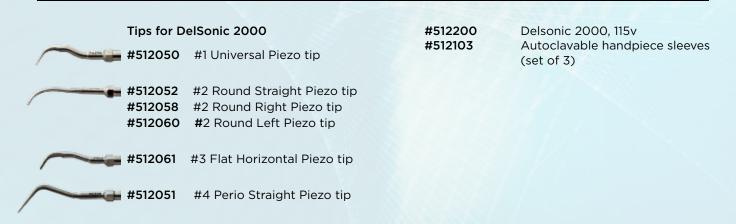

### Ultrasonic Inserts

A complete selection of ultrasonic scaling inserts for all stack-type (magnetostrictive) ultrasonic scaler units such as Cavitron\* scalers and others.

- Fully autoclavable
- Choice of styles available
- The better choice of inserts at highly competitive prices
- · Resin or metal handles
- Available in 25k and 30k series
- 3 month warranty

	DELDEI	NT			
Item Description	Item #	Best For	Dentsply		
RESIN-HANDLE INSERTS "DF-SERIES"					
DF #3 Flat 25K	508600	- Heavy to moderate supragingival deposits	TFI-3 Beavertail/#6121702		
DF #3 Flat 30K	508700	- Cement removal	TFI-3 Beavertail/#6374102		
DF #10 Universal 25K	508601	- Heavy to moderate supragingival deposits - Light to moderate	TFI-10 Straight/#6121708		
DF #10 Universal 30K	508701	subgingival deposits	TFI-10 Straight/#6374106		
DF #1000 3-Bend 25K	508607	- Light to moderate	TFI-1000 (Triple-Bend)/#6122001		
DF #1000 3-Bend 30K	508707	supragingival deposits	TFI-1000 (Triple-Bend)/#6374110		
Thinline Subging	gival Inse	rts (Low Power Only)			
DF #100 Straight 25K	508604	- Light subgingival	FSI-100/#80795, FSI-SLI 10/#80392		
DF #100 Straight 30K 508704 deposits	deposits	FSI-100/#80798, FSI-SLI 10/#80395			
ALL METAL INSERTS "DP-SERIES"					
DP #10 Universal 25K	508801	Heavy to moderate     supragingival deposits     Light to moderate	P10/#6120108		
DP #10 Universal 30K	508901	subgingival deposits.	-		
Subgingival Inserts (Low Power Only)					
DP #100 Straight 25K	508804	- Light subgingival	SLI-10S Straight/#79290		
DP #100 Straight 30K <b>508904</b> deposits		deposits	SLI-10S Straight/#79611		

### Ultrasonic Inserts

(Continued)

Longer plates vibrate at 25k frequency

Resin (plastic) handles have water spray that comes from inside the tip

Shorter plates vibrate at 30k frequency

Metal handles have external water spray

SIMILAR TO					
Hu-Freidy	Coltene/Whaledent	Parkell	Premier	Tony Riso	<b>Henry Schein</b>
#3 Streamline Beavertail/ #UI325K	-	#3 Beavertail 25K	-	-	-
#3 Streamline Beavertail/ #UI330K	-	#3 Beavertail 30K	-	-	-
Original Prophy Plus/ #UI1025K	#10/#US1025K	-	#1005300	IF 50 25K	#102-3077
Original Prophy Plus/ #UI1030K	#10/#US1030K	-	#1005310	IF 50 30K	#102-8413
#1000 Streamline Triple-Bend/ #UI100025K	#1000 Triple Bend/ #US100025K	-	#1005305	-	#102-7537
#1000 Streamline Triple-Bend/ #UI100030K	#1000 Triple Bend/ #UI100030K	-	#1005315	-	#102-3078
#100 Streamline/ #UI25K100S	#10 Slim/#US1025KSP	-	-	IF 100 25K	-
#100 Streamline/ #UI30K100S	#10 Slim/#US1030KSP	-	-	IF 100 30K	-
#10 After Five/ #UI25KSL10S	-	DUE-25	#1006803	P50 25K	#102-8168
#10 After Five/ #UI30KSL10S	-	DUE-30	#1006903	P50 30K	#102-4144
-	-	DPE-25	#1006801	P100 25K	#102-8404
-	-	DPE-30	#1006902	P100 30K	#102-3079

7

### Miniblaster™

Fully autoclavable clinical sandblaster designed to enhance adhesive dentistry procedures.



- Produces microretentive surfaces to enhance bonding to all dental materials.
- Index-finger controlled valve gives maximum ergonomic control and requires less operating pressure.
- Designed for left- and right-handed operators.
- Unique rachet system for simplified nozzle rotation and angulation;
   no need to waste time loosening and tightening collars.
- FDA approved for intra-oral use.
- Tungsten carbide nozzle for wear resistance.
- Optional nozzle extension with 360° rotation for difficult to reach areas.
- In-line high pressure stop-cock with built-in on/off toggle included.
- · Simple installation kit included.
- Can also use at 35psi on existing handpiece tubing or handpiece quick disconnect (Kavo, W&H, Sirona).
- 1 year warranty.

#505000	Miniblaster
#505130	Miniblaster Extension Nozzle
#505050	Aluminum Oxide Powder 50 micron
	(2 pound bottle)
#505090	Aluminum Oxide Powder 90 micron
	(2 pound bottle)

#### **Miniblaster Clinical Applications**

- Can be used for surface preparation prior to bonding to all dental materials, including: precious, semi-precious, and non-precious metals; porcelain; composite; and resin materials.
- Macro-retention such as undercuts should be employed where possible in conjunction with the microretention produced by the Miniblaster.
- Use 50 micron aluminum oxide for routine work and 90 micron for more aggressive work and cement removal.



Dust-Inn 2000™

Sleek, high-tech dust collector with an all-around see-through dome.

No more flying grit or dust!

- Self-contained filtration system (no large vacuum systems to set up).
- For use with Miniblaster or other clinical sandblasters and for other grinding procedures.
- Compact, lightweight and portable.
- Attractive design and easier to use than all other dust collectors.
- Quiet, efficient fan extractor system with easily exchangeable filters.
- Hygienic wipe-down design.
- Large, transparent, replaceable lid.
- Used abrasive collects beneath the perforated worktable for simple disposal - just lift and empty.
- · Ideal for chairside use in a dental office or dental laboratory.
- 1 year warranty.

#### **Miniblaster Clinical Applications (Cont'd)**

Increase retention of all bonding surfaces





BEFORE AFTER 400% increase of bond strength with sandblasting!

#502900 #502602

#502608

Dust-Inn 2000 110v Dust-Inn Replacement Dome with Handguards

Dust-Inn Replacement Filters (5 Pack)

Intraoral Surface Prep



Chairside convenience with handpiece tubing connections

Cement Removal



90 micron Aluminum Oxide for aggressive work

Orthodontic Preparations



-Improve retention of bands -Remove composite from detached brackets

### JetStream Powder

Fresh-Tasting Sodium
Bicarbonate Powder with a
Hint of Mint



- Convenient 2.2 pound boxes have 4 foil safety-sealed, snaplid plastic jars.
- No more pouring from a big clumsy jar!
- Cost effective, stay-dry snap lids protect powder.
- Suitable for use in Deldent's JetPolisher and JetSonic 2000M and in all standard prophy units.

#500004

Jetstream mint-flavored prophy powder (4x9 oz. tubs per carton)

### Aluminum Oxide

Available in 50 and 90 micron



- 50 micron for general etching and blasting.
- 90 micron for heavier work.
- Economical 2-pound bottles.
- For use with the Deldent Miniblaster and all clinical sandblasters.

#505050

Aluminum Oxide Powder 50 micron (2 pound bottle)

#505090

Aluminum Oxide Powder 90 micron

(2 pound bottle)

Cost effective & convenient packaging!

# Parts & Accessories

IETDOLICHED 2000 AND IETCONIC 2000

JETPOLISH	ER 2000 AND JETSONIC 2000M
#509120	Regular Handpiece
#509111	Booster Handpiece
#509550	Spray Head Tubes (bag of 50)
#500004	Jetstream Mint Powder (4 x 9 oz. tubs per carton)
DELSONIC	2000
#512103	Autoclavable Handpiece Sleeves (set of 3)
#524000	Stream Selector 2000
ULTRASON	IC INSERTS
#508006	O Rings (Package of 12)
MINIBLAST	ER
#505133	Miniblaster Holder
#505110	Miniblaster Nozzle Replacement Kit
#505119	Miniblaster Powder Bottle Complete
#505121	Miniblaster Pressure Stop-Cock
#505130	Miniblaster Extension Nozzle
#505123	Miniblaster Quick Disconnect Hook-Up Kit w/Metal 1/4 'T'
#505124	Miniblaster Quick Disconnect Hook-Up Kit w/Metal ¾ 'T'
#505126	Miniblaster Quick Disconnect Hook-Up Kit w/Plastic 'T'
#505128	Miniblaster Male Quick Disconnect
#505129	Miniblaster Standard Hook-Up Kit
#505116	Miniblaster Pick-Up Tube
#505111	Miniblaster Adaptor for 4 Hole High Speed Connector
#505112	Miniblaster Adaptor for 2-3 Hole High Speed Connector
#505113	Miniblaster Adaptor for Kavo-Style High Speed Connector
#505114	Miniblaster Adaptor for Sirona-Style High Speed Connector
#505117	Miniblaster Adaptor for W&H-Style High Speed Connector
#505050	Aluminum Oxide Powder 50 Micron (2 pound bottle)
#505090	Aluminum Oxide Powder 90 Micron (2 pound bottle)
DUST-INN :	
#502602	Dust-Inn 2000 Replacement Dome with Handguards
#502608	Dust-Inn 2000 Replacement Filters (5 pack)
#157009	Dust-Inn 2000 115v Transformer
MINIDI ATE	<b>R</b> (Used to increase bond strength after sandblasting)
#504000 #504009	Miniplater Kit  Miniplater Tip Plating Solution (10cc)
#504009	Miniplater Tip Plating Solution (10cc)  Miniplater Applicator Tips (20 pcs)
#504007	
	Miniplator Sot of Patteries (2 pcs)
#504003	Miniplater Set of Batteries (3 pcs)
#504012	Cotton Bud Holder

### DiaGlaze



- Creates high lustre and natural enamel finish when polishing porcelain, gold base metals, and composite restorations.
- Advanced natural diamond particle technology for superb natural appearance and lustre.
- 3 controlled grades of diamond paste provide unmatched esthetics.
  - Extra fine grade contains ultrafine diamond particles of O-1/2 micron in size to achieve surface smoothness never possible before.
- Food grade components have high viscosity to prevent splattering and have a fresh mint scent.
- Optimum results are achieved with a progressive method that begins with the regular grade, followed by fine, and finished with extra fine.

#520005	Diaglaze 2gm EXTRA FINE
#520006	Diaglaze 2gm FINE
#520007	Diaglaze 2gm REGULAR
#520004	Diaglaze 3-Pack of 2gm Syringes:
	REGULAR, FINE, and
	EXTRA FINE





When it comes to dental polishing, etching, blasting, and cleaning, our high quality Deldent equipment and supplies provide the most attractive features at the best prices anywhere. With production and engineering capabilities in the U.S., Japan, China, and Israel, Johnson-Promident is now your onestop shop for all small equipment requirements, including operatory and lab handpieces, ultrasonic and sonic scalers and inserts, curing lights, air polishers, sand blasters, and many other essential related accessories.