

Fill, pack and go

5mm



+



one bulk fill
placement up
to 5mm

20 sec. LED
light cure



“

I have been using **aura bulk fill** for three years now. The resin comes out smooth with no pull back and cures extremely strong. The polish is amazing and my posterior composites have never been better. The characteristics are ideal for any restoration, primary or permanent teeth. I have recommended this product to many associates and I plan to continue using the **aura bulk fill** as my go to composite.

Dr Liel Grinbaum, DMD

SDI

Your Smile. Our Vision.
www.sdi.com.au

aura bulk fill
ultra universal restorative material

aura bulk fill

One Shade. One Fill.

Speedy and simple posterior composite restorations.

- Natural-looking finish with just one universal shade
- High depth of cure without the need for a capping layer
- Easy placement in a matter of seconds
- The waxy consistency is packable and non-sticky to optimize handling



The Dental Advisor, Vol. 33,
No. 07, September 2016



2.21%

Low volumetric shrinkage

aura bulk fill delivers reduced stress to provide long-lasting restorations.

115 MPa

Excellent flexural strength

aura bulk fill resists major deformations without fracturing.

375 MPa

High compressive strength

aura bulk fill resists the high forces of compression in the mouth.


Depth of Cure

The opacity of **aura bulk fill** is a function of the difference in the refractive index of its two components, the filler and the resin.

The curing process alters the refractive index of the resin marginally, to match the refractive index of the filler.

This lowers the opacity temporarily, allowing deeper light penetration for a high depth of cure.

Upon curing, the indices move apart again to arrive at the ideal opacity for **aura bulk fill's** brilliant chameleon-like finish.

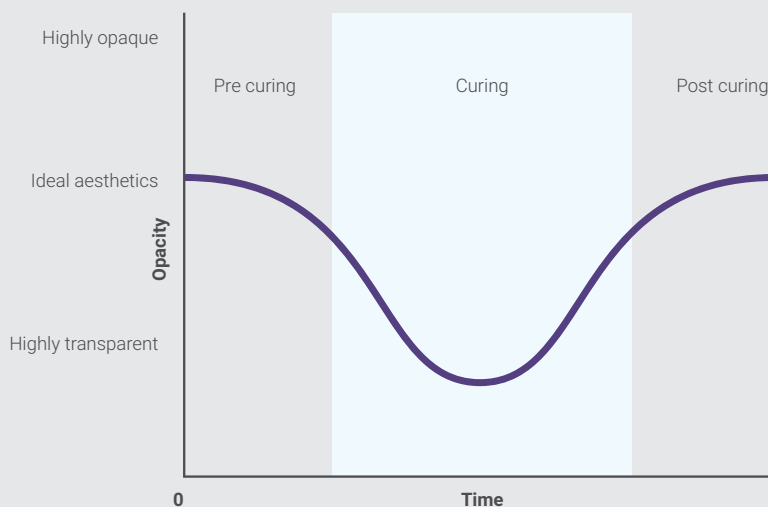

BEFORE



AFTER

— Dr. Bob Lowe (USA)



BEFORE



AFTER

— Dr. Gregor Thomas (Germany)

// 8565012 Bulk Fill 20 x 0.25g Complet Refills // 8560012 Bulk Fill 1 x 4g Syringe Refills



Your Smile. Our Vision.
www.sdi.com.au

Australia 1 800 337 003
Austria 0800 0225 5734
Brazil 0800 770 1735
France 0800 0225 5734
Germany 0800 1005759
Ireland 0800 0225 5734

Italy 00800 0225 5734
New Zealand 0800 734 034
Spain 00800 0225 5734
United Kingdom 0800 0225 5734
USA & Canada 1 800 228 5166

