


# Filtek™ Bulk Fill

Flowable Restorative

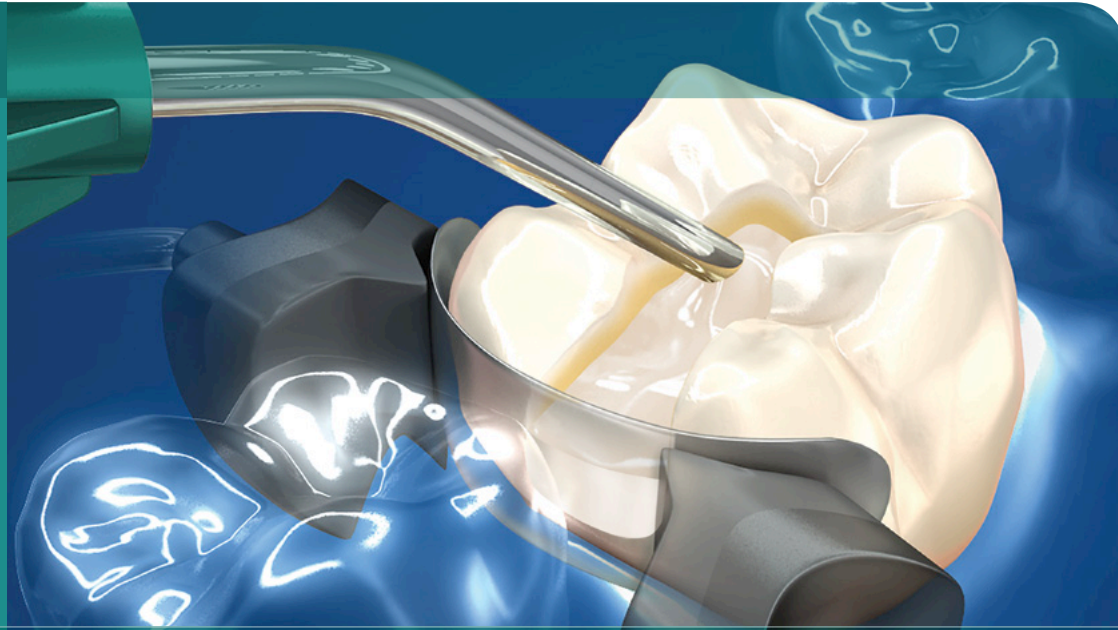

The bulk fill flowable  
that's fast and strong.

# Filtek™ Bulk Fill

## Flowable Restorative

Reduce procedure time  
and increase your  
confidence with this  
bulk fill flowable.

With the trend toward faster and easier procedures, you have more choices than ever for a flowable composite. So how do you choose the one that stands up to your needs and expectations for efficiency, strength and all-around reliable performance?

Rely on the reputation of a leading innovator of dental restorative solutions ... 3M ESPE.




**Excellent**  
handling ✓


## Excellent as a liner

Flowable composites are designed to adapt easily, reducing the chance for voids and providing a good marginal seal to reduce the risk of post-operative sensitivity. These benefits are likely why 75% of dentists surveyed use a flowable composite as a liner.<sup>1</sup>

**Very**  
**Low** ↓  
polymerization stress  
compared to traditional flowables  
placed and cured in 2mm increments<sup>2</sup>

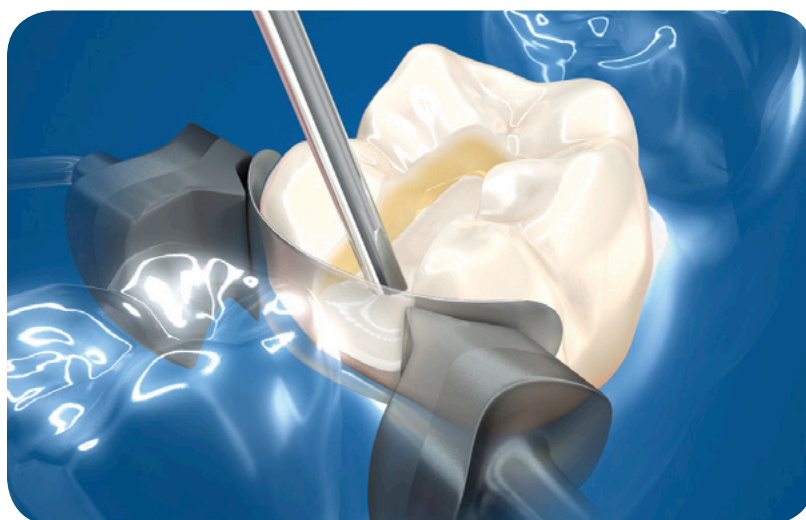

**Easy**  
cavity  
adaptation


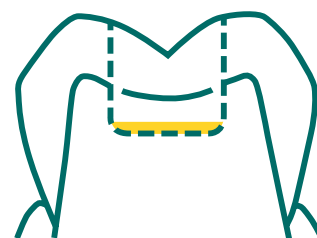


1. Clinicians Report, July 2013 Vol. 6, Issue 7.

2. Compared to Grandio Flow, TPH3 Flow, Premise Flowable, Revolution Formula 2, EsthetX Flow and Tetric EvoFlow. Source: 3M ESPE internal data.


## Excellent as a time-saving bulk fill base in deep restorations

Filtek Bulk Fill Flowable Restorative is formulated for low stress and an increased depth of cure, so it can be placed with confidence in one convenient 4 mm increment. The restoration is then completed with a minimum 2 mm thickness of a posterior or universal composite such as Filtek™ Supreme Ultra Universal Restorative.

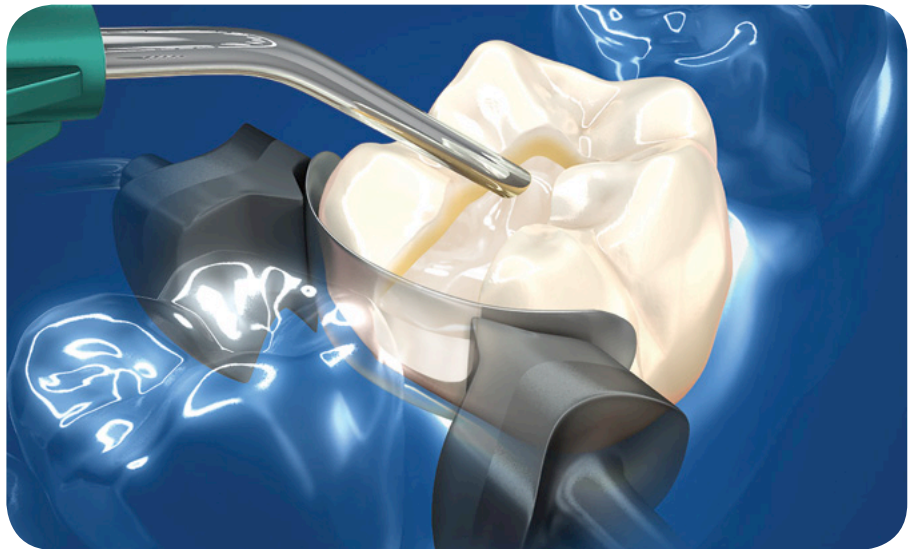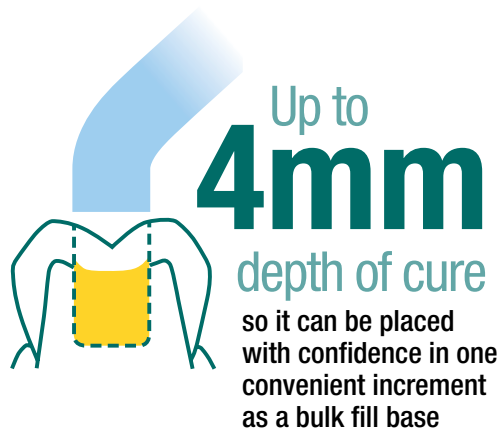

## Less shrinkage stress

than more  
highly filled,  
sculptable composites<sup>4</sup>

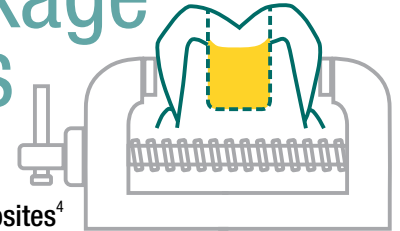

## Highly Radiopaque

compared to traditional  
flowables and other  
bulk fill flowables<sup>3</sup>



## Excellent for surface restorations

High strength and excellent wear resistance make it a versatile solution for surface restorations including small, non-stress bearing occlusal restorations.

Over  
**3X** the wear  
resistance  
of the leading bulk fill flowable<sup>5</sup>

3. 3M ESPE internal data

4. Source: J Dent Res 92 (Spec Issue A):2439 El-Damanhoury.

5. In-vitro data on file. Source: J Dent Res 92 (Spec Iss A):abs. #2441 K. SAWLANI, P. BECK, L.C. RAMP, D. CAKIR-USTUN, and J. BURGESS, University of Alabama at Birmingham, Birmingham, AL.

6. 3M ESPE internal data. Compared to Surefill SDR Flow.

# Choose Filtek™ Bulk Fill Flowable Restorative for fast and simple posterior restorations.



**1**  
Prepare the tooth and apply adhesive.



**2**  
Dispense up to 4 mm of Filtek™ Bulk Fill Flowable Restorative into the prepared cavity and cure.



**3**  
Place a 2 mm (minimum) capping layer of an esthetic composite, such as a Filtek™ universal restorative, and cure.

Clinical photos courtesy of Dr. Reich, Seefeld, Germany

## Ordering Information



Filtek™ Bulk Fill Flowable Restorative  
Refill—Syringe



Filtek™ Bulk Fill Flowable Restorative  
Refill—Capsule

### Filtek™ Bulk Fill Flowable Restorative

Refill	Item #	Item #
SHADE	SYRINGE	CAPSULE
A1	4862A1	4861A1
A2	4862A2	4861A2
A3	4862A3	4861A3
Universal	4862U	4861U
Item #	Product Description	
3740T	Dispensing Tips	
	20 – 19-gauge dispensing tips	

**SYRINGE SHADE REFILLS** include 2 – 2g syringes; 20 dispensing tips; Technique Guide; Instructions for Use.

**CAPSULE SHADE REFILLS** include 15 – 0.2g capsules; Technique Guide; Instructions for Use.

Customer Care Center: 1-800-634-2249 [www.3MESPE.com/FiltekBulkFill](http://www.3MESPE.com/FiltekBulkFill)

★★★★★ Review these products at [www.3MESPE.com/review](http://www.3MESPE.com/review)

**3M ESPE**

#### Dental Products

2510 Conway Avenue  
St. Paul, MN 55144-1000 USA  
1-800-634-2249

#### 3M Canada

Post Office Box 5757  
London, Ontario N6A 4T1  
Canada  
1-888-363-3685

Please recycle. Printed in U.S.A.  
70-xxxx-xxxx-x

3M, ESPE and Filtek are trademarks of 3M or 3M Deutschland GmbH. Used under license in Canada. © 3M 2014. All rights reserved. EsthetX, Grandio, Premise, Revolution Formula 2, SureFil SDR, Tetric EvoFlow, TPH and Venus are not trademarks of 3M.