


Specifications

- Power adapter: 90-250V AC, 50-60 Hz
- Class: Class II (ICE601-1)
- Power consumption: less than 10W
- Wavelength range: 430-490 nm
- Light Intensity: 1,250 mW / cm²
- Dimensions: 7"L x 1"W x 1"H
- Weight: 160 g
- Warranty: 2 years

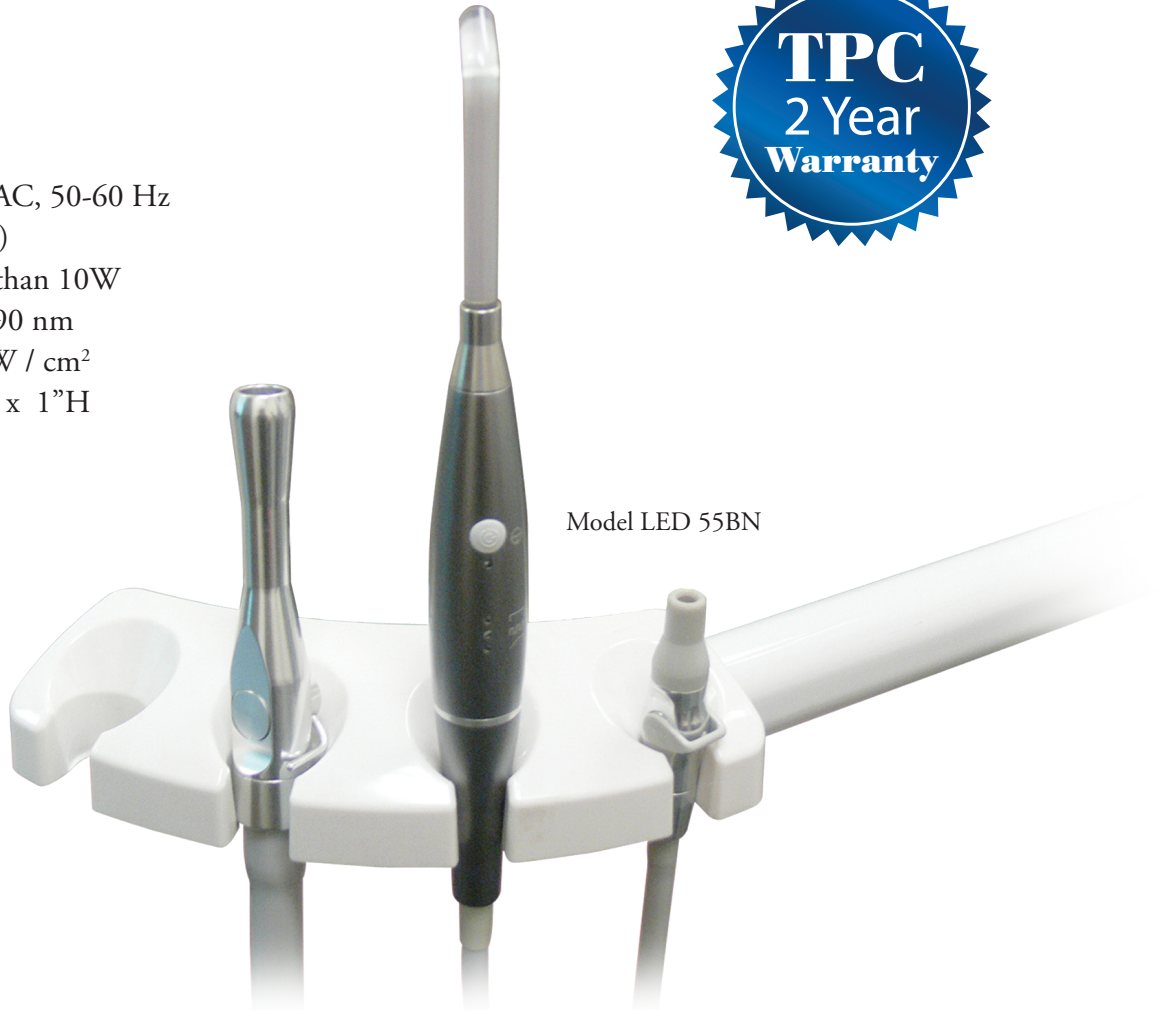




Model LED 55BN

Curing Mode Light Indicators

LED color	Working Mode
Green	Fast curing: Time setting 20 seconds (every 10 second, it sounds once)
Red	Pulse curing: interrupt radiating (every 10 seconds, it sounds once)
Orange	Step curing: soft start for 10 second from 0 to full power, then full power during the later 10 second

Curing Status Light Indicators

LED color	Status Indication
Green	Normal
Red Flashing	LED curing light handpiece is over temperature

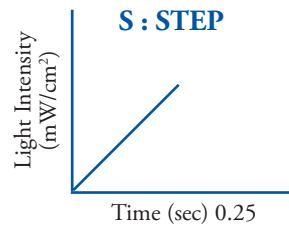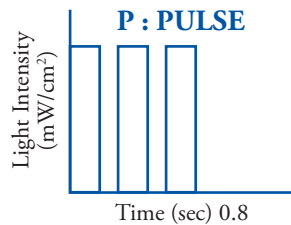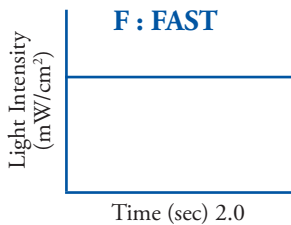

Shipping Dimensions: 10"L x 9"W x 5"H / 3 lbs.