



AdvanceCAM Intraoral Camera

Corded Direct USB Camera System

New Standards for Image Quality and Value...

Model #AIC5899



Model #AIC5111



Simplicity • Portability • Versatility • Clarity



Introducing the AdvanceCAM USB Direct: TPC continues our long tradition of providing dental professionals with quality imaging equipment with our brand new, lightweight intraoral camera. Combining the latest video technologies, the AdvanceCAM USB Direct is an exceptional tool for patient education, higher case acceptance, and improved patient communication.

Portability

- High-speed direct USB connection with no docking station required
- Features a lightweight handpiece and reduced cabling
- Disconnects effortlessly for quick transfer to another operator

Specifications

Imager: Sony 1/4" CCD
Capture Size: 640 x 480
Light Source: 6 white LEDs
Power: USB 5 VDC
Frame Rate: 30 fps max
PC Connection: USB 2.0 Port

Technology

- Virtually eliminates image distortion by utilizing true optics
- Optional Mini LCD screen
- State-of-the-art 4-LED lighting system

System Requirements

Intel Pentium 4.15 GHz or Higher
512 MB RAM or higher
Windows XP or higher
USB 2.0 port

Shipping Dimensions: 12"L x 9"W x 4"H / 2 lbs.
(10/cs.) 27"L x 12"W x 8"H / 13 lbs.





AdvanceCAM Intraoral Camera

Mini SD Card Camera System

New Standards for Image Quality and Value...


Model #AIC5855+AIC5901+AIC5111



Simplicity • Portability • Versatility • Clarity

Model #AIC5855

Mini SD memory card (included)

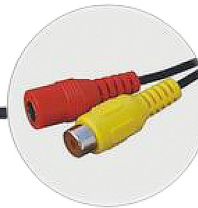

Model #AIC500

Introducing the AdvanceCAM Mini SD: TPC continues our long tradition of providing dental professionals with quality imaging equipment with our brand new, lightweight intraoral camera. Combining the latest memory card technologies, the AdvanceCAM Mini SD is an exceptional tool for patient education, higher case acceptance, improved patient communication and maximum portability.

Maximum Portability

- Mini SD memory storage
- Optional battery pack for complete cordless operation
- Features a lightweight handpiece and no cabling
- Connects directly to LCD monitor via composite video cable

Advanced Technology

- State-of-the-art 4-LED lighting system
- Mini SD card technology
- Optional Mini LCD

Specifications

Imager:	Sony 1/4" CCD
Capture Size:	640 x 480
Light Source:	6 white LEDs
Power:	5 VDC
Supported Mini SD:	Up to 16 GB

Shipping Dimensions: 12"L x 9"W x 4"H / 2 lbs.
(10/cs.) 27"L x 12"W x 8"H / 13 lbs.

